



your trusted partner, roadmapping  
the future of behavioral research



# CATALOGUE 2015




**ugo basile**

TRANSFORMING IDEAS INTO INSTRUMENTS

Phone +39 0332 744574  
sales@ugobasile.com

[www.ugobasile.com](http://www.ugobasile.com)

# what's new

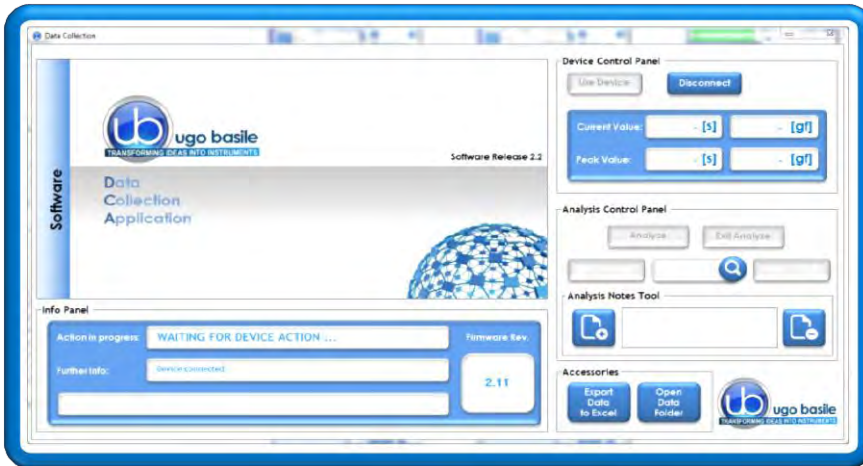

**new DCA Software  
for PAM, eVF,  
GSM**

**new conditioning  
cage system, new  
software**

**one touch-screen unit to  
manage different tests**

- fear conditioning
- active avoidance
- passive avoidance
- learned helplessness



**mouse multi-maze system:  
T, Y and radial maze with a  
single equipment**



**your trusted partner...**



**roadmapping the future  
of behavioral research!**



## PAIN AND INFLAMMATION

|   |              |   |    |
|---|--------------|---|----|
| ● | 37140        | Plethysmometer .....  | 1  |
| ● | 37215        | Analgesy-Meter .....  | 3  |
| ● | 35100        | Hot/Cold Plate .....  | 5  |
| ● | 37370        | Plantar Test (Hargreaves's Apparatus) .....                       | 7  |
| ● | 37360        | Tail-Flick Unit .....   | 9  |
| ● | 37300        | I.R. Heat-Flux Radiometer .....                                   | 11 |
| ● | 37450        | Dynamic Plantar Aesthesiometer .....                              | 13 |
| ● | <b>38500</b> | P.A.M. Pressure Application Measurement ..... <b>NEW SOFTWARE</b> | 15 |
| ● | <b>38450</b> | e-VF Electronic Von Frey ..... <b>NEW SOFTWARE</b>                | 17 |
| ● | 37450-275    | Von Frey Hairs .....  | 19 |
| ● | 31300        | Orofacial Stimulation Test (Hargreaves Method) .....              | 21 |
| ● | 37100        | Durham Animal Holders (for trigeminal stimulation) .....          | 23 |

## MOTORY COORDINATION, GRIP STRENGTH, ACTIVITY

|   |              |   |    |
|---|--------------|---|----|
| ● | <b>47200</b> | Grip Strength Meter ..... <b>NEW SOFTWARE</b> | 25 |
| ● | 47420        | Activity Cage .....                           | 27 |
| ● | 47700        | Rat Rota-Rod .....                            | 29 |
| ● | 47600        | Mouse Rota-Rod .....                          | 31 |
| ● | 43000        | Rotometer .....                               | 33 |
| ● | 6650         | Hole Board .....                              | 35 |
| ● | 1800         | Activity Wheels for Rats .....                | 37 |
| ● | 1850         | Activity Wheels for Mice .....                | 37 |

## VENTILATORS and GAS ANESTHESIA

|   |         |  |    |
|---|---------|--|----|
| ● | 28025   | Mouse Ventilator .....                                     | 39 |
| ● | 7025    | Rodent Ventilator .....                                    | 41 |
| ● | 6025    | Cat-Rabbit Ventilator .....                                | 43 |
| ● | 7020    | Bronchospasm Transducer .....                              | 45 |
| ● | 21000   | Gas Anaesthesia Systems .....                              | 47 |
| ● | 7900/10 | Anaesthetizing Boxes and Gas Anaesthesia Accessories ..... | 50 |





## BEHAVIOUR, CONDITIONING, REWARD

|   |                  |   |                            |    |
|---|------------------|---|----------------------------|----|
| • | <b>40500-001</b> | Touch-Screen Controller .....                         | <b>NEW FUNCTIONS</b> ..... | 51 |
| • | <b>46000</b>     | Fear Conditioning Systems .....                       | <b>NEW SOFTWARE</b> .....  | 53 |
| • | <b>40532</b>     | Active Avoidance Set-Up (Shuttle-Cage) for Rats ..... | <b>NEW MODEL!</b> .....    | 55 |
| • | <b>40533</b>     | Active Avoidance Set-Up (Shuttle-Cage) for Mice ..... | <b>NEW MODEL!</b> .....    | 55 |
| • | <b>40552</b>     | Passive Avoidance (Step-through) for Rats .....       | <b>NEW MODEL!</b> .....    | 57 |
| • | <b>40553</b>     | Passive Avoidance (Step-through) for Mice .....       | <b>NEW MODEL!</b> .....    | 57 |
| • | <b>40570</b>     | Passive Avoidance (Step-down) for Mice .....          | <b>NEW MODEL!</b> .....    | 59 |
| • | <b>47502</b>     | Learned Helplessness for Rats .....                   | <b>NEW MODEL!</b> .....    | 61 |
| • | <b>47503</b>     | Learned Helplessness for Mice .....                   | <b>NEW MODEL!</b> .....    | 61 |
| • | <b>42552</b>     | Conditioned Place Preference (CPP) for Rats .....     |                            | 63 |
| • | <b>42553</b>     | Conditioned Place Preference (CPP) for Mice .....     |                            | 63 |
| • | <b>45100</b>     | Lickometer – Vogel Test for Rats .....                |                            | 65 |
| • | <b>45150</b>     | Lickometer – Vogel Test for Mice .....                |                            | 65 |

## BEHAVIOUR, MAZES, TRACKING

|   |                |   |    |
|---|----------------|---|----|
| • | <b>46503</b>   | Sociability Apparatus (3-chambered social test) .....   | 67 |
| • | <b>40100</b>   | Hydraulic Atlantis Platforms for Morris water maze .....  | 69 |
| • | <b>60000</b>   | ANY-maze VideoTracking System, <i>see also</i> <a href="http://www.ub.anymaze.com">www.ub.anymaze.com</a> ..... | 71 |
| • | <b>40125</b>   | Morris Water Maze .....   | 73 |
| • | <b>40132/3</b> | T- Maze for Rats/Mice .....   | 73 |
| • | <b>40142/3</b> | Plus-Maze for Rats/Mice .....   | 73 |
| • | <b>40153</b>   | Radial-Maze for Maze .....  | 73 |
| • | <b>40162/3</b> | Elevated Zero-Maze for Rats/Mice .....  | 73 |
| • | <b>40172/3</b> | Y-Maze for Rats/Mice .....  | 73 |
| • | <b>40192/3</b> | Barnes Maze for Rats/Mice .....   | 73 |
| • | <b>40150</b>   | Open Fields for Videotracking, for Rats/Mice .....  | 73 |
| • | <b>41500</b>   | Modular Multi-Maze System for Mice (Y-, T-, Radial Maze) ... <b>NEW!</b> .....                                  | 75 |

## TISSUE BATHS, TRANSDUCERS, STIMULATORS

|   |             |   |    |
|---|-------------|---|----|
| • | <b>4000</b> | Isolated Organ Baths single chamber / 2 chambers / 4 chambers ..... | 77 |
| • | <b>3165</b> | Multiplexing Pulse Booster .....                                    | 79 |



|   |       |  |    |
|---|-------|--|----|
| • | 14900 | Superfusion System .....                                 | 81 |
| • | 7003  | Isometric Force Transducers DY1 .....                    | 83 |
| • | 7004  | Isometric Force Transducers DY2 .....                    | 83 |
| • | 7005  | Isometric Force Transducers DY3 .....                    | 83 |
| • | 7010  | Isometric Force Transducers DY0 .....                    | 83 |
| • | 7006  | Isotonic Transducer .....                                | 85 |
| • | 17308 | DataCapsule-Evo Digital Recorder ..... <b>NEW MODEL!</b> | 87 |

## MISCELLANEOUS, ECT, LESION MAKING DEVICE

|   |       |   |    |
|---|-------|---|----|
| • | 57800 | ECT Unit .....                          | 89 |
| • | 53500 | Lesion Making Device .....              | 91 |
| • | 51600 | Stoelting Stereotaxic Instruments ..... | 93 |
| • | 5000  | KDS Infusion Pumps .....                | 95 |

## BLOOD PRESSURE, VITAL FUNCTIONS

|   |         |   |     |
|---|---------|---|-----|
| • | 58500   | Blood Pressure Recorder, non Invasive, for Rats ..... | 97  |
| • | 58600   | Blood Pressure Recorder, non Invasive, for Mice ..... | 97  |
| • | 17844   | Blood Pressure Transducer Invasive .....              | 99  |
| • | MouseOx | Starr MouseOx Pulse Oximeter .....                    | 101 |

## METABOLISM, FEEDING BEHAVIOUR

|   |       |  |     |
|---|-------|--|-----|
| • | 41800 | Metabolic Cages and Feeding Analyser ..... | 103 |
| • | 47552 | Rat Feeding & Activity Analyser .....      | 105 |
| • | 47553 | Mouse Feeding & Activity Analyser .....    | 105 |

## MUROMACHI MICROWAVE FIXATION

|   |        |   |     |
|---|--------|---|-----|
| • | MMW-05 | Muromachi Microwave Brain Fixation System ..... | 107 |
|---|--------|---|-----|

## BROWSE PRODUCTS BY APPLICATION (each caption is a link to the related web page)

### NEUROPATHIC PAIN, HYPERALGESIA, INFLAMMATION,

Analgesy-Meter  
Hot/Cold Plate  
P.A.M. Pressure Application Measurement  
Plantar Test (Hargreaves Apparatus)  
Plethysmometer  
Tail-Flick Unit  
Orofacial Stimulation Test (Fehrenbacher, Henry, Hargreaves method)

### ALLODYNIA, HYPERSENSITIVITY, SOMATOSENSATION

Dynamic Plantar Aesthesiometer  
Von Frey Hairs (with grid)  
Hot/Cold Plate  
P.A.M. Pressure Application Measurement  
e-VF Electronic Von Frey  
Durham Animal Holders

### MOTOR FUNCTION, PARKINSON'S, STRENGTH, EXERCISE

Mouse / Rat Rota-Rod  
Rotometer  
Activity Cage  
Isolated Organ Baths  
Mouse Open Field, 44cm  
Grip-Strength Meter (mice and rats)  
Hole Board

### VENTILATION, ANESTHESIA, SURGICAL MONITORING

Blood Pressure Recorder, non Invasive  
Blood Pressure Transducer, Invasive  
Cat/Rabbit Ventilator  
Gas Anesthesia System  
Mouse Ventilator  
Pulse Oximeter for Mice & Rats

### MEMORY, LEARNING, ALZHEIMER

NEW Fear Conditioning System  
Atlantis Platforms for Water Maze  
Passive Avoidance - step down NEW MODEL  
Any-maze Video-Tracking Software  
Water Maze Pool  
Barnes Maze  
Mouse Open Field, 44cm  
Passive Avoidance - Step through - New Model  
Multi-Maze System for Mouse  
Active Avoidance Set-Up (Shuttle-Box)

### ADDICTION AND REWARD, SOCIAL BEHAVIOUR AND AUTISM

Sociability Apparatus (3-chambered social test)  
NEW Place Preference (CPP)  
KDS Infusion Pumps  
Any-maze Video-Tracking Software  
Lickometer Vogel Test  
Mouse Open Field, 44cm

### ANXIETY, DEPRESSION, FEAR, STRESS

Learned Helplessness  
Activity Cage  
Any-maze Video-Tracking Software  
Lickometer Vogel Test  
Elevated Plus Maze  
Mouse Open Field, 44cm  
Elevated Zero-Maze

### ELECTROLYTIC LESIONS AND INFUSION

Lesion Making Device  
Stoelting Stereotaxic Instrument  
KDS Infusion Pumps  
NEW DataCapsule-Evo Digital Recorder

### EPILEPSY, SEIZURES, CONVULSIONS

Rotometer  
ECT Unit

### BRAIN CHEMISTRY, PHOSPHORYLATION,

Microwave Brain Fixation NEW 5KW

# Plethysmometer

Cat. No. 37140

## General

In research on rheumatoid arthritis, the central development of oedema, and its modifications by pharmacological processes, it has proved of great value to measure inflammatory processes in the rat paw.

Our **Plethysmometer 37140** displays the exact paw volume on the graphic LCD read-out. Small differences are detected by a transducer of original design.

The 37140 is provided with a pedal holding-command which freezes the reading, enabling the operator to concentrate its attention to the paw dipping.

The paw volume is shown on the multifunction graphic display in four digits, with 0.01 ml resolution. A zero key is provided to zero the meter before each measurement.



Now supplied with  
measuring cell for  
both RAT & MOUSE  
paw!!

**FOR ACCURATE  
MEASUREMENT OF:**

- RAT paw oedema
- MOUSE paw oedema



## MICROPROCESSOR Controlled Instrument. Main Features:

- Computer compatibility : direct connection to PC (via the 52050 Software included)
- Read-out : multifunction graphic display
- Print-out : by optional thermal MiniPrinters 57145

## Volume Measuring Water Cell

The measuring cell consists of two vertical interconnected Perspex tubes; the animal paw is dipped in the larger tube (1.8cm diam) to measure water displacement. A tube of smaller diameter (1.3cm) is also included for measuring the mouse paw.

The smaller diameter side tube contains the transducer which measures the conductance between two vertical wire electrodes; conductance is linearly proportional to the water level, hence to the displaced volume.

## Data Acquisition

The 37140 Plethysmometer is microprocessor controlled, featuring direct PC output. Internally stored data can be routed to the PC serial (RS232) or USB port (via adaptor).

Communication is managed by the dedicated Software Cat. 52050-02, a Windows® based Data Acquisition Software Package.

The 52050-02 enables the data storage into individual files, ready to be easily managed by most statistical analysis packages available on the market.

**Ask for details!**

### Also Available

- 37140-25** **Plethysmometer**, complete with water cell diam. **2.5cm** & standard accessories
- 37140-35** **Plethysmometer**, complete with water cell diam. **3.5cm** & standard accessories

### Other Available Water Cells

- 7157** Special Water Cell, diam. 2.5cm, complete with Transducer 7153-L
- 7159** Special Water Cell, diam. 3.5 cm, complete with Transducer 7153-L

### Optional

- 57145** Thermal Mini-Printer
- 37400-305** Thermal Paper Roll for 57145

### Specifications

|                   |  |
|-------------------|--|
| Power Requirement | Universal input 85-264 VAC, 50-60Hz, 40 W max.                               |
| Data Read-out     | multifunction graphic display  |
| Data Format       | 4 digits (2 integers, 2 decimals)  |
| Resolution        | 0.01 ml  |
| Commands          | via soft-buttons   |
| Connection to PC  | direct connection to PC USB port, via serial cable and serial to USB adaptor |
| Data Print-Out    | via the optional MiniPrinter 57145   |

### Physical

|                    |                |
|--------------------|----------------|
| Weight             | 4.8 Kg         |
| Shipping Weight    | 8.1 Kg approx. |
| Shipping Dimension | 67x42x53cm     |

## Bibliography

- D. Piomelli et alia: "Anandamide suppresses pain initiation through a pe-ripheral endocannabinoid mechanism". Nature NSC, 2010
- F. Vincenzi et alia: "A2A Adenosine Receptors Are Differentially Modulated by Pharmacological Treatments in Rheumatoid Arthritis Patients and Their Stimulation Ameliorates Adjuvant-Induced Arthritis in Rats" PLoS ONE 8(1): e54195, 2013
- T. Bertaim et alia: "Dose and Administration Schedule Effect of Tiludro-nate on Joint Damage in the Model of Complete Freund Adjuvant Induced Monoarthritis in Rats" Open Journal of Rheumatology and Autoimmune Diseases 3: 18-25, 2013
- E. Borbély et alia: "Role of Tachykinin 1 and 4 Gene-Derived Neuropeptides and the Neurokinin 1 Receptor in Adjuvant-Induced Chronic Arthritis of the Mouse" PLoS ONE 8(4): e61684, 2013
- R. Korhonen et alia: "Attenuation of TNF production and experimentally induced inflammation by PDE4 inhibitor rolipram is mediated by MAPK phosphatase-1" Br. J. Pharmacol. 169 (7): 1525-1536, 2013
- A. Finlay et alia: "Sphingosine 1-Phosphate Mediates Hyperalgesia via a Neutrophil-Dependent Mechanism" PLoS ONE 8(1): e55255, 2013

## Ordering Information

- 37140** **PLETHYSMOMETER**, standard package including:-
- 7141** Electronic Block
- 7152-S** Standard Water Cell, diam. 1.8cm, including Mouse paw tube **7186**, diam. 1.3cm
- 7153-L** Conductance Transducer
- 7140-154** Water Reservoir
- 7155** Calibration Probes (0.1, 0.2, 0.5, 1, 2, 4ml)
- 7160** Wetting Compound, 100ml bottle
- 7165** Connection tube (cell-reservoir & drain vessel)
- 37215-303** "Hold" Pedal Switch
- 52050-02** CUB Dedicated Software (on USB drive)
- 37140-302** Instruction Manual (on USB drive)
- 52010-320** USB to serial port converter
- 52010-322** Connecting cable 9 to 9 pin
- 4210** Three Claw Stand, 10mm diam. Upright, complete with **4003** Open Side Boss-Head
- E-WP 008** Mains Cord



# Analgesy-Meter

Cat. No. 37215

## General

The 37215 is the up to date version of the classical 7200 paw pressure test which, since 1965, in a number of academic and industrial laboratories, is helping to perform a rapid precise screening of analgesic drugs

The force is applied to the animal's paw, which is placed on a small plinth under a cone-shaped pusher with a rounded tip.

The 37215 features a low voltage synchronous motor and conforms the CE rules.

The operator depresses a pedal switch to start the mechanism which exerts the force.

When the rat struggles, the operator releases the pedal and reads off the scale the force at which the animal felt pain.

The standard 37215 can be conveniently used with mice. However, a dedicated model is also available, when lower pressure (50%) is desirable, model **37216**, which includes a special chisel-shaped pusher



- **Randall-Selitto Paw Pressure Test**

- **Rapid Precise Screening of Analgesic Drugs**

## Main Features

- Same instrument, three force ranges (from 0 to 250, 500, 750 g)
- Simple and reliable: no calibration needed!
- Classic method since the 1960s : hundreds of papers published!

## Principle of Operation

The force applied to the paw by the plinth, increases at a constant rate, thus enabling perfect reproducible measurements to be made. The motor stops immediately the pedal is released.

After each test the slide should be returned to its starting point by lifting it and pushing it to the left.

The force is measured on the scale calibrated in 10-gram steps, by a pointer riveted to the slide. The scale can be multiplied by 2 or 3, by placing on the slide one or two discs provided with the standard package.


## Data Acquisition

The 37215 incorporates an optical switch, to enable the connection to the optional Multifunction Printer Cat. 2600, via the dedicated cable 2610-D.

The Multifunction Printer is a microprocessor controlled device, designed to acquire data from 6 independent channels (each Analgesy-Meter requires 1 channel).

The data, stored in the 2600 internal memory and shown on its graphic display, can be printed out in real time and/or routed to the PC via the 52050 CUB Software Package included with the 2600.

### Ask for details!

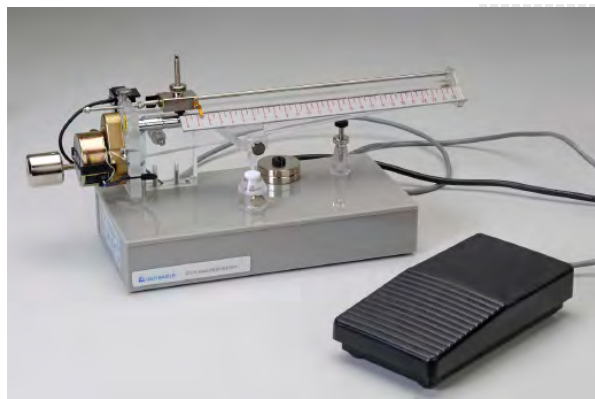

## Specifications:

### Power

- Requirements : 115 or 230V, 50/60Hz, 15W max.  
 Start / Stop : by pedal switch  
 Force Range : 37215: 0 to 250, 500, 750 grams  
 37216: 0 to 125, 250, 375 grams

### Physical:

- Dimensions : cm 40 x 16 x 14  
 Packing : cm 55 x 45 x 36  
 Weight : 2.1Kg  
 Shipping Weight : 5.0Kg approx

## Ordering Information

**37215** **ANALGESY-METER**, complete with following standard accessories:-

- 37215-302** Instruction Manual (on USB key)
- 37215-303** Pedal Switch, complete with cable
- 37215-323** Set of discs for additional weight
- 37215-321** Plinth
- 37215-322** Standard Pusher \*
- E-WP008** Mains Cord

\* The pusher can be ordered in special material and/or shapes, according to customer's requirements.

**37216** **ANALGESY-METER**, low-pressure model, suitable for mice, complete with chisel-shaped pusher 37215-326

## Bibliography

### METHOD PAPER

- L.O. Randall and J.J. Selitto: "A Method for Measurement of Analgesic Activity on Inflamed Tissue" *Arch. Int. Pharmacodyn. CXI*, No. 4: 409-419, 1957.

### REFERENCE TO UB ANALGESY-METER (RAT)

- E.K. Joseph et alia: "Vascular Endothelial Cells Mediate Mechanical Stimulation-Induced Enhancement of Endothelin Hyperalgesia via Activation of P2X2/3 Receptors on Nociceptors" *J. Neuroscience* 33 (7): 2849-2859, 2013
- L. Ferrari et alia: "Role of Nociceptor  $\alpha$  CaMKII in Transition from Acute to Chronic Pain (Hyperalgesic Priming) in Male and Female Rats" *J. Neuro-science* 33 (27): 11002-11011, 2013
- D.A. Andersson et alia: "TRPA1 Has a Key Role in the Somatic Pro-Nociceptive Actions of Hydrogen Sulfide" *PLoS ONE* 7(10): e46917, 2012
- K. Király et alia: "The Dipeptidyl Peptidase IV (CD26, EC 3.4.14.5) Inhibitor Vildagliptin is a Potent Antihyperalgesic in Rats by Promoting Endomorphin-2 Generation in the Spinal Cord" *Eur. J. Pharmacol.* 650: 195-199, 2011
- Zs. Helyes et alia: "Involvement of Transient Receptor Potential Vanilloid 1 Receptors in Protease-Activated Receptor-2-induced Joint Inflammation and Nociception" *Eur. J. of Pain* 14 (4): 351-358, 2010

### REFERENCE TO UB ANALGESY-METER (MOUSE)

- K. Sugimoto et alia: "The Impact of Low-Dose Insulin on Peripheral Nerve Insulin Receptor Signaling in Streptozotocin-Induced Diabetic Rats" *PLoS ONE*: 8(8): e74247, 2013
- M.J. Hussey et alia: "Deletion of the Adenosine A2A Receptor in Mice Enhances Spinal Cord Neurochemical Responses to an Inflammatory Nociceptive Stimulus" *Neuroscience Letters* 506(2): 198-202, 2012
- M.S. Nash et alia: "7-tert-Butyl-6-(4-Chloro-Phenyl)-2-Thioxo-2,3-Dihydro-1H-Pyrido[2,3-d]Pyrimidin-4-One, a Classic Polymodal Inhibitor of Transient Receptor Potential Vanilloid Type 1 with a Reduced Liability for Hyperthermia, Is Analgesic and Ameliorates Visceral Hypersensitivity" *J. Pharmacol. Exper. Therap.* 342 (2): 389-398, 2012

# Hot / Cold Plate

Cat. No. 35100

## General

This new instrument, with temperature presettable in the range 2°C to 66°C, can be used as:

- A **conventional HOT PLATE**, to carry out a rapid precise screening of narcotic type analgesic drugs according to the well known Hot Plate Test devised by N.B. Eddy and D. Leinbach.
- As a **COLD PLATE**; the **Cold Plate Test** is useful in studying cold receptors and cold allodynia, a phenomenon very frequently observed in chronic pain on humans.

The **two operating modes** allow for testing at fixed temperature or at increasing/decreasing temperature. The latter experimental scheme is obtained by simply setting on the keypad starting and final temperature.

An optional **auxiliary Hot Plate** (see figure 1, overleaf) permits running two tests simultaneously.

Model **35200 "Thermal Escape"** device is also available, encompassing a Hot/Cold Plate, an auxiliary Hot Plate and an interconnecting corridor, to perform the **Temperature Preference Tests**; see figure 2, overleaf).



**For Rats**

**For Mice**

**IT CAN BE USED AS:**

- Hot Plate
- Cold Plate
- Temperature Preference Tests

## Main Features

- Operating Temperature: 2°C to 66°C presettable by function keys
- Two operating modes: fixed or ramping temperature, for dynamic experiments
- PC Interface: USB
- Computer compatibility: direct connection to a PC, via the dedicated software included as standard
- Data Portability: via the USB Memory-Key, included in the standard package
- Print-out: by optional thermal MiniPrinter

## Instrument Description

The Instrument features:

- a cabinet of original design, incorporating the Heating/Cooling Plate and a Command/Display Module
- a convenient Plexiglas restrainer (suitable to restrain either mice or rats).

The multifunction liquid-crystal graphic display monitors the plate temperature in 0.1°C steps and shows the reaction time in 0.1 s increments on the graphic display; the display also presents all available commands: the operator sets the experiment configuration via the command keyboard located on the right of the display.

The plate temperature can be set by the operator in the range 2-66°C. The extremes of this ample range can be reached, provided the room temperature remains in the interval 18-23°C.

## Data Acquisition

The 35100 is a microprocessor controlled unit. The experimental data, stored in its internal memory can be directly exported to the PC USB or serial ports.

Communication is managed by the dedicated CUB Data Acquisition Software Package, **Cat. 52050-11**, included as standard; this Windows®-based Software enables the user to route to the PC the data originated by UB instruments and store them into individual files, ready to be easily managed by most statistical analysis packages available on the market.

The 35100 is provided with a **memory key**, to record all the experimental data of one or more sessions and to program the experiment layout from a remote PC.

## Options

An **"auxiliary" conventional Hot Plate 35100-002** is available as optional for connection to the 35100, with the same operational features of a complete Hot Plate, powered and controlled by the main unit.

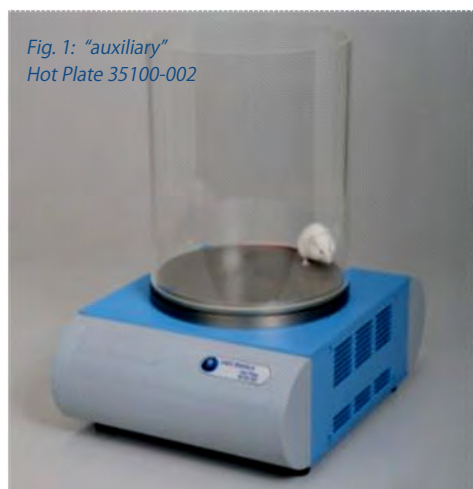


Fig. 1: "auxiliary"  
Hot Plate 35100-002

Model **35200 "Thermal Escape"** device is also available; the 35200 encompasses a Hot/Cold Plate, an auxiliary

Hot Plate and an interconnecting corridor, to perform **Temperature Preference Tests**, see figure below



Fig. 2: "Thermal  
Escape" Device

## Ordering Information

|                  |   |
|------------------|---|
| <b>35100</b>     | <b>HOT / COLD PLATE</b> , standard package, including:-                   |
| <b>35100-001</b> | Cabinet (controller/display and Plate assembly)                           |
| <b>35100-286</b> | Perspex Animal Restrainer, for Mice and Rats), 25cm height                |
| <b>35100-302</b> | Instruction Manual (on USB key)   |
| <b>37215-303</b> | Pedal Switch  |
| <b>E-AU 041</b>  | Memory Key  |
| <b>E-WP 008</b>  | Mains Cord  |
| <b>52050-11</b>  | CUB Data Acquisition Software Package, complete with USB Connection Cable |

### Optional

|                  |   |
|------------------|---|
| <b>35100-002</b> | Auxiliary Hot Plate                       |
| <b>35200</b>     | Thermal Escape Device, complete           |
| <b>57145</b>     | Thermal MiniPrinter with connection cable |

### Physical

|                 |                               |
|-----------------|-------------------------------|
| Universal input | 85-264 VAC, 50-60Hz           |
| Dimensions      | 25x37x47(h)cm with restrainer |
| Weight          | 7.0Kg                         |
| Shipping Weight | 9.6Kg approx.                 |
| Packing         | 65x34x28cm                    |

## Bibliography

- M. Sakurai *et alia*: "Oxaliplatin-induced neuropathy in the rat: involvement of oxalate in cold hyperalgesia but not mechanical allodynia". *Pain* 147 (2009) 165-174
- L. Yu *et alia*: "Effects of calcitonin gene-related peptide-(8-37) on withdrawal responses in rats with inflammation" *EJP* 347 (1998) 275-282
- D. Piomelli *et alia*: "Anandamide suppresses pain initiation through a peripheral endocannabinoid mechanism". *Nature NSC* (2010)



# Plantar Test (Hargreaves Apparatus)

Cat. No. 37370

For Rats

For Mice

AUTOMATIC  
MEASUREMENT OF THE  
ANIMAL RESPONSE

## General

Determination of acute nociceptive thermal threshold in laboratory animals has primarily relied upon the tail flick and hot plate methods.

Although both methods are used frequently in pharmacological studies, they are not without limitation. In addition, neither method has been extended to investigating behavioural responses to hyperalgesia.

The Plantar Test represents a remarkable advance in methodology, as it combines the best features of all other methods of measuring pain sensitivity. Unique to the Plantar Test, **the animal is unrestrained and unhandled during experiments.**



## Main Features

- Automatic detection of paw withdrawal (no visual score needed!)
- I.R. intensity adjustable in the interval 01-99 (in one digit steps)
- Software included
- Modular animal enclosure, from 3 to 12 spaces, conveniently designed to restrain mice or rats
- Optional 37300 Radiometer for calibration
- Data portability via the included memory key
- NEW: orofacial stimulation by optional holders

## Instrument Description

The Instrument basically consists of:-

- a Movable I.R. (infra-red) Source
- a Controller (the picture shows the optional printer 37000-145 mounted on the top panel)



- a Glass Pane supported by columns on a base platform onto which the movable source glides.
- a modular enclosure of new design, in which the 3 spaces can be further divided into 2 or 4 by removable partitions, obtaining up to 12 spaces

After the acclimation period, the I.R. source placed under the glass floor (see the picture) is positioned by the operator directly beneath the hind paw. A trial is started by depressing a key on the I.R. source.

When the animal feels pain and withdraws its paw, the I.R. source switches off and the reaction time counter stops. The withdrawal latency to the nearest 0.1s is automatically determined and recorded.

## Data Acquisition

The 37370 is a microprocessor controlled unit. The experimental data, stored in its internal memory can be directly exported to the PC USB or serial ports.

Communication is managed by the dedicated CUB Data Acquisition Windows®-based Software Package **52050-10**, included as standard, which enables the user to route the experimental data to the PC and store them into individual files, to be managed by most statistical analysis packages available on the market.

The 37370 is provided with a **memory key**, to record all the experimental data of one or more sessions and to program the experiment parameters from a remote PC.

## Calibration Radiometer

Each Plantar Test Unit is accurately calibrated via an **Heat-Flux I.R. Radiometer Cat. 37300**.

The end user should consider this extremely useful optional accessory, which enables the experimenter to:

- Make sure that two or more units deliver thermal nociceptive stimuli (expressed in mW per square cm) of **exactly the same intensity**.
- Measure the I.R. energy (1mW for the duration of 1s corresponds to 1mJ) **in absolute terms**


## Ordering Information

**37370 Plantar Test (Hargreaves' test)**, complete with following standard accessories:

|                  |                                       |
|------------------|---------------------------------------|
| <b>37370-001</b> | Plantar Test Controller               |
| <b>37370-002</b> | Emitter/Detector Vessel, with cable   |
| <b>37000-003</b> | Platform                              |
| <b>37370-327</b> | Supporting columns                    |
| <b>37000-006</b> | Modular Animal Enclosure              |
| <b>37370-005</b> | Framed Glass Pane                     |
| <b>37370-302</b> | Instruction manual (on the USB key)   |
| <b>52050-10</b>  | CUB Software (USB key) with USB cable |
| <b>E-WP 008</b>  | Mains Cord                            |

## Optional Spares & Accessories

|                  |   |
|------------------|---|
| <b>37000-145</b> | Panel-Mount Printer   |
| <b>37300</b>     | Heat-Flux I.R. Radiometer   |
| <b>E-HR 002</b>  | Replacement Bulb  |
| <b>37370-278</b> | Additional stimulation base, complete with glass pane and animal enclosure                            |
| <b>37100</b>     | Set of two <b>Durham Holders</b> for orofacial stimulation ( <a href="#">see separate datasheet</a> ) |

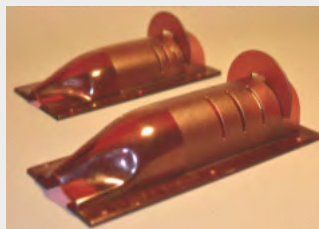

## Physical

|                 |                                  |
|-----------------|----------------------------------|
| Universal Mains | 85-264 VAC - 50-60Hz - 20 W max. |
| Dimensions      | 85 x 40 x 35 cm (assembled)      |
| Weight          | 13.00 Kg                         |
| Packing         | 98 x 49 x 47 cm                  |
| Shipping Weight | 27.50 Kg approx                  |

## Bibliography

### Method Paper:

- K.M.Hargreaves, R.Dubner, F.Brown, C.Flores & J.Joris: "A New and Sensitive Method for Measuring Thermal Nociception in Cutaneous Hyperalgesia" *Pain* 32: 77-88, 1988.
- D.C. Yeomans & H.K. Proudfit: "Characterization of the Foot Withdrawal Response to Noxious Radiant Heat in the Rat" *Pain* 59: 85-97, 1994.

### Papers mentioning UB model:

- D. Piomelli et alia: "Anandamide suppresses pain initiation through a peripheral endocannabinoid mechanism" *Nature NSC*, 2010
- A. Finley et alia: "Sphingosine 1-Phosphate Mediates Hyperalgesia via a Neutrophil-Dependent Mechanism" *PLoS ONE* 8(1): e55255, 2013
- M. Arai et alia: "The miRNA and mRNA Changes in Rat Hippocampi after Chronic Constriction Injury" *Pain Medicine* 14 (5): 720-729, 2013
- P.J. Austin et alia: "G. Chronic Constriction of the Sciatic Nerve and Pain Hypersensitivity Testing in Rats" *JoVE* 61, e3393, doi:10.3791/3393, 2012

# Tail-Flick Unit

Cat. No. 37360

Dedicated Software

Memory Key included

**RAPID and PRECISE  
SCREENING OF  
ANALGESIC DRUGS  
ON THE RAT TAIL**

## General

This new style Tail Flick Unit has been designed to perform rapid precise screening of analgesic drugs via heat stimulation on the rat tail, **according to D'Amour & Smith**, see bibliography. It basically consists of an I.R. source, whose radiant energy of adjustable intensity is focused on the rat tail by an embodied parabolic mirror.

The rat is held by the operator on the instrument unobstructed upper panel (see picture) in such a way that its tail, placed over a flush mounted window, receives the I.R. energy.

The operator starts the stimulus and the related solid state second counter. When the rat feels pain and **flicks** its tail, a sensor detects it, stops the second counter and switches off the bulb. The **reaction time** of the animal is thus determined and automatically recorded.



## Main Features

- Automatic detection of the animal response
- Data portable to USB pen-drive stick or to PC (USB)
- Comfortable, unobstructed working surface (no protruding elements)
- Excellent reproducibility thanks to optics lodged in a rigid structure & electronically controlled I.R. flux

## Instrument Description

The instrument components are neatly arranged in a box of new design, which contains the I.R. source, the sensor, the microcontroller and the electronic circuit.

When the counter stops, the **display** remains frozen on the indicated time. Latency time is thus automatically recorderd.

An inclined **Mouse Restrainer** is supplied as **optional**, to be used with the mouse to compensate for its tendency to hold its tail at 45 degrees up and therefore away from the heat source.

In fact, the availability of **mice** with specific gene(s) knock-outs is driving a substantial shift from rats to mice as a research animal of first choice.



## Data Acquisition

The 37360 is a microprocessor controlled unit. The experimental data, stored in its internal memory can be directly exported to the PC USB or serial ports.

Communication is managed by the dedicated CUB Data Acquisition Software Package, **Cat. 52050-09**, included as standard.

The CUB Windows®-based Software Package enables the user to route the experimental data to the PC and store them into individual files, to be managed by most statistical analysis packages available on the market.

The 37360 is provided with a **memory key**, to record all the experimental data of one or more sessions and to program the experiment layouts from a remote PC.

## Calibration Radiometer

Each Tail Flick Unit is accurately calibrated via an **Heat-Flow I.R. Radiometer Cat. 37300**.

The end user should consider this extremely useful accessory, which enables the experimenter to:

- Make sure that two or more units deliver thermal nociceptive stimuli (expressed in mW per square cm) of **exactly the same intensity**.
- Know the I.R. energy (1mW for the duration of 1s corresponds to 1mJ) in **absolute terms**

## Ordering Information

**37360** **TAIL-FLICK UNIT**, complete with following standard accessories:-

**37215-303** Pedal Switch, complete with cable  
**37360-302** Instruction Manual (on USB key)  
**52050-09** CUB Software (on USB key)  
**52010-323** USB cable  
**E-WP008** Mains Cord

### Accessories and Optionals

**57145** MiniPrinter  
**37300** Heat-Flux I.R. Radiometer  
**E-HR 002** Replacement Bulb  
**37360-325** Mouse Holder, 25mm diam.  
**37360-330** Mouse Holder, 30mm diam.

## Basic Specifications

I.R. Intensity adjustable in the interval  
 01-99 (in one digit steps)  
 Reaction Time three digits, 0.1s steps  
 Calibration via appropriate I.R. Radiometer  
 Universal Mains 85-264 VAC - 50-60Hz - 20 W max.

### Physical

Dimensions 43x22x10cm  
 Weight 4.0 Kg  
 Packing 65x34x28cm  
 Shipping Weight 5.8 Kg approx.

## Bibliography

### Method Paper:

- F.E. D'Amour & D.L. Smith: "A Method for Determining Loss of Pain Sensation" *J. Pharmacol. Exp. Therap.* 72: 74-79, **1941**

### Papers mentioning UB model:

- T.O. Lilius et alia: "The Mineralocorticoid Receptor Antagonist Spironolactone Enhances Morphine Antinociception" *Eur. J. of Pain* online view, **2013**
- J.W. Little et alia: "Spinal Mitochondrial-Derived Peroxynitrite Enhances Neuroimmune Activation During Morphine Hyperalgesia and Antinociceptive Tolerance" *Pain* 154 (7): 978-986, **2013**
- P.J. McLaughlin et alia: "The Mineralocorticoid Receptor Antagonist Spironolactone Enhances Morphine Antinociception" *Eur. J. of Pain* online, **2013**
- T.A. Kosten et alia: "A Morphine Conjugate Vaccine Attenuates the Behavioral Effects of Morphine in Rats" *Progr. in Neuro-Psychopharmacol. and Biol. Psychiatry* 45: 223-229, **2013**
- J. Walsh et alia: "Disruption of Thermal Nociceptive Behaviour in Mice Mutant for the Schizophrenia-Associated Genes NRG1, COMT and DISC1" *Brain Res.* 1348: 114-119, **2012**



## I.R. Heat-Flux Radiometer

Cat. No. 37300

### General

The Heat-Flux I.R. Radiometer Cat. 37300 has been designed to **calibrate** I.R. sources, in particular the classic Tail-Flick 37360 and Plantar Test 37370 of our make.

The purpose of this extremely useful accessory is to make sure different I.R. sources deliver the same **power flux** (expressed in mW per square cm), hence a nociceptive stimulus of the **same intensity**.

The standard package of this portable self-sufficient instrument includes an I.R. Probe, a Digital Meter, and Adaptors for Tail-Flick and Plantar Test (see picture), all parts of neatly lodged in a sturdy plastic case with punched foam lining.



- For Precise Calibration of Infrared Analgesia Meters

- To calibrate the I.R. emission of Ugo Basile Plantar Test & Tail Flick



### Main Features

- Provides a measure of stimulus intensity in mW/cm<sup>2</sup>
- Assures that all infrared instruments are emitting the same level of stimulus intensity

The I.R. output of a I.R. Tail-Flick or Plantar Test may, over the course of one-two years, undergo to 2-3% reduction, due to dust gathered on the optics, darkening of the I.R. bulb, accidental knocks, aging of components due to thermal cycles, etc.


Moreover, if the bulb is replaced or the electronics serviced, output alteration of more significant magnitude, say, 8-10%, may take place.

**The 37300 Radiometer enables the experimenter to:**

- **Check** (and adjust if necessary) **the actual emission of an I.R. source**
- **Ensure** that two or more Tail-Flick/Plantar Test Units deliver thermal nociceptive stimuli of exactly the **same intensity**. Balance them, if necessary.
- **Know the I.R. energy** in absolute terms: 1mW for the duration of 1s corresponds to 1 mJ. A useful datum to compare with any equal or different method/instrument described in the literature.

## Principle of Operation

This simple and reliable I.R. Radiometer uses miniature flat "temperature gradient sensors", whose output signal is proportional to the temperature difference between their top and bottom surface.

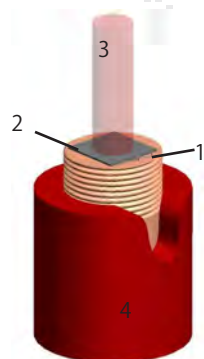

1 Heat-Sink  
2 Temperature Gradient Sensor  
3 I.R. Beam  
4 Plastic Guard

In fact, the temperature reached by the top surface of the sensor attains few degrees Celsius over the heat-sink temperature and hence involves negligible convection and radiation losses.

At the equilibrium, the I.R. power flux  $p$  (mW per square cm) is given by the formula:

$$p = \frac{\Delta T}{\rho d}$$

Where  $\Delta T$  is the temperature difference between top and bottom surfaces of the sensor,  $\rho$  is its thermal resistivity and  $d$  its thickness.

It is notable that the determination of  $p$  is not affected by the heat-sink temperature.  $\Delta T$  only comes into play. The time constant of the system  $\zeta$  (zeta), i.e., the time to reach the equilibrium is given by the formula:

$$\zeta = \rho d C$$

where  $C$  is the thermal capacity \* of the sensor.

$\rho d$  and  $C$  are very small, which leads to the equilibrium and hence to the exact determination of the I.R. power flux in a matter of 3-4 seconds.

Note : \* thermal capacity = mass by specific heat  
\*\* the heat propagates by radiation - conduction - convection

## Practical Clues

The measure, as previously mentioned, requires only a few seconds. The I.R. probe is positioned on the Tail-Flick/Plantar Test, after the suitable adaptor is fitted on the threaded head of its heat sink.

The reading on the digital display gives the I.R. power output in mW per square centimetre.

The calibration (if necessary) of the I.R. radiation source is carried out by adjusting the supply current of the I.R. bulb, see the instruction manuals of the Tail Flick and, respectively, the Plantar Test.

## Ordering Information

|                  |  |
|------------------|--|
| <b>37300</b>     | <b>I.R. HEAT-FLUX RADIOMETER,</b><br>standard package, including:- |
| <b>37300-001</b> | Heat-Flux Meter (complete with<br>cable/connector & 9V battery)    |
| <b>37300-002</b> | Heat-Flux Probe  |
| <b>37300-302</b> | Instruction Manual (on CD)   |
| <b>37300-320</b> | Probe Front Cover  |
| <b>37300-321</b> | Adaptor for Tail-Flick   |
| <b>37300-322</b> | Adaptor for Plantar Test   |
| <b>I-A 073</b>   | Instrument case  |

## Physical

### 37300 complete standard package, lodged in its case:

|                 |            |
|-----------------|------------|
| Dimensions      | 37x32x11cm |
| Weight          | 1.5Kg      |
| Packing         | 46x38x27cm |
| Shipping Weight | 3.6Kg      |

# Dynamic Plantar Aesthesiometer

Cat. No. 37450

- Mechanical Stimulation
- With large platform
- Modular animal cage for Mice & Rats

**ASSESSMENT OF ANIMAL SENSITIVITY TO LIGHT TOUCH OF THE PAW**

## General

The Dynamic Plantar Aesthesiometer has been designed to assess **"touch sensitivity"** on the plantar surface of the rodents.

Somaesthetic (mechanical) stimulation has a long history of effective clinical use to diagnose pathologies of hyper- or hypo-aesthesia, brought about by drugs, neural pathology or experimental lesions, etc., in model and experimental systems using laboratory animals.



## Main Features

- Automatic detection of animal response
- Consistent application of force at an adjustable rate (force ramp)
- Software included as standard
- Data Portability via the Memory-Key provided with the standard package
- Print-out: by optional panel mount or independent thermal mini-printer

The **37450** encompasses:-

- a movable **touch-stimulator unit**, complete with filament actuator and adjustable angle mirror
- a microprocessor controlled **electronic unit**, of new design provided with graphic display, internal memory for data storage, memory stick and optional printer.
- a large **testing surface**
- a modular **animal enclosure**, in which the 3 spaces can be further divided into 2 or 4 by removable partition, thus obtaining up to 12 spaces.

## Operation

The animal moves freely in one of the enclosure compartments, positioned on the testing surface.

After cessation of exploratory behaviour, the operator places the touch-stimulator below the target area of the animal paw, using the adjustable angled mirror to position the filament.

The **START** key provided at both sides of the touch-stimulator handle, invokes the following automatic sequence:

- an electrodynamic actuator of proprietary design lifts a straight metal (NiTi alloy) filament
- the small diameter rod touches the plantar surface and begins to exert an upward force below the threshold of feeling
- the force increases at the preset application rate, until a stop signal is attained, either when the animal removes its paw or when the preset force is reached.

The filament (0.5mm diameter) transmits force over the entire range of typical aesthesiometers. Paw withdrawal reflex is automatically recorded using two metrics: the latency until withdrawal, in seconds, and the force at which paw was withdrawn, in grams.

## Basic Specifications

|                       |  |
|-----------------------|--|
| Starting              | via keys on the touch-stimulator vessel                  |
| Force range           | 0 to 50.0 grams, in 0.5g steps                           |
| Force increasing rate | adjustable in the interval 1 to 20 seconds, in 1 s steps |
| Filament travel       | 12 mm  |
| Latency time          | on graphic display, in 0.1s steps                        |
| Connection to PC      | through DELTA 9-pin connector                            |

## Data Acquisition

The 37450 is a microprocessor controlled unit. The experimental data, stored in its internal memory can be directly exported to the PC USB or serial ports.

Communication is managed by the dedicated CUB Data Acquisition Software Package, **Cat. 52050-12**, included as standard. The CUB Windows®-based Software Package enables the user to route the experimental data to the PC and store them into individual files, to be managed by most statistical analysis packages available on the market.

The 37450 is provided with a **memory key**, to record all the experimental data of one or more sessions and to program the experiment layouts from a remote PC.

## Ordering Information

|                  |  |
|------------------|--|
| <b>37450</b>     | <b>DYNAMIC PLANTAR AESTHESIOMETER</b> , complete with following standard accessories:                              |
| <b>37450-001</b> | Microprocessor controlled electronic unit, with USB key  |
| <b>37400-002</b> | Touch stimulator   |
| <b>37000-003</b> | Large platform   |
| <b>37400-327</b> | Supporting Columns   |
| <b>37450-005</b> | Framed testing surface (perforated platform)   |
| <b>37000-006</b> | Modular animal enclosure (3 to 12 spaces)  |
| <b>37450-302</b> | Instruction manual (on USB key)  |
| <b>37400-321</b> | Set of two 0.5mm diam. NiTi alloy filaments, two calibration weights (5 & 50 g) and accessories, in a plastic case |
| <b>E-WP 008</b>  | Mains Cord   |
| <b>52050-12</b>  | CUB Data Acquisition Software Package, with USB Connection Cable   |

### Optional

|                  |  |
|------------------|--|
| <b>37000-145</b> | Panel-Mount Thermal Printer  |
| <b>57145</b>     | Thermal MiniPrinter  |
| <b>37450-278</b> | Additional stimulation base, with perforated platform and animal enclosure                                 |
| <b>37100</b>     | Set of two <b>Durham Holders</b> for <b>orofacial</b> stimulation ( <a href="#">see separate leaflet</a> ) |

### Physical

|                 |                                  |
|-----------------|----------------------------------|
| Universal Mains | 85-264 VAC - 50-60Hz - 20 W max. |
| Total Weight    | Kg 12.5                          |
| Packing         | 98x49x47cm                       |
| Shipping Weight | Kg 21 approx.                    |

## Bibliography

- R. Lu, A. Schmidt: "Direct Intrathecal Drug Delivery in Mice for Detecting In Vivo Effects of cGMP on Pain Processing" *Methods in Molecular Biology* 1020: 215-221, **2013**
- I.Q. Russe et alia: "Activation of the AMP-Activated Protein Kinase Reduces Inflammatory Nociception" *Journal of Pain* 2, **2013**
- J. Btesh et alia: "Mapping the Binding Site of TRPV1 on AKAP79: Implications for Inflammatory Hyperalgesia" *J. Neuroscience* 33 (21): 9184-9193, **2013**
- V. Brázda et alia: "Dynamic Response to Peripheral Nerve Injury Detected by In Situ Hybridization of IL-6 and its Receptor mRNAs in the Dorsal Root Ganglia is not Strictly Correlated With Signs of Neuropathic Pain" *Molecular Pain* 9(42), **2013**
- D. Piomelli et alia: "Anandamide Suppresses Pain Initiation Through a Peripheral Endocannabinoid Mechanism" *Nature NSC*, **2010**
- P.J. Austin et alia: "G. Chronic Constriction of the Sciatic Nerve and Pain Hypersensitivity Testing in Rats" *JoVE* 61, e3393, doi:10.3791/3393, 2012 <http://www.jove.com/video/3393/chronic-constriction-sciatic-nerve-pain-hypersensitivity-testing>



# PAM

## PRESSURE APPLICATION MEASUREMENT

Cat. No. 38500

### General

The new P.A.M. (Pressure Application Measurement) from Ugo Basile is a novel, easy-to-use tool for measuring mechanical pain threshold in experimental **joint hypersensitivity models in rodents**.

The PAM device has been designed and validated specifically for the mechanical stimulation and assessment of **joint pain**, and therefore is especially useful in studying **arthritis**.

The PAM applies a quantifiable force for **direct stimulation of the joint** and automatic readout of the animal response.

The operator simply wears on his/her thumb a special force sensor, specially designed to apply force to **rat and mouse joints**, and measures the force which elicits the animal response (normally, limb withdrawal).

Each PAM device comes standard with two force sensors, a **large one** useful for stimulating rat joints, a **smaller sensor** recommended to test mice; an optional **paw transducer/applicator** is also available, to stimulate the animal paw.



**Joint Pain**

**Arthritis**

**MECHANICAL PAIN  
THRESHOLD IN:**

- Joint Hypersensitivity
- Chronic Joint Inflammation

### Main Features

- Software included - **NEW 2014 release**
- Maximum Applicable Force: 1500g
- Resolution: 0.1g
- Automatic recording of Limb Withdrawal
- User-controlled application of pressure directly to the joint

## Rationale of the Technique

Arthritis is associated with chronic, debilitating pain in the joints. Current metrics of arthritic pain in animal models are indirect, by scoring the level of motor activity or the animal weight distribution (Barton et al. 2007); while correlating well with the level of joint pain, their metric is a composite picture of complex pain responses, and provides little direct information about local stimulation and locally-evoked responses.

The quantification of localized joint hypersensitivity is not common in animal experiments; in this sense the PAM device represents a step forward toward multifactorial measurement of pain-related behavior in animal research; the **PAM** is the **first instrument designed specifically to apply force to the joint** and automatically detect the animal response.

## Instrument Configuration

**Pressure transducers:** the PAM device comes with 2 transducers, each tested and validated. Both flat and round, the **large transducer** is suitable for rat, the **small one** is ideal for mouse.



Fig. 1: "Joint Transducer"

An optional **paw transducer/appliator** is also available, rapidly transforming the PAM into a Digital Randall-Selitto for pressure application on paws, muscles, tail.



Fig. 2: "Paw Transducer"

**Electronic Unit:** the compact PAM controller connects to the mains or can be battery-operated. A foot pedal switch is provided for manual score of the peak force.



Fig. 3: "PAM device standard package (38500), shown with pedal switch, small and large joint transducer and Usb cable".

## Data Monitoring and Storage

The device includes as standard both a control unit with internal memory and the NEW DCA software for signal monitoring, data transfer and analysis.

Once saved, data can be browsed on the control unit and/or transferred to a PC in proprietary, Excel (.xls) or text (.txt) format, to be managed by most statistical analysis packages available on the market.



## Acknowledgements

The PAM was invented and validated in the University of Edinburgh by the team of Prof. Daniel McQueen, Susan Bond and colleagues and Dr. Harry Brash, who built the first prototypes.

## Ordering Information

|                  |  |
|------------------|--|
| <b>38500</b>     | <b>PAM</b> , standard package, including the following components: |
| <b>38500-001</b> | Electronic Unit  |
| <b>38500-002</b> | <b>Large</b> Joint Transducer                                      |
| <b>38500-003</b> | <b>Small</b> Joint Transducer                                      |
| <b>38500-011</b> | DCA Software (on USB Key)  |
| <b>38500-302</b> | Instruction Manual (on USB Key)                                    |
| <b>38500-303</b> | Pedal Switch   |

All components lodged in a dedicated plastic case

### Options

|                  |   |
|------------------|---|
| <b>38500-006</b> | Paw Transducer                              |
| <b>38550</b>     | PAM, high-pressure model for large animals* |

### Physical

|                 |                              |
|-----------------|------------------------------|
| Weight          | 1.4 Kg (in the plastic case) |
| Shipping weight | 2.7 Kg                       |
| Packing         | 46x38x27cm                   |
| Shipping Weight | 27.50 Kg approx              |

## Bibliography

- **Method Paper:** N. J. Barton et al.: "A novel behavioural technique for measuring hypersensitivity in a rat model of joint pain". *J. Neurosc. Methods*, 163, 67-75, 2007.
- J. Leuchtweis et al.: "Validation of the Digital Pressure Application Measurement (PAM) Device for Detection of Primary Mechanical Hyperalgesia in Rat and Mouse Antigen-Induced Knee Joint Arthritis..." *Methods & Findings in Exp. & Clinical Pharmacol.*, 32(8): 581-589, 2010
- T. Schwagarus et alia: "A New Method for Measuring CFA-induced Mechanical Hyperalgesia in the Rat" *Evotec* 2012
- D. Amorim et alia: "Amitriptyline reverses hyperalgesia and improves associated mood-like disorders in a model of experimental monoarthritis" *Behav. Brain Res* 265: 12-21, 2014
- **38550 (\*)**: P. Di Giminiani et alia: "Nociceptive responses to thermal and mechanical stimulations in awake pigs" *Eur. J. Pain* 17(5): 638-648, 2013

# e-VF

ELECTRONIC VON FREY

Cat. No. 38450

## General

Ugo Basile introduces an electronic apparatus for applying light touch to the rodent foot, the **e-VF, Electronic Von Frey**.

A touch stimulator transducer is mounted on a Perspex bar so that routine procedures may be employed to examine and test the animal skin sensitivity. A **prism** of proprietary design is a useful tool to locate and aim the stimulation area.

The completion of each test may be indicated either by the sudden release of the paw or by pressing the external foot-pedal. The display then gives the operator a summary of the results of the test (i.e. force and time corresponding to the animal response).

The operator may choose to reject the results or to accept them, in which case they are recorded in the e-VF internal memory. The results of several hundred tests may be stored in the e-VF for transfer them to a PC when convenient.


The rate of application of the force is set by the operator and the **NEW** e-VF includes software tools that help in consistently applying the force at the desired rate.



**Sensitivity**

**Allodynia**

**ASSESSMENT OF  
HYPERSENSITIVITY  
IN RATS & MICE**

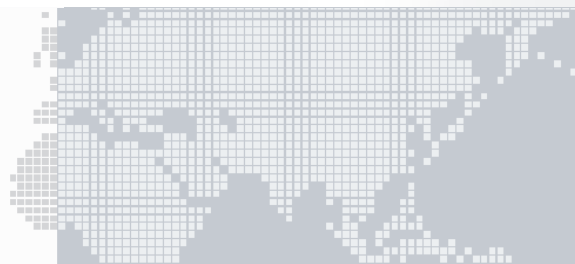

## Main Features

- Software included - **NEW 2014 release**
- Maximum Applicable Force: 1000g
- Resolution: 0.1g
- Automatic recording of animal response
- User-controlled application of force rate
- Location of the target via the original prism-design



## Rationale of the technique

Impaired cutaneous sensation is usually first made evident as a loss of light-touch detection. The Electronic Von Frey was developed to quantify the sensitivity to light touch in the laboratory animal.

The classic instrument for test of touch sensitivity is the **Semmes-Weinstein set of Von Frey Hairs**, i.e., 20 monofilaments in a linear scale of physical force. The Semmes-Weinstein set can be used on rodents, which respond to light touch of the paw, when they feel it, by a paw withdrawal reflex.

However, the involved procedure is tedious and time-consuming because several stimulations must be performed for a single test (i.e., a different filament for each force level).

Compared to the classic Von Frey Hairs, the **Electronic Von Frey (e-VF)** has the advantage of ensuring a continuous force application along the whole force range of the sensor (0-1000g), by using a single rigid metal tip.

Moreover, the metal tip used in the e-VF is the same as the one used in the classic **Ugo Basile Dynamic Plantar Aesthesiometer** (PN 37450), allowing consistent comparison of results among the two instruments.



Fig. 1: "touch stimulator" transducer, including prism  
Grid mesh not included (optional, see ordering information)

## Data Monitoring and Storage

The device comes standard with both a control unit with internal memory and the **new DCA software** for signal monitoring, data transfer and analysis.

Once saved, data can be browsed on the control unit and/or transferred to a PC in proprietary, Excel (.xls) or text (.txt) format, to be managed by most statistical analysis packages available on the market.

## Ease of use

The e-VF device has been designed to make sensitivity experiments easy and consistent, thanks to its:

- Effective **peak detector**, for a reliable and automated detection of the animal response
- **Ratometer** and **Slope** feature, ensuring the desired force is applied at a consistent rate



- **NEW Software**, acting as a quality control tool, by showing the applied pulling force (red line), the desired target force rate (blue line), and the peak detection in real time, see picture above

## Instrument configuration

The e-VF comes as a complete package including **touch stimulator transducer** with **prism**, **electronic unit** with power supply, foot pedal, **software** & **USB cable**. The mesh grid with platform, and animal enclosure are optional.



Fig. 2: electronic unit, usb cable and foot pedal

## Ordering Information

**38450 e-VF, Electronic Von Frey**, complete with following standard parts

**38450-001** Electronic Unit, with power supply

**38450-004** Touch-Stimulator Transducer with

**38450-310** Prism

**38500-011** DCA Software (on USB Key)

**38450-302** Instruction Manual (on USB key)

*All components lodged in a dedicated plastic case*

### Options

**37450-005** Perforated Metal Sheet for plantar stimulation

**37450-278** Base assembly for plantar stimulation, with perforated metal sheet & animal enclosure

### Physical

|                 |            |
|-----------------|------------|
| Weight          | 1.4Kg      |
| Shipping Weight | 2.7Kg      |
| Packing         | 46x38x27cm |



# Von Frey Hairs

Cat. No. 37450-275

## General

Von Frey hairs (named after the German physiologist Max von Frey, 1852–1932) were been originally produced from animal and human hairs of different diameter; nowadays they are nylon monofilaments; the diameter determines the resistance of the monofilament to bending. A filament is placed perpendicularly to the skin with slowly increasing force until it bends, thereby determining the amount of force applied.

The **Aesthesio®** set of 20 monofilaments is based on the Semmes Weinstein monofilament set, **but now features retractable filaments** to protect the filament and allow the evaluator to carry a few around in a pocket.

The set of monofilaments provides an approximately logarithmic scale of actual force, and a linear scale of perceived intensity.

They have a long history of effective use in clinical settings, and can be used to diagnose pathologies of hyper- or hypo-aesthesia.

Subsets within the set of 20 probes distinguish pathologies on different parts of the body (foot, hand, lip, cheek, etc.).

Individual filaments are also sold separately individually.



Hypersensitivity

Touch Threshold

Semmes Weinstein  
Von Frey Filaments  
for Touch  
Assessment

## Main Features

- 20 Filament Kit
- Graded Series of Nylon Monofilament, color-coded
- Rotating sleeve protects precision filament while in closed position

Von Frey Filaments have a long history of effective use in clinical settings, and can be used to diagnose pathologies of hyper- or hypo-aesthesia.

The operating principle remains the same: when the tip of a fiber of given length and diameter is pressed against the skin at right angles, the force of application increases as long as the researcher continues to advance the probe, until the fiber bends. After the fiber bends, continued advance creates more bend, but not more force of application.


**This principle makes it possible for the researcher using a hand held probe to apply a reproducible force, within a wide tolerance, to the skin surface.**

Rodents exhibit a paw withdrawal reflex when the paw is unexpectedly touched. The Touch Test™ Sensory Evaluator can be used on the Plantar surfaces of the foot of a rat or mouse, and the animal will indicate sensation by pulling back its paw.

**Replacement filaments available.** Subsets within the set of 20 probes distinguish pathologies on different parts of the body (foot, hand, lip, cheek, etc.). **Rotating sleeve** protects precision filament while in closed position.



Grip sleeve and rotate handle

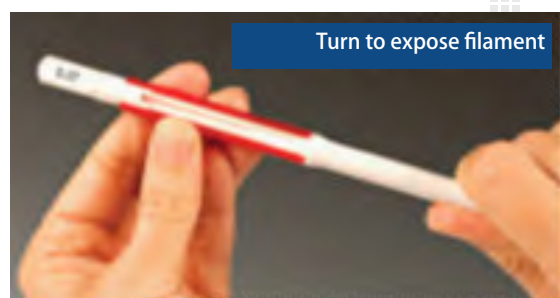

Turn to expose filament



Pivot filament head to begin testing

## Accessories

For easy and quick stimulation of the plantar surface with Von Frey filaments, we offer a 90x38cm **perforated metal platform**, cat. 37450-005. Laser-cut perforations form a mesh-like open grid of square holes ~5X5 mm; intervening metal grid is ~1mm wide, comfortable to the animal and easy to view the target area of the paw.

The shelf is coated with a polymer resin that is easy to clean and which will not be spoiled by fluids or waste materials. Mount the shelf on the wall.

In alternative we offer a **shelf with 40 or 80cm legs**, 37450-045 & 37450-085 respectively, which can be completed with our standard animal enclosure 37000-006; the latter is the **modular enclosure**, used with our Plantar Test & Dynamic Plantar Aesthesiometer, in which the 3 spaces can be further divided by partitions into 2 or 4, thus lodging up to 12 rats or mice.



You might also consider the **complete stimulation base** 37450-278, including supporting columns, shelf, and animal enclosure.

## Ordering Information

**37450-275** **Aesthesio®** Sensory Evaluator, Kit of 20 Von Frey filaments in a carrying case

### Physical

Weight 0.4 Kg  
Shipping Weight 0.9 Kg  
Packing 24x22x5cm

### Options

**37450-005** Large Perforated Metal Platform (testing shelf) for plantar stimulation  
**37450-045** Platform 37450-005, with 40cm legs  
**37450-085** Platform 37450-005, with 80cm legs  
**37000-006** Multiple-configuration animal-enclosure, from 3 to 12 spaces  
**37450-277** Set of 20 VonFrey Filaments 37450-275 & complete base assembly 37450-278  
**37450-278** Base Assembly for plantar stimulation, incl. supporting columns, perforated metal sheet and multiple-configuration animal-enclosure, from 3 to 12 spaces

## Orofacial Stimulation Test

*Fehrenbacher, Henry and Hargreaves Method*

Cat. No. 31300

Mechanical Nociception

Thermal Nociception

Trigeminal  
hyperalgesia

### General

The **Orofacial Stimulation Test** by **Ugo Basile** measures hypersensitivity to thermal or mechanical stimulation of the trigeminal area.

Rats voluntarily contact a thermal or a mechanical stimulator with their *unshaved vibrissal pad* in order to access a food reward. Metrics obtained are the **duration** of feeding and the **number of feeding** attempts, measured by interruption of an infrared barrier traversing the opening to the reward.

Feeding duration and number of attempts are strongly dependent on changes in the applied thermal or mechanical stimulus.



### Main Features

- Mechanical and thermal nociception assays within the same experiment
- High throughput: up to 16 animals can be tested simultaneously
- Intact vibrissal pad, as the test does not require any shaving
- The ORO-Software, included as standard, manages up to 16 cages



## Instrumentation and Methodology

Orofacial pain problems are common and involve structures and mechanisms unique to the trigeminal nerve. Few methods are currently available for orofacial pre-clinical research, and none incorporates parallel measurement of mechanical or thermal stimulation within the same experiment.

Moreover, while most of the current assays measure unlearned behaviors, such as flinching or withdrawal reflexes, the new **Orofacial Stimulation Test**, developed by Fehrenbacher, Henry and Hargreaves, integrates higher-order brain functions into measurements of orofacial nociception.


This innovative approach permits highly integrated nociceptive responses to thermal or mechanical stimulation.

Animals are trained & tested in standard home cages.

The snout is inserted through an opening to lick the reward bottle. Tests are performed in the presence of thermal or mechanical stimuli contacting the vibrissal pad.

Following treatment to induce hypersensitivity, (e.g., trigeminal ligation or injection) trials are repeated to determine the effect of treatment on feeding behavior/reward. Assay sensitivity (inflammation-induced decreases in feeding behavior and reversal of hypersensitivity by local and systemic administration of analgesics) has been proven (Hargreaves et alia, ms in prep.); the feeding behavior is strongly correlated to mechanical or thermal orofacial nociception, as the animal must contact the stimulator in order to access the food reward.

The **Ugo Basile Orofacial Stimulation Test** quantifies feeding behavior by measuring and recording the beam-break number and duration (including min, max and mean), via the **ORO-Software** included; the software acquires data from up to 16 cages simultaneously.

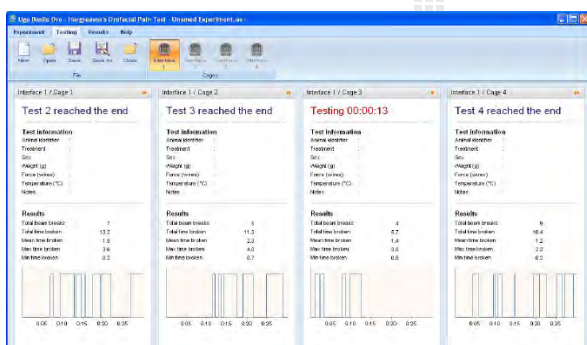


Orofacial Software: testing window


The Data are shown in real-time both as numeric summary results and in a graphic format. Data are automatically analyzed across time according to an adjustable time window, independently viewable for each of the 16 cages. The results of all the tests are available in a spreadsheet format which can easily be copied to other programs for further analysis.

Either the thermal or the mechanical stimulator is mounted onto a **stimulus/detection "wall"**, which

also incorporates a drinking bottle and fits inside standard rat home cages (e.g. Tecniplast or Allentown).

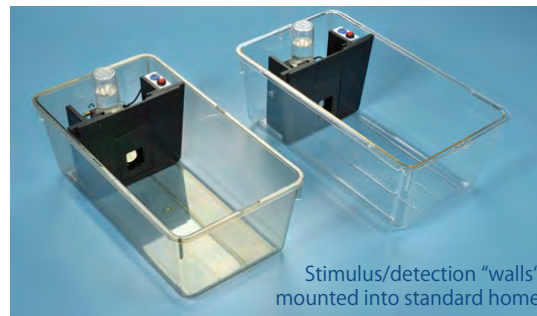

The **thermal stimulator** relies on a copper tubing loop and a circulating water bath, whose temperature can be adjusted from ambient to 70°C, to reach hot nociceptive thresholds. Chin inserts are included to test animals of different size.

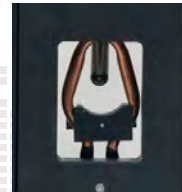

The **mechanical stimulator** relies on thin wires attached to a mounting plate. The system comes with several plates, each with a different number of wires in order to apply different force levels to the animal vibrissal pad.



A kit of Mouse adaptors for both thermal and mechanical stimulation is available, see ordering information.

*The "System and Method for Assessing Hypersensitivity to Orofacial, Thermal and Mechanical Stimulation" (U.S. Provisional Patent Application 61/235,590) was invented by J. Fehrenbacher, M. Henry and K. Hargreaves, in the Laboratory of Dr. Hargreaves at UT San Antonio and developed commercially by Ugo Basile R&D. Dr. Fehrenbacher is now at IUPUI, Indianapolis.*

## Ordering Information

- |                  |   |
|------------------|---|
| <b>31300</b>     | Complete system for one animal  |
| <b>31320</b>     | Complete system for two animals   |
| <b>31340</b>     | Complete system for four animals  |
| <b>31300-001</b> | Electronic unit (four channels)   |
| <b>31300-002</b> | Additional cage assembly (includes thermal and mechanical stimulators and feeding detector) |
| <b>31300-003</b> | Circulating water bath  |
| <b>31300-010</b> | ORO-Software, for data acquisition and analysis from up to 16 cages                         |
| <b>31300-323</b> | Optional Kit of Mouse adaptors for thermal and mechanical stimulation (for 1 cage)          |

## Bibliography

- M. Prochazkova et alia: "Activation of cyclin-dependent kinase 5 mediates orofacial mechanical hyperalgesia" *Molecular Pain* 9:66: 1-12, 2013
- M. Cha et alia: "Assessment of chronic trigeminal neuropathic pain by the orofacial operant test in rats" *Behav. Brain Research* 234: 82-90, 2012
- Fehrenbacher, J.C. et al. 2010. "Characterization of a novel orofacial behavioral assay to assess hyperalgesia to thermal and mechanical stimulation". (submitted).



## Durham Animal Holders

*New animal holders for trigeminal stimulation*

Cat. No. 37100

- Orofacial Pain assessment
- Mechanical and Thermal Nociception

Trigeminal  
hyperalgesia

### General

The **Durham Animal Holders** are the newest accessory for use with the **Plantar Test / Hargreaves Test**, and **Dynamic Plantar Aesthesiometer**, manufactured by Ugo Basile.

These animal holders complete the scope of the infrared (IR) thermal stimulus of the Plantar Test, or the mechanical stimulus of the Dynamic Plantar Aesthesiometer, used for assessing hind paw withdrawal. This new invention allows the application of the same stimulus to the region innervated by the trigeminal nerve.

The 37100 includes two holders, form molded for testing specific size ranges of animals; the two sizes have been optimized for young adult rats as well as for bigger rats.



*"Very nicely done - easy to use and it greatly facilitates consistent handling of animals"*

*Dr. Ken Hargreaves, UT Texas*

### Main Features

- Correlation thresholds in submandibular (trigeminal) region and hindpaw plantar surface
- Test orofacial nociception using a standard Plantar Test (Hargreaves) device, a Dynamic Plantar Aesthesiometer, or an eVF Electronic Von Frey

## Innovative design and material

The Durham Holders are designed to hold an animal comfortably and effectively. They are made of a proprietary polymer with a deep-red color which appears dark to the animal.

The holders conformation is optimized to two specific animal size ranges; the smaller holder will accommodate rats from 175 grams to 250 grams, and the larger holder will accommodate animals from 225 grams to over 400 grams.

In practice, the rat crawls in happily and becomes snugly nestled within the holder. Normally the rats don't back out, but inserting the vertical back plate ensures that the animal stays in place.

The position of the removable back panel insert can be adjusted from slot to slot, which allows the animal to be securely held in place, without being crowded.

The rat crawling towards the front helps quite a lot and the subject is almost self-positioning for applying the IR stimulus to the submandibular region of the rat face.

## Access Panels

There are two different windows through which the stimulus may be presented:

- **Submandibular access panel:**

The opening under the chin is a perfectly sized rectangular aperture just below the animal's chin. It allows the IR or mechanical stimulus to be aimed precisely and to stimulate the area innervated by the mandibular branch of the trigeminal nerve.

The aperture is large enough that both right and/or left side may be individually stimulated!

- **Plantar access panel:**

The holder allows the animal to be positioned in such a way to use the classic Plantar Test instruments for stimulating the hindpaw, as well as the areas innervated by the trigeminal nerve.



The picture above shows a Durham Holder positioned on a classic Ugo Basile Plantar Test (Hargreaves) device.

## Rationale of the technique

The Durham Holders have distinct advantages which make them ideal as accessories to the classical Hargreaves test and they represent a step forward toward a multifactorial measurement of pain-related sensitivity in animal research.

They may appear similar to the classic Broome style animal holder; however, those animal restrainers are clunky, made of clear acrylic, and do not have stimulus apertures, so they could never be used for this stimulation.

Quantification of localized hypersensitivity is common in the clinic, but not in animal experiments.



## Acknowledgements

The Durham Holders were invented and validated at the Center of Biomedical and Life Sciences at Missouri State University; specifically, in the laboratory of Dr. Paul Durham, director of Biomedical & Life Sciences and Professor of Cell Biology at Missouri State University.

Filip Garrett and Allison Overmyer performed the validations. Prototypes were put together by Larry Vause.

## Ordering Information

**37100 Set of two Durham Holders for rats:**  
37102 medium size  
37103 large size

## Physical

Weight 0.4 Kg (two holders)  
Gross weight 1.0 Kg  
Packing 39x27x21cm

## Bibliography

- F.G. Garrett, et alia: "Validation of a Novel Rat-Holding Device for studying heat- and mechanical-Evoked Trigeminal Nocifensive Behavioral Responses" *J. Orofacial Pain*, 26 No. 4, 336- 344, 2012
- F.G. Garrett, A.E. Overmyer, L.A. Vause, J.L. Hawkins, J.B. Hayden, and P.L. Durham "Development of a novel device for measuring withdrawal latency by thermal stimulation in rodent facial pain models using the Hargreaves Plantar Apparatus" Poster presented at SFN 2010

## Grip Strength Meter

Cat. No. 47200


### General

The Ugo Basile Grip Strength Meter automatically measures grip-strength (*i.e.* peak force and time resistance) of forelimb or hindlimb (via the optional grid) in rats and mice.

The Grip Strength test is a perfect complement to the gold standard Ugo Basile Rota-Rod device for motor coordination and motor function experiments. The effects of drugs, toxins, muscle relaxants, disease, ageing or neural damage on muscle strength may be assessed.

The animal is placed over a base plate, in front of a grasping tool (either T-shaped, trapeze-shaped or grid), whose height is adjustable.

The bar is fitted to a force sensor connected to the control unit, which can be used as a stand-alone or connected to a PC via the USB port, for monitoring and data recording, via the **NEW** software provided as standard

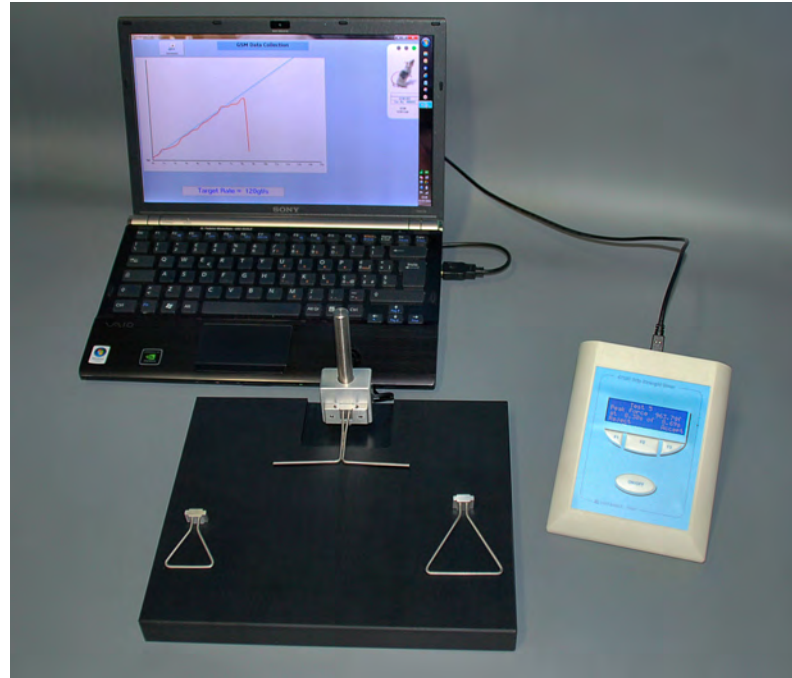

High Consistency  
with force-rate  
monitoring tool

for Rats

for Mice

### Features and Benefits

- Software included - **NEW 2014 Release**
- Grasping tools included for rats and mice (grid is optional)
- No calibration needed
- Force-rate monitoring (via software or LCD display)
- Grasping bar / grasping trapeze positioned at adjustable height
- Maximum applicable force 1500g; resolution 0.1g



## Rationale of the Grip Strength test

When pulled by the tail, the animal grasps at the bar. Rodents instinctively grab anything they can, to try to stop this involuntary backward movement, until the pulling force overcomes their grip strength. After the animal loses its grip on the grasping bar, the peak amplifier **automatically stores the peak pull-force achieved by the limbs** and shows it on the display.

The instrument basically consists of a base plate of black sand-blasted Perspex, complete with a force transducer and a grasping device (bar, trapeze or the optional grid), which can be positioned at an adjustable height.

The force transducer has a maximum applicable force of 1500g, with a resolution 0.1g.

The transducer incorporates a proprietary memory chip to store all calibration parameters, so that no further calibration is required for normal use; moreover, the controller will prompt to auto-zeroing routine at every measurement to automatically adjust any offset.

## Data Monitoring and Storage

The device comes standard with both a control unit with internal memory and the **new DCA software** for signal monitoring, data transfer and analysis.

Once saved, data can be browsed on the control unit and/or transferred to a PC in proprietary, Excel (.xls) or text (.txt) format, to be managed by most statistical analysis packages available on the market.


## Ease of use

The GSM device has been designed to make sensitivity experiments easy and consistent, thanks to its:

- Effective **peak detector**, for a reliable and automated detection of the animal response
- **Ratometer** and **Slope** features, ensuring the desired force is applied at a consistent rate
- **NEW Software**, acting as a quality control tool, by showing the applied pulling force (red line), the desired target force rate (blue line), and the peak detection in real time.

The experimenter can consistently apply the force (i.e. pull the animal) at the desired rate, by simply making sure that the red trace lays on the blue line, see figure 1

## Optional Grids

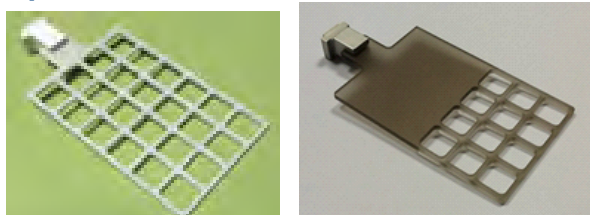

Optional grids are available for integrated measurement of the four limbs (left) or hindlimbs (right)



Figure 1: Screenshot of the GSM software showing the force trace (in red) and the desired target force rate (in blue) - slope function

## Ordering Information

|                  |  |
|------------------|--|
| <b>47200</b>     | <b>Grip-Strength Meter</b> , new model for rats & mice, complete with following standard accessories |
| <b>47200-001</b> | Control Unit, with Power Supply  |
| <b>47200-002</b> | Force Sensor   |
| <b>47200-004</b> | Baseplate and upright  |
| <b>38500-011</b> | DCA Software (on USB Key)  |
| <b>M-LM 589</b>  | T-shaped Grip-Bar  |
| <b>M-LM 590</b>  | Grip-Trapeze for Rat   |
| <b>M-LM 588</b>  | Grip-Trapeze for Mouse   |
| <b>38500-303</b> | Pedal Switch   |
| <b>52010-325</b> | USB Cable  |

All components lodged in a dedicated plastic case

## Optional

|                  |                                   |
|------------------|-----------------------------------|
| <b>47200-325</b> | Mouse Grasping Grid               |
| <b>47200-326</b> | Mouse Grasping Grid ("blind" top) |

## Physical

|                 |            |
|-----------------|------------|
| Weight          | 4.8kg      |
| Shipping weight | 6.5Kg      |
| Packing         | 46x38x27cm |

## Bibliography

- G.J. Huang et alia: "Ectopic Cerebellar Cell Migration Causes Maldevelopment of Purkinje Cells and Abnormal Motor Behaviour in Cxcr4 Null Mice". *PLoS ONE* 9 (2): e86471, 2014 (Mouse)
- R. Barone et alia: "Endurance Exercise and Conjugated Linoleic Acid (CLA) Supplementation Up-Regulate CYP17A1 and Stimulate Testosterone Biosynthesis". *PLoS ONE* 8 (11): e79686, 2013 (Mouse)
- N. Lange et alia: "Behavioural and Pharmacological Examinations in a Transgenic Mouse Model of 2 early-onset torsion dystonia". *Pharmacology, Biochemistry and Behavior* 97 (4): 647-655, 2011 (Mouse)
- M. Savic et alia: "Behavioural Characterization of Four Endemic Stachys Taxa". *Phytother. Res.*, 2010 (Rat)
- A. Coluccia et alia: "Developmental Omega-3 Supplementation Improves Motor Skills in Juvenile-Adult Rats". *Int. J. Devl Neuroscience* 27: 599-605, 2009 (Rat)
- A. Coluccia et alia: "Late Embryonic Exposure to All-Trans Retinoic Acid Induces a Pattern of Motor". *NeuroToxicology* 30: 1120-1126, 2009 (Rat)



## Multiple Activity Cage

Cat. No. 47420

### General

An animal level of general activity or locomotion is an indicator of drug action, toxic substances, neurological damage, or daily rhythms in activity.


Activity data may be automatically and unobtrusively collected by many methods.

The Ugo Basile Activity Cage has proved to be of great value to record spontaneous co-ordinate activity in rats and mice (individual or groups) and variation of this activity in time.

As the animal moves about a clear acrylic cage, it interrupts one or more infrared beams. The beams are arranged in an array of emitters on one side of the cage, detectors on another.

The lower IR array monitors horizontal movement while the upper IR array monitors vertical or rearing activity.

The number of beam breaks is correlated with the amount of movement about the cage.

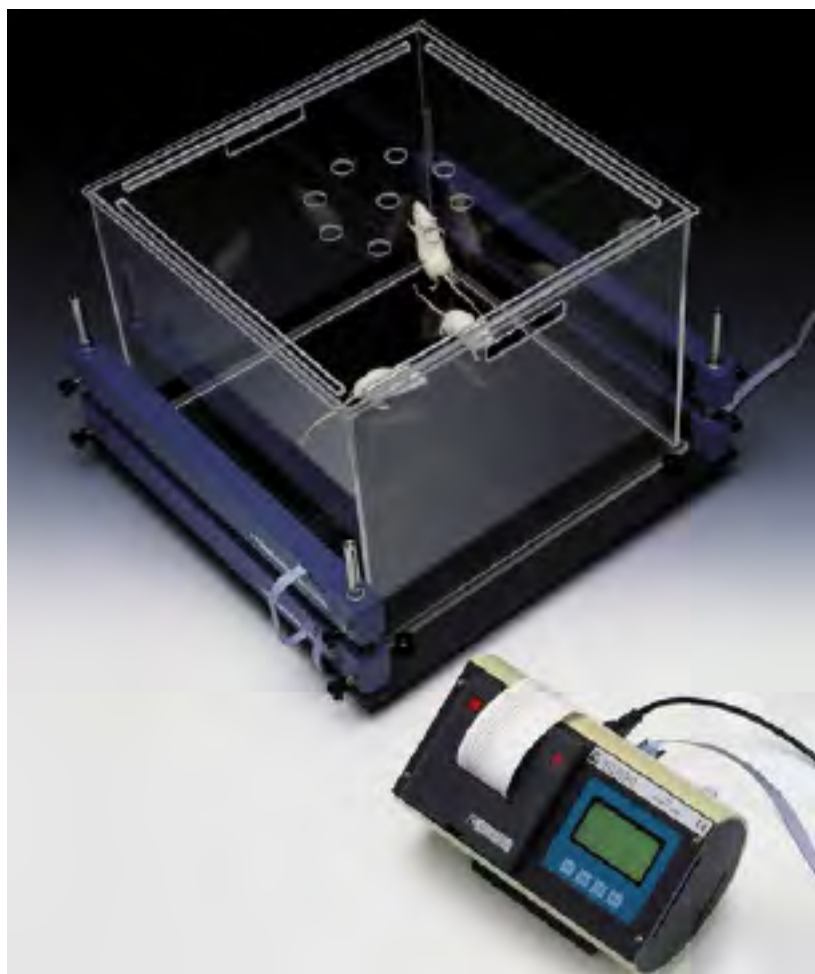

**With dedicated software included**

**STAND-ALONE CONTROLLER WITH EMBEDDED PRINTER FOR GLP AND DATA SAFETY**

### MAIN FEATURES

- Measures **horizontal and vertical activity in rats and mice**, useful in the following types of investigation:-
- **General Toxicology**, ascertaining the action of a drug on the animal's activity
- **Psychopharmacology**, screening drugs which are potentially active on the CNS
- **Behavioural Sciences**, in evaluating the variations of spontaneous activity after changes in environmental conditions

## Instrument Description

The **47420 MULTIPLE ACTIVITY CAGE** package comprises:

- an **Electronic Unit**, Cat. 7441
- an **I.R. Beam Cage**

This set-up can accept up to 5 additional cages, for a total of 6.

### Electronic Unit

The **7441**, designed to process the data originated by **up to 6 Cages**, incorporates a graphic display, a thermal printer and a serial port RS232 for direct connection to the PC via the software Cat. 52050 included. A serial to USB adaptor is also included.

The graphic display presents all available commands. The operator sets the experiment configuration via the keyboard located below the display.

The activity data are displayed at preset intervals and printed/routed to the computer according to the selected configuration. The data can be customized by adding animal & experiment numbers, gender, etc.

Its internal memory is capable to store the data of several experiments, to be unloaded to the PC later.

### Cage

The **7433** Cage consists of a cubicle, entirely made of clear Perspex, dimensioned 41x41x33(h)cm. Upper lid and bottom catch pan detachable for cleaning.

The cubicle rests on a sturdy base, provided with four vertical notched bars of stainless steel to which the horizontal/vertical detecting systems 7435 and/or 7436 can be fastened.

The **7435** consists of two facing blocks containing an I.R. array of emitters and, respectively, sensors, to record the **horizontal activity**. A similar system, Cat. **7436**, whose height can be adjusted, assesses the **vertical activity (rearing)**.

**Open-field cages** are also available, in different dimensions and colors: ask for additional details

## Data Acquisition

The electronic unit is microprocessor controlled and features direct PC output. Internally-stored data can be routed via a 9-pin D-type connector to the PC serial port (RS232).

Data output is managed by **52050-04** Data Acquisition Software Package (Windows® based), which enables the research worker to store the data into individual files, ready to be easily managed by most statistical analysis packages available on the market.

Combination with ANY-maze videotracking software is also possible, to integrate the quantitative measure of general locomotor activity, collected by our Activity cage, with more detailed information about the animal activity.

Moreover, the 47420 will add vertical activity (rearing) to videotracking data. **Ask for additional information!**

### Ordering Information

**47420 MULTIPLE ACTIVITY CAGE**, standard package, including following parts:

|             |   |
|-------------|---|
| <b>7441</b> | Electronic Unit   |
| <b>7433</b> | Animal Cage   |
| <b>7435</b> | Set of emitter/receiver sensor arrays for horizontal activity |
| <b>7436</b> | Set of emitter/receiver sensor arrays for vertical activity   |

**47420-302** Instruction manual (on CD)

**37400-305** Package of 10 Heat Sensitive Paper Rolls

**E-WP008** Mains Cord

**52050-04** Dedicated Software Package CUB

**52010-320** USB to serial port converter

**52010-322** Serial cable 9 to 9 pin

### Physical

|        |             |                              |
|--------|-------------|------------------------------|
| Weight | <b>7441</b> | 2.7Kg                        |
|        | <b>7433</b> | 11.8Kg (including 7435/7436) |

### Dimensions

|             |            |
|-------------|------------|
| <b>7441</b> | 27x16x19cm |
| <b>7433</b> | 54x50x37cm |

Shipping weight 26Kg (whole set-up)

Packing 80x60x44cm

### Bibliography

- V. Labrie et alia: "**Genetic loss of D-amino acid oxidase activity reverses schizophrenia-like phenotypes in mice**" *Genes, Brain and Behavior*, 9: 11–25, 2010
- J. Vlainic, et alia: "**Zolpidem is a potent anticonvulsant in adult and aged mice**" *Brain Res.*, 1310 181–188, 2010
- A. Betourne et alia: "**Central locomotor and cognitive effects of a NPFF receptor agonist in mouse**" *Peptides* 31, 221–226, 2010
- A. Marazioti et alia: "**Somatostatin Receptors in the Ventral Pallidum/Substantia Innominata Modulate Rat Locomotor Activity**" *Psychopharmacol.*, 181:2, 319–326, 2005
- W. Ponti et alia: "**In vivo Model for the Evaluation of Molecules Active Towards Transmissible Spongiform Encephalopathies**" *Veter. Res. Communicat.*, 28:1, 307–310, 2004
- T. Dolezal et alia: "**Guaifenesin Enhances The Analgesic Potency of Paracetamol in Mice**" *Arch. Pharmacol.*, 366:6, 551–554, 2002

# Rat Rota-Rod

Cat. No. 47700

## General

The "Rota-Rod" technique has been originated by a 1957 paper of N.W Dunham and T.S Miya and has proved to be of great value in research involving screening of drugs which are potentially active on motory coordination.

The **Ugo Basile Rota-Rods**, the original Rota-Rods, are the result of many years of research in cooperation with the latest development in behavioral and pharmacological research.

The "Rota-Rod" technique is essential in screening drugs which have side effects on motor coordination.

The new Rota-Rod for Rats has four lanes, and operates in constant speed, accelerating, and reverse/rocking modes! The tiltable digital panel shows accrued rotations and time on the rod.



- THE ORIGINAL ROTA-ROD:  
THOUSANDS OF CITATIONS SINCE 1960s!
- THREE OPERATION MODES:  
CONSTANT SPEED, ACCELERATING OR ROCKING

## Main Features

- Adjustable speed (2-80 rpm) and acceleration ramp (6" - 600")
- Tiltable graphic display for optimal reading
- PC Interface: serial and USB (via the adaptor provided)
- Computer compatibility: direct connection to PC (via the **dedicated software included** as standard)

## Basic Features

The new Ugo Basile Rota-Rod replaces both previously available constant speed and accelerating models. They basically consist of four 6cm diam. drums, suitably machined to provide grip.

Five flanges divide the four 8.7cm lanes, enabling four rats to be on the treadmill simultaneously.

When a rat falls off its cylinder section on to the plate below, the plate trips thereby recording the animal's endurance time in seconds. Height to fall is 30cm.

A large, very readable backlit graphic display shows the actual angular speed (RPM). At the end of a run, the display shows for each animal the running time and the instrument rotation speed at the time that animal fell off.

The panel can be oriented to select the most comfortable angle for the operator, to avoid glare, etc.

The main features, digitally preset by function keys are the following:-

- the angular speed can be preset in the range 2-80 RPM (revolutions per minute)
- in the acceleration mode, the change of speed can be preset in 6 second-10 minute interval in 6 second steps.
- reverse rotation can be selected, which takes place at minimum speed, at the end of a programmable acceleration-deceleration sequence.
- a rocking motion is also presettable, with adjustable angular amplitude, speed and acceleration.

## Data Acquisition

The 47700 is microprocessor controlled and features direct PC output. Internally stored data can be routed to the PC serial port (RS232) or USB (via adaptor provided).

Data output is managed by 52050-07 Data Acquisition Software Package (Windows® based), which enables the research worker to store the data into individual files, ready to be easily managed by most statistical analysis packages available on the market. **Ask for details!**

### 47800 COMBO-PACKAGE FOR MOUSE & RAT

If you plan to work with both rats and mice, you should consider the Combination Package 47600 (Mouse) + 47700 (Rat) Rota-Rod, offered at a special price.



47600 Mouse Rota-Rod

## Ordering Information

**47700** RAT ROTA-ROD, standard package, including:

- 47700-320** Trip Plate, complete (4 pieces)
- 47700-321** Transmission Belt (2 pieces)
- 47700-302** Instruction Manual (on USB drive)
- 52050-07** Dedicated Software Package (on USB drive)
- 52050-322** Serial Cable
- 52010-320** Serial to USB adaptor
- E-WP008** Mains Cord

## Optional

- 57145** Thermal Mini-Printer
- 47800** Combination Package 47700 Rat Rota-Rod and 47600 Mouse Rota-Rod

## Physical

|                 |                      |
|-----------------|----------------------|
| Universal input | 85-264 VAC, 50/60 Hz |
| Dimensions      | 50(w)x49(d)x63(h)cm  |
| Weight          | Kg 10.0              |
| Shipping Weight | Kg 18.0 (approx.)    |
| Packing         | 80x60x44cm           |

## Bibliography

### Method Papers

- N.W. Dunham & T.S. Miya: "A Note on a Simple Apparatus for Detecting Neurological Deficit in Rats & Mice" *J. Am. Pharmaceut. Assoc., Scientific Edit., XLVI: No. 3, 1957*
- B.J. Jones & D.J. Roberts: "The Quantitative Measurement of Motor Incoordination in Naive Mice Using an Accelerating Rotarod" *J. Pharm. Pharmacol.* 20: 302-304, 1968

### Papers Dealing With Rota-Rod Technique

- L.D.C. Mannelli et alia: "Involvement of  $\alpha 7$  nAChR Subtype in Rat Oxali-platin-induced Neuropathy: Effects of selective activation" *Neuropharmacology* 79: 37-48, 2014
- F. Barthel alia: "Long-term Application of Glycine Transporter Inhibitors Acts Antineuropathic and Modulates Spinal N-methyl-D-aspartate Receptor Subunit NR-1 Expression in Rats" *Anesthesiology*, 2014
- A. De Visser et alia: "The Adjuvant Effect of Hypertension Upon Diabetic Peripheral Neuropathy in Experimental Type 2 Diabetes" *Neurobiol. Of Disease* 62: 18-30, 2014
- J. Wang et alia: "Comparison of Anticoagulation and Thrombolysis Treatments in a Rat Model of Superior Sagittal Sinus Thrombosis" *Intl. J. Neuroscience* posted online 2013
- F. E. Padovan-Neto et alia: "Anti-Dyskinetic Effect of the Neuronal Nitric Oxide Synthase Inhibitor is Linked to Decrease of Fosb/Deltafosb Expression" 541: 126-131, 2013
- C.D. Heldermon et alia: "Development of Sensory, Motor and Behavioral Deficits in the Murine Model of Sanfilippo Syndrome Type B": *PLoS ONE*: 8 (e772): 2007 (*rocking*)



# Mouse Rota-Rod

Cat. No. 47600

## General

The "Rota-Rod" technique has been originated by a 1957 paper of N.W Dunham and T.S Miya and has proved to be of great value in research involving screening of drugs which are potentially active on motory coordination.

The **Ugo Basile Rota-Rods**, the original Rota-Rods, are the result of many years of research in cooperation with the latest development in behavioral and pharmacological research.

The "Rota-Rod" technique is essential in screening drugs which have side effects on motor coordination.

The new Rota-Rod for Mice has five lanes, and operates in constant speed, accelerating, and reverse/rocking modes! The tiltable digital panel shows accrued rotations and time on the rod.



- THE ORIGINAL ROTA-ROD:  
THOUSANDS OF CITATIONS SINCE THE 1960s!
- THREE OPERATION MODES:  
CONSTANT SPEED, ACCELERATING OR ROCKING



## Main Features

- Adjustable speed (2-80 rpm) and acceleration ramp (6" - 600")
- Tiltable graphic display for optimal reading
- PC Interface: serial and USB (via the adaptor provided)
- Computer compatibility: direct connection to PC (via the **dedicated software included** as standard)

## Basic Features

The new Ugo Basile Rota-Rod replaces both previously available constant speed and accelerating models. It basically consists of five 3cm diam. drums, suitably machined to provide grip. Six flanges divide five 5.7cm lanes, enabling **five** mice to be on the treadmill simultaneously.

When a mouse falls off its cylinder section on to the plate below, the plate trips thereby recording the animal's endurance time in seconds. Height to fall is 16cm.

A large, very readable backlit graphic display shows the actual angular speed (RPM). At the end of a run, the display shows for each animal the running time and the instrument rotation speed at the time that animal fell off.

The panel can be oriented to select the most comfortable angle for the operator, to avoid glare, etc.

The main features, digitally preset by function keys are the following:-

- the angular speed can be preset in the range 2-80 RPM (revolutions per minute)
- in the acceleration mode, the change of speed can be preset in 6 second-10 minute interval in 6 second steps.
- reverse rotation can be selected, which takes place at minimum speed, at the end of a programmable acceleration-deceleration sequence.
- a rocking motion is also presettable, with adjustable angular amplitude, speed and acceleration.

## Data Acquisition

The 47600 is microprocessor controlled and features direct PC output. Internally stored data can be routed to the PC serial port (RS232) or USB (via adaptor provided).

Data output is managed by 52050-07 Data Acquisition Software Package (Windows® based), which enables the research worker to store the data into individual files, ready to be easily managed by most statistical analysis packages available on the market. **Ask for details!**

## 47800 COMBO-PACKAGE FOR MOUSE & RAT

If you plan to work with both rats and mice, you should consider the Combination Package 47600 (Mouse) + 47700 (Rat) Rota-Rod, offered at a special price.



47700 Rat Rota-Rod

## Ordering Information

**47600** **MOUSE ROTA-ROD**, standard package, including:

**47600-320** Trip Plate, complete (5 pieces)

**47600-321** Transmission Belt (2 pieces)

**47600-302** Instruction Manual (on USB key)

**52050-07** Dedicated Software Package (on USB)

**52050-322** Serial Cable

**52010-320** Serial to USB adaptor

**E-WP008** Mains Cord

## Optional

**57145** Thermal Mini-Printer

**47800** Combination Package 47600 Mouse Rota-Rod and 47700 Rat Rota-Rod

## Physical

Universal input 85-264 VAC, 50/60 Hz

Dimensions 40(w)x30(d)x38(h)cm

Weight Kg 6.0

Shipping Weight Kg 12.0 (approx.)

Packing 66x50x63cm

## Bibliography

### Method Papers

- N.W. Dunham & T.S. Miya: "A Note on a Simple Apparatus for Detecting Neurological Deficit in Rats & Mice" J. Am. Pharmaceut. Assoc., Scientific Edit., XLVI: No. 3, 1957
- B.J. Jones & D.J. Roberts: "The Quantitative Measurement of Motor Incoordination in Naive Mice Using an Accelerating Rotarod" J. Pharm. Pharmacol.: 20: 302-304, 1968

### Papers Dealing With Rota-Rod Technique

- B.J. Turner et alia: "Overexpression of Survival Motor Neuron Improves Neuromuscular Function and Motor Neuron Survival in Mutant SOD1 Mice" Neurobiol. of Aging 35 (4): 906-915, 2014
- M. Milanese et alia: "Knocking Down Metabotropic Glutamate Receptor 1 Improves Survival And Disease Progression in the SOD1<sup>G93A</sup> Mouse Model of Amyotrophic Lateral Sclerosis" Neurobiol. of Disease 64: 48-59, 2014
- J.E. Lorenz: "Oxidant-Induced Activation of cGMP-Dependent Protein Kinase I  $\alpha$  Mediates Neuropathic Pain After Peripheral Nerve Injury" Antioxidants & Redox Signaling Jan. 2014
- J.N. Justice et alia: "Battery of Behavioral Tests in Mice that Models Age-Associated Changes in Human Motor Function" Age Oct 2013
- R. Günther et alia: "Clinical Testing and Spinal Cord Removal in a Mouse Model for Amyotrophic Lateral Sclerosis (ALS)" JoVE J. Vis. Exp. (61), e3936, doi:10.3791/3936, 2012
- C.D. Heldermon et alia: "Development of Sensory, Motor & Behavioral Deficits in the Murine Model of Sanfilippo Syndrome Type B": PLoS ONE: 8 (e772): 2007 (*rocking*)

# Rotometer

Cat. No. 43000

## General

The Rotometer is widely used in research on motor assessment tests, in traumatic and acquired brain injury research and spinal cord injury research.

There are several well-characterized causes for animals to exhibit rotational behavior:

- Uneven/unilateral higher expression of levels of neurotransmitters (such as GABA or dopamine). Some brain tumors can cause aberrant expression levels to occur. Injury may also interfere with proper neurotransmitter expression and/or cause some localized change in neurotransmitter expression.
- Developmental anomalies can also cause rotational behavior.
- Anxiety/stress may cause this aberrant behavior.
- Exposure to some drugs, or drug abuse, or withdrawal from some drugs; all may cause rotational sequences.
- Physical lesions also can cause rotational behavior in an animal



**No Tether !**

**No Jacket !**

**TRULY  
UNRESTRAINED  
MICE**

## Main Features

- No jacket or tether is necessary: the animal is completely free
- Stand-alone, with internal memory
- Quick and simple to use: no training, turn-key system with software included

## Freely Moving Animals

To quantify rotational behavior in a freely moving mouse is a significant development.

This **new Rotometer** accomplishes this task precisely, using new and clever technology to count clockwise (CW) and counterclockwise (CCW) rotations in an open field.

The animal just carries a small magnet (not much larger than a grain of rice) on its nape or on its tail.

The magnet can be surgically implanted or injected subcutaneously; however, a convenient method is to attach it to the base of the mouse tail by using standard laboratory tape. This easy and efficient method, involves minimal stress for the animal, and has the advantage of requiring no anesthesia procedure.



Fig. 1: "2x15mm magnet, attached to the mouse tail"

Our **magnets** are encapsulated within a proven **bio-compatible material** (Paryline), to be implanted or injected subcutaneously, and fit within syringes normally used for the injection of identification transponders.



Fig. 2: "four Rotometers set up for high throughput screening, for testing several animals at the same time"

## Principle of Operation

The animal is placed in the open field (20cm diam. circular arena, enclosed in a 25cm tall acrylic cylinder). Our Rotometer is dimensioned for mice, but small rats can also be tested conveniently.

The design of this detecting system is very advanced, to enable the arena to be quite large whilst the magnet aboard is very small.

When the mouse circles within the open field, or rotates in place, the magnet (carried by the mouse) also rotates.

Sensors below the open field pick up these rotations, and the electronics record their number over time, discriminating Clockwise from Counterclockwise rotation.

As CW and CCW rotations accrue, they are displayed on the front panel and stored in the instrument internal memory; experiments may be qualified with animal data, date, time, and other diagnostic data.

## Data Acquisition

The 43000 is a microprocessor controlled unit. The experimental data, stored in its internal memory can be directly exported to the PC USB port, or to a flash drive (included).

Communication is managed by the dedicated CUB Data Acquisition Software Package, **Cat. 52050-13**, included as standard. The CUB Windows®-based Software Package enables the user to route the experimental data to the PC and store them into individual files, to be managed by most statistical analysis packages available on the market.

### Ordering Information

**43000 ROTOMETER**, standard package, including:

**43000-001 Main Unit with display**

**35100-286** Perspex Animal Restraint (25cm h)

**43000-302** Instruction Manual (on USB key)

**E-E 018** Paryline-coated Magnet, 2x12mm (2pcs)

**E-E 019** Paryline-coated Magnet, 2x15mm (2pcs)

**E-AU 041** Memory Key

**52050-13** CUB Data Acquisition Software Package and USB cable

#### Optional:

**57145** Thermal MiniPrinter

**43000-321** Syringe Kit, incl. implant, replacement needle & injectable magnets, 2x12 & 2x15 mm, 10 each

43000-012 Set of 10 Paryline-coated Magnets (2x12mm)

43000-015 Set of 10 Paryline-coated Magnets (2x15mm)

43000-052 Set of 50 Paryline-coated Magnets (2x12mm)

43000-055 Set of 50 Paryline-coated Magnets (2x15mm)

#### Specifications:

Read-out multifunction graphic display  
Print-out by optional thermal MiniPrinter  
Universal Mains 85-264 VAC - 50-60Hz - 30 W max.

Dimensions 25(w)x37(d)x16(h)cm, plus restraint  
Animal Restraint 20 (diam.) x 25 (h) cm

Weight 3.5Kg  
Shipping Weight 7.0Kg approx.  
Packing 65x34x28cm



# Hole Board

Cat. No. 6650

## General

The Hole-Board 6650 has been conceived to study the innate **exploratory behavior** of the mouse confronted with a new environment (head plunging stereotype), according to the classic method devised by Boissier-Simon.

The normal mouse of either gender, when confronted with a new environment, will explore holes in the substrate of its environment by **poking its nose** in and out of the hole a few times, then moving on to the next hole.

The initial exploration activity of the animal and its variations brought about by psychotropic drugs are unmistakably assessed. The nose poke frequency provides an indicator of exploratory behavior.

The test lasts few minutes and does not require any previous training/conditioning of the animal.



- Quick Test for Exploratory Behavior in Mice

- The classical "Planche à Trous" Test by Boissier & Simon

## Main Features

- The recording of the "nose poking" stereotype takes place automatically
- A few minute test is sufficient for most screening tests
- No previous training/conditioning required

## Instrument Descriptions

The “Méthode de la Planche à Trous” devised by Boissier & Simon (see bibliography) can be performed under optimum conditions: the **recording** of the “**head plunging**” or “**nose poking**” **stereotype** takes place automatically, via miniature I.R. emitters/receivers embodied in the “holes”.

The instrument consists of a “Board” and a Control Unit.

### Control Unit 6651

The control unit is lodged into a resilient cabinet whose front panel features the ACTIVITY display, the RESET and TEST keys, the LED visual indicators.

At every head plunging, the ACT (activity) LED blinks and the read-out increases by one digit.

A time-constant has been provided to inhibit the circuit to record a rapid up & down nose poking as it were a multiple event.

The figure remains frozen until the operator depresses the reset key again, when placing a fresh mouse on the board.

### Board 6652

The 40x40 cm board, 2.2 cm thick, is made of grey Perspex. The matt finishing avoids reflections which may alter the behavior of the animal.

The board embodies 16 “head-plunging detectors”, each comprising an I.R. emitter and a diametrically opposed receiver, flush mounted 1 cm below the upper panel.

The dimensioning of the board and holes has been optimized for mice in the 15-30 g range, to provide negligible false recordings.

## Data Acquisition

The 6650 Hole Board is provided with a connector for branching it to the **MULTIFUNCTION PRINTER Cat. 2600**, a microprocessor controlled device designed to acquire data from 6 independent channels.

The data, stored in the 2600 internal memory and shown on its graphic display, can be printed out in real time and/or routed to the PC, via the **52050-01 DAS Software Pack-age** provided with the 2600 package.

The **52050** is a Windows® based Data Acquisition Software Package, which enables the research worker to store the data into individual files, ready to be easily managed by most statistical analysis packages available on the market.

## Special Model for Videotracking

A special model of Mouse Hole-Board is also available, with no electronics, ideal for Videotracking.

The **46653** is a simple open field, dimensioned 40x40cm, with 16 holes diam 3cm, spaced 10cm apart (from center to center), enclosed in transparent (or opaque) walls.

A similar model, the **46652**, is also available, dimensioned 1mx1m, 35cm high, 16 holes diameter 3.8cm, to test rat exploratory behavior.

In alternative, the standard **6650** can be connected to a PC parallel port, via the optional cable **46650-330**, to acquire the number of pokes via **ANY-maze** and combine them with data collected by the videotracking software.

### Ordering Information

|                |   |
|----------------|---|
| <b>6650</b>    | <b>HOLE BOARD, standard package</b> , including:- |
| <b>6651</b>    | Control Unit                                      |
| <b>6652</b>    | Board   |
| <b>6655</b>    | Instruction Manual (on CD)                        |
| <b>E-WP008</b> | Mains Cable                                       |

### Physical

|                 |  |
|-----------------|--|
| Power           | 115 or 230V, 50/60Hz, 15W max.                       |
| Dimensions      | 40x40x2.2(h)cm (board)<br>26x15x25(h)cm (controller) |
| Weight          | 5.5Kg  |
| Shipping Weight | 10Kg approx.   |
| Packing         | 67x42x53cm   |

### Bibliography

#### Method Paper

- J.R. Boissier et P. Simon: “**Dissociation de deux composantes dans le comportement d’investigation de la souris**” *Arch Int. Pharmacodyn* 147, No. 3-4, 1964.
- J.R. Boissier et P. Simon: “**L’utilisation d’une réaction particulière de la souris (Méthode de la planche à trous) pour l’étude des médicaments psychotropes**” *Thérapie XIX*, 571-589, 1964.

#### Papers Mentioning 6650

- O.D. Can et alia: “**Anti-depressant-like effect of vitexin in BALB/c mice and evidence for the involvement of monoaminergic mechanisms**” *Eur. J. Pharmacol* 699 (1-3): 250-257, 2013
- M.G. Blake et alia: “**Scopolamine prevents retrograde memory interference between two different learning tasks**” *Physiology & Behavior* 102 (3-4): 332-337, 2011
- P.B. Gomes et alia: “**Anxiolytic-like effect of the monoterpene 1,4-cineole in mice**” *Pharmacology, Biochemistry & Behavior* 96: 287-293, 2010
- S.W. Zhu et alia: “**Influence of environmental manipulation on exploratory behaviour in male BDNF knockout mice**” *Behavioral Brain Res.* 197, 339-346, 2009

# Rotating Wheels for Rodent Activity

Cat. No. 1800 / 1850


**EASY MONITORING OF  
RODENT MOTOR  
ACTIVITY**

Data Acquisition  
available as optional  
(2600 Multifunction Printer)

## General

The Activity Wheels are designed to provide an easy and convenient method for measuring motor activity over long periods of time in laboratory rodents.

Especially useful for research on circadian rhythms or motor function, when connected to the 2600 Multifunction Printer or to any other data acquisition systems.

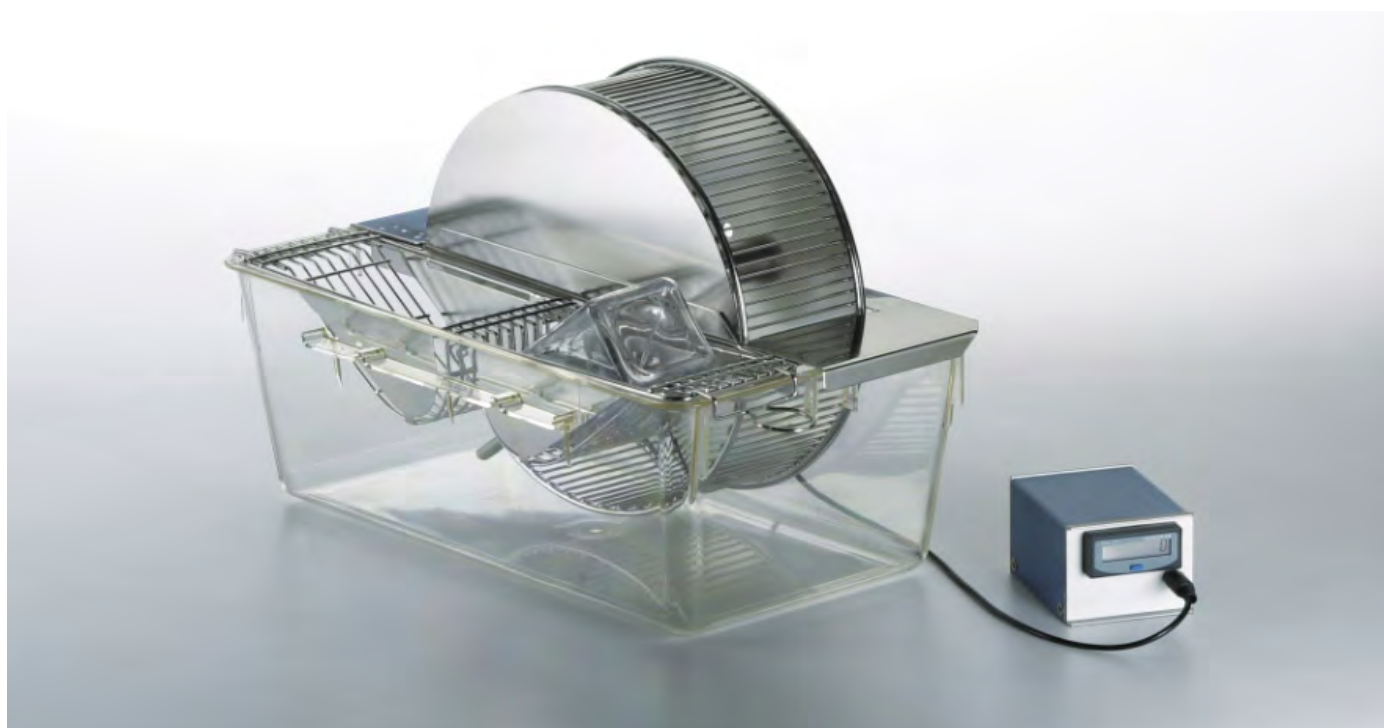

## Main Features

- Flexibility: version for rats or mice
- Easy monitoring (compatible with any Data Acquisition System)
- All stainless-steel wheel construction
- Clear polycarbonate cage for total visibility



## 1850 Mouse Cage

The 11850 is the classic **25 cm diameter running-wheel** made of stainless steel, provided with low friction Teflon bushing, for quite smooth action. The mouse runs on 2mm bars, placed 7 mm apart.

The wheel is housed in a clear polycarbonate cage. A stainless steel wire lid with exclusive lid locks, incorporates a 500 ml water bottle and a U-shaped food hopper for pellets.

The **Mouse cage is dimensioned 37(h)x26(w)x358d) cm.**



## 1800 Rat Cage

The Rat Cage is similar to the mouse model; the **running wheel has 35 cm diameter**. The 2 mm bars are placed 8.8 mm apart.

Dimensions of the **Rat Cage are 48(h)x32(w)x47(d) cm.**

## Revolution Counter

Each cage is complete with magnetic switch and LCD counter. The switch counts whole revolutions of the activity wheel and operates on an extended-life battery (included).

Cages without counter, models 1800-S and 1850-S, are also available, for data collection via PC, see paragraph below.

## Data Acquisition

For data acquisition a Multifunction Printer is required.

This is a microprocessor controlled device, designed to acquire data from 6 Cat. **2600**) independent channels (each Activity Wheel requires 1 channel).

The data, stored in the 2600 internal memory and shown on its graphic display, can be printed out in real time and/or routed to the PC, via the CUB software provided as standard.

When working with the Multifunction Printer, the counter is not required, so you may consider models **1800-S** or **1850-S**.



The picture above features a Multifunction Printer, with the necessary multi-connection cable 2610-F to connect up to 6 activity wheels.

## Ordering Information

- 1800 Rat Activity Wheel**, complete with polycarbonate cage, magnetic switch and LCD revolution counter
- 1850 Mouse Activity Wheel**, complete with polycarbonate cage, magnetic switch and LCD revolution counter
- 1800-S Rat Activity Wheel**, complete with polycarbonate cage & magnetic switch, without counter
- 1850-S Mouse Activity Wheel**, complete with polycarbonate cage & magnetic switch, without counter

## Multifunction Printers

- 2600 Multifunction Printer, 6 input channels**, with microprocessor for direct connection to the PC. Complete with dedicated software 52050-01, serial cable & USB adaptor

## 2610-F Multi-Connection Cable

### Physical

|                 |             |                      |
|-----------------|-------------|----------------------|
| Dimensions      | <b>1800</b> | 48(h)x32(w)x47(d) cm |
|                 | <b>1850</b> | 37(h)x26(w)x358d) cm |
| Weight          | <b>1800</b> | 7Kg                  |
|                 | <b>1850</b> | 5Kg                  |
| Shipping weight | <b>1800</b> | 11Kg                 |
|                 | <b>1850</b> | 7Kg                  |

# Mouse Ventilator

Cat. No. 28025

## General

This new Respirator, which completes the well known Ugo Basile line of Ventilators, features:-

- The **tidal volume**, in the range 0.1-1 ml (or 0.05-0.5 with the smaller piston installed), can be selected via its knob either while the pump is running or at a standstill. The stroke volume scale is ample, provided with precise engraved marks.
- The **rate**, selected by a knob, is indicated by a 3-digit solid state display, in the range 60-300 strokes per minute.
- Suitable channels and ports provide the witching of the air flow, with practically **no dead space**.
- A unique **variable stroke linkage** mechanism operates the piston.

The reciprocating motion is adjusted and transmitted to the piston by rods and articulated joints only, which leads to minimal wear, no backlash, silent operation and exact stroke reproducibility.



Unique Design

Reliable

Compact


Silent

## Main Features

- Ideal for use with mice, small birds and perinatal rats
- Optional 0.5 ml cylinder/piston assembly
- Purely mechanical, with impeccable finishing: lifetime lasting
- Quiet operation and negligible electrical noise

The instrument is compact and light, cm 20x13x18.5 and 2.5 Kg, and it is self-contained: in other words, it embodies its power supply which feeds the geared motor, its feedback controller and the rate display.

## The Connection Square

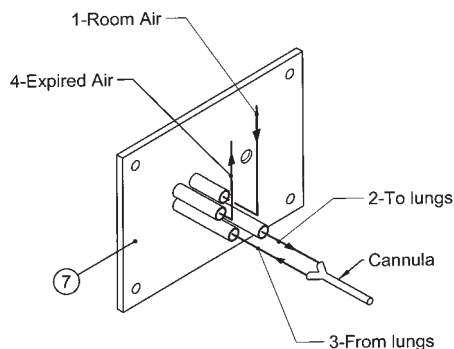

As illustrated in the drawing above, and pictured below, a connection square of four ports include:-

1. intake of air or other non-explosive gas mixture
2. delivery of air to the animal lungs
3. return air from animal
4. exhaust, for sampling, partial recycling, testing positive expiration pressure, etc.

so closely packed, that the connection tubes are cut in different lengths, to ease the insertion of the tubing.



## Start / Stop Model

A Mouse Ventilator version is available, Cat. **28125**, which embodies a controlled pause feature.

**The synchronised START/STOP function** gives the operator a means to stop and restart the respirator at "full lungs" point, via an external trigger pulse, when it is beneficial if not essential to minimize any extraneous movement of the anesthetized animal during electrophysiological recording, X-ray and imaging, etc.

## Specifications

|                    |  |
|--------------------|--|
| Rate               | 60 to 300 strokes for minute   |
| Rate Read-out      | on digital display   |
| Stroke Volume      | 0.1 to 1ml (with standard 1 ml piston)<br>0.05 to 0.5ml (optional 0.5ml piston)<br>Reproducibility $\pm 2\%$ |
| Volume Scale       | precision engraved, 0.05ml divisions   |
| Start-Stop         | by synchronised command<br>(model 28125 only)  |
| Power Requirements | 115 or 230V, 50/60Hz, 10W max.   |

## Physical

|                 |               |
|-----------------|---------------|
| Dimensions      | 20x13x18.5cm  |
| Net weight      | 2.2Kg         |
| Shipping Weight | 4.6Kg approx. |
| Packing         | 40x39x30cm    |

## Ordering Information

|                  |  |
|------------------|--|
| <b>28025</b>     | <b>MOUSE VENTILATOR</b> , complete with following standard accessories :-  |
| <b>28025-010</b> | 1ml Cylinder/piston assembly   |
| <b>28025-302</b> | Instruction Manual (on CD)   |
| <b>28025-321</b> | Perspex Vertical Lid   |
| <b>28025-323</b> | Cannula/Y-connection assembly (0.7mm & 1mm ID), tube, etc., in a plastic case  |
| <b>E-WP008</b>   | Mains Cord   |
| <b>Options</b>   |  |
| <b>28025-5</b>   | <b>Mouse Ventilator</b> , with 0.5ml cylinder/piston assy. & standard accessories  |
| <b>28025-005</b> | 0.5ml Cylinder/piston assembly   |
| <b>28125</b>     | <b>Mouse Ventilator</b> , with <b>synchronised START/STOP</b> feature, with 1ml cylinder/piston assy. & standard accessories   |
| <b>28125-5</b>   | <b>Mouse Ventilator</b> , with <b>synchronised START/STOP</b> feature, with 0.5ml cylinder/piston assy. & standard accessories |

## Bibliography

- M.S.Karbalaei et alia: "Impaired contractility and detrusor hypertrophy in cavin-1-deficient mice" *Eur.J.Pharmacol.* 689 (1-3): 179-185, 2012
- D. Chen et alia: "Angiotensin Type 1A Receptors in C1 Neurons of the Rostral Ventrolateral Medulla Modulate the Pressor Response to Aversive Stress" *J.Neurosc.* 32 (6): 2051-2061, 2012
- G. Grände et alia: "Unaltered size selectivity of the glomerular filtration barrier in caveolin-1 knockout mice" *Am.J.Physiol.* 297: F257-262, 2009
- S. Finotto et alia: "Asthmatic Changes in Mice lacking T-beta mediated by IL-13" *Intl. Immunology* 17, No. 8: 993-1007, 2005



# Rodent Ventilator

Cat. No. 7025

## General

The 7025 Rodent Ventilator is a volume-controlled mechanical ventilator (according to Starling's ventilator method), designed for use with rats, guinea pigs, mice and small birds.

The 7025 drive consists of a variable speed geared motor linked by a novel variable stroke mechanism to easily interchangeable cylinder/piston assemblies.

In particular, the **7025 can be equipped with 5, 10 or 30ml** cylinder/piston assembly.

Its precisely regulated geared-motor speed provides the most accurate and reliable stroke rate control of any respirator available.

The operation of the 7025 may be "paused" by an external TTL logic signal.

*The picture features a Rodent Ventilator 7025, together with the 6025 for Cat/Rabbit*



**Best available  
Starling  
Pumps**

**THE CHOICE OF  
THE CRITICS!**

*"We have four of your respirators in our extended lab and they are wonderful - as is your service"*

*Dr. Nicholas Price, Monash University*

## Main Features

- Interchangeable cylinder/piston assemblies (5, 10, 30ml)
- Quiet operation, both acoustically and electrically (negligible R.F. broadcasting)
- Reliable mechanics and impeccable finishing: lifelong lasting
- Synchronised START/STOP function available as optional

## Instrument description

The **unique linkage mechanism** insures that:

- 1) The piston almost touches the cylinder end with each stroke, regardless of the pre-set volume, thus insuring all air taken into the pump is expelled with each stroke.
- 2) The volume, clearly indicated on a **stationary dial**, is adjustable by means of a knob while the pump is either running or at standstill.
- 3) The reciprocating motion is generated, adjusted and transmitted to the piston by rods and articulated joints only.

The **lack of sliding friction** leads to:

- a) practically no wear
- b) no backlash and hence silent operation and exact stroke reproducibility.

## Hook-up to animal

Four ports (*Intake, To Animal, From Animal and Exhaust*) allow flexibility in air channelling.

The input may be room air or any non-explosive gas mixture. The exhaust air may be partially or totally recycled or collected for analysis.

## Ventilator Controls

The speed control knob adjusts the geared motor to the desired speed, which is indicated on the 3-digit LED display labelled STROKES P.M.

The operation of Ugo Basile Ventilators may be "paused" by an external TTL logic signal.

## Start / Stop Model

For more demanding electrophysiological-pharmacological investigations, in particular when the operation of the Ventilator is software controlled, a **synchronised command** is available to START-STOP the Ventilator at completed forced inspiration.

Ask for special models 7125.

## Specifications

|                     |   |
|---------------------|---|
| Rate                | 10 to 180 strokes for minute                                  |
| Rate Read-out       | digital display   |
| Stroke Volume       | 0.5 to 5; 1 to 10 or 3 to 30 ml, depending on cylinder/piston |
| Stroke Vol. Scale   | 1-10 ml   |
| Stroke Vol. Reprod. | ±2%   |
| Universal input     | 85-264 VAC, 50-60Hz, 40 VA max.                               |

## Physical

|                 |              |
|-----------------|--------------|
| Dimensions      | 27x26x19cm   |
| Net weight      | 9.5Kg        |
| Shipping Weight | 16Kg approx. |
| Packing         | 67x42x53cm   |

## Ordering Information

|                 |  |
|-----------------|--|
| <b>7025</b>     | <b>RODENT VENTILATOR</b> , complete with following standard accessories: |
| <b>7026</b>     | 10ml Cylinder/piston assembly, complete                                  |
| <b>7032</b>     | Perspex Lid  |
| <b>7033</b>     | Lithium-Grease Tube  |
| <b>7044</b>     | Y-Canula   |
| <b>7025-302</b> | Instruction Manual (on CD)   |
| <b>E-WP 008</b> | Mains Cord   |

## Other available models and accessories

|                 |   |
|-----------------|---|
| <b>7025-5</b>   | <b>RODENT VENTILATOR</b> , as above, 5ml  |
| <b>7025-30</b>  | <b>RODENT VENTILATOR</b> , as above, 30ml |
| <b>7128</b>     | 5ml Cylinder/piston assembly, complete    |
| <b>7027</b>     | 30ml Cylinder/piston assembly, complete   |
| <b>7025-150</b> | Anesthesia Kit                            |

## Models with synchronised START/STOP feature

|                |                                 |
|----------------|---------------------------------|
| <b>7125</b>    | <b>Rodent Ventilator</b> , 10ml |
| <b>7125-5</b>  | <b>Rodent Ventilator</b> , 5ml  |
| <b>7125-30</b> | <b>Rodent Ventilator</b> , 30ml |

See also our Anesthesia Systems, series 21100, the ideal match to our Ventilators!



## Bibliography

- R. Rudz et alia: "Acute Myocardial Infarction Inhibits the Neurogenic Tachycardic and Vasopressor Response in Rats via Presynaptic Cannabinoid Type 1 Receptor" *J.Pharmacol.Exp.Therap.* 343: 198-205, 2012
- A. Andersen et alia: "Ischemic Preconditioning Reduces Right Ventricular Infarct Size through Opening of Mitochondrial Potassium Channels" *Cardiology* 123 (3): 177-180, 2012
- Y. Kalam et alia: "Levosimendan Does Not Improve Cardiac Output or Blood Pressure in a Rodent Model of Propranolol Toxicity When Administered Using Various Dosing Regimens" *Intl. J. Toxicol.* 31 (2): 166-174, 2012
- D. Bilbija et alia: "Protecting the heart through delivering DNA encoding for heme oxygenase-1 into skeletal muscle" *Life Sciences* 91 (17-18): 828-836, 2012
- G.E. Peoples et alia: "Autologous pump-perfused rat hindlimb preparation for investigating muscle function and metabolism in vivo" *Microcirculation* (accepted article), 2013

## Cat/Rabbit Ventilator

Cat. No. 6025

### General

The 6025 Cat/Rabbit Ventilator is a volume-controlled mechanical ventilator (according to Starling's ventilation method), designed for use with cats, rabbits and animals of similar size.

The 6025 drive consists of a variable speed geared motor linked by a novel variable stroke mechanism to easily interchangeable cylinder/piston assemblies.

In particular, **the 6025 can be equipped with 50 or 100ml cylinder/piston assembly.**

Its precisely regulated geared-motor speed provides the most accurate and reliable stroke rate control of any respirator available

The operation of the 6025 may be "paused" by an external TTL logic signal.

*The picture features a Rodent Ventilator 7025, together with the 6025 for Cat/Rabbit*



**Best available  
Starling  
Pumps**

**THE CHOICE OF  
THE CRITICS!**

### Main Features

- Interchangeable cylinder/piston assemblies (50 and 100ml)
- Quiet operation, both acoustically and electrically (negligible R.F. broadcasting)
- Reliable mechanics and impeccable finishing: lifelong lasting
- Synchronised START/STOP function available as optional



## Instrument description

The **unique linkage mechanism** insures that:

- 1) The piston almost touches the cylinder end with each stroke, regardless of the pre-set volume, thus insuring all air taken into the pump is expelled with each stroke.
- 2) The volume, clearly indicated on a **stationary dial**, is adjustable by means of a knob while the pump is either running or at standstill.
- 3) The reciprocating motion is generated, adjusted and transmitted to the piston by rods and articulated joints only.

The **lack of sliding friction** leads to:

- a) practically no wear
- b) no backlash and hence silent operation and exact stroke reproducibility.

## Hook-up to animal

Four ports (*Intake, To Animal, From Animal* and *Exhaust*) allow flexibility in air channelling.

The input may be room air or any non-explosive gas mixture. The exhaust air may be partially or totally recycled or collected for analysis.

## Ventilator Controls

The speed control knob adjusts the geared motor to the desired speed, which is indicated on the 3-digit LED display labelled STROKES P.M.

The operation of Ugo Basile Ventilators may be "paused" by an external TTL logic signal.

## Start / Stop Model

For more demanding electrophysiological-pharmacological investigations, in particular when the operation of the Ventilator is software controlled, a **synchronised command** is available to START-STOP the Ventilator at completed forced inspiration.

Ask for special models 6125.

## Specifications

|                     |   |
|---------------------|---|
| Rate                | 10 to 100 strokes for minute                                |
| Rate Read-out       | digital display   |
| Stroke Volume       | 10 to 50; 20 to 100, depending on cylinder/piston installed |
| Stroke Vol. Scale   | 10-50 ml  |
| Stroke Vol. Reprod. | ±2%   |
| Universal input     | 85-264 VAC, 50-60Hz, 40 VA max.                             |

## Physical

|                 |              |
|-----------------|--------------|
| Dimensions      | 27x26x19cm   |
| Net weight      | 10.5Kg       |
| Shipping Weight | 16Kg approx. |
| Packing         | 67x42x53cm   |

## Ordering Information

|                 |  |
|-----------------|--|
| <b>6025</b>     | <b>CAT/RABBIT VENTILATOR</b> , complete with following standard accessories: |
| <b>6026</b>     | 50ml Cylinder/piston assembly, complete                                      |
| <b>6027</b>     | Set of 2 Lip-Seal Rings for 50ml piston                                      |
| <b>7032</b>     | Perspex Lid  |
| <b>7033</b>     | Lithium-Grease Tube  |
| <b>7034</b>     | Set of 3 Hex. Wrenches (2, 2.5, 3 mm)  |
| <b>6044</b>     | Y-Canula   |
| <b>6025-302</b> | Instruction Manual (on CD)   |
| <b>E-WP 008</b> | Mains Cord   |

## Other available models and accessories

|                 |  |
|-----------------|--|
| <b>6025-100</b> | <b>Cat/Rabbit Ventilator</b> , as above, 100ml |
| <b>6029</b>     | Set of 2 Lip-Seal Rings for 100ml piston       |
| <b>6025-150</b> | Anesthesia Kit                                 |

## Models with synchronised START/STOP feature

|                 |                                      |
|-----------------|--------------------------------------|
| <b>6125</b>     | <b>Cat/Rabbit Ventilator</b> , 50ml  |
| <b>6125-100</b> | <b>Cat/Rabbit Ventilator</b> , 100ml |

See also our **Anesthesia Systems, series 21100**, featured in the picture together with a 6026 Ventilator.



The ideal match to our Ventilators!

## Bibliography

- A. Ahmed et alia: "Development of an In Vitro Model to Assess Deposition of Aerosol Particles in a Representative Replica of the Rat's Respiratory Tract" *J. of Aerosol Med.* 25 (3): 169-178, 2012
- L. Monassier et alia: "Prevention by NMDA receptor antagonists of the centrally-evoked increases of cardiac inotropic responses in rabbits" *Br. J. Pharmacol.* 111 (4): 1347-1354, 2012
- T. Hoch et alia: "Modulation of the amplitude of  $\gamma$ -band activity by stimulus phase enhances signal encoding" *Eur. J. Neuroscience* 33 (7): 1223-1239, 2011
- T. Tchumatchenko et alia: "Ultrafast Population Encoding by Cortical Neurons" *J. Neuroscience* 31 (34): 12171-12179, 2011

## Bronchospasm Transducer

Cat. No. 7020

### General

This transducer is designed to perform the bronchospasm test on the guinea pig and is particularly suitable for connecting to UGO BASILE DataCapsule-Evo Recorder.

It enables the research worker to evaluate the spasm-inducing effect of drugs having a very wide range of action, not necessarily intended to act on respiratory dynamics.

The Bronchospasm Transducer 7020 is also a useful research tool for screening substances inducing the opposite effect, both those causing active bronchodilation in basal conditions and those which antagonize test drugs such as histamine, bradykinin, etc.

It is basically an air flow meter provided with a water input valve with adjustable pressure threshold.

The whole device is a compact unit made entirely of Perspex, mounted on a base along with its own power supply and controls.

The picture shows a complete set-up for bronchodynamics studies which includes the Rodent ventilator 7025 and a pen recorder



- Evaluates the bronchospasm inducing effect of drugs


### Main Features

- Simple and reliable method to assess airflow resistance
- The effect of bronchodilators agents is quickly assessed (by simply connecting to an animal ventilator and to a data acquisition system or chart recorder)

## Experimental Layout

The experimental layout follows the well-known Konzett-Roessler arrangement (see BIBLIOGRAPHY) with the anaesthetized guinea pig breathing via a reciprocating pump, according to Starling's mode of operation.

See sketch below:

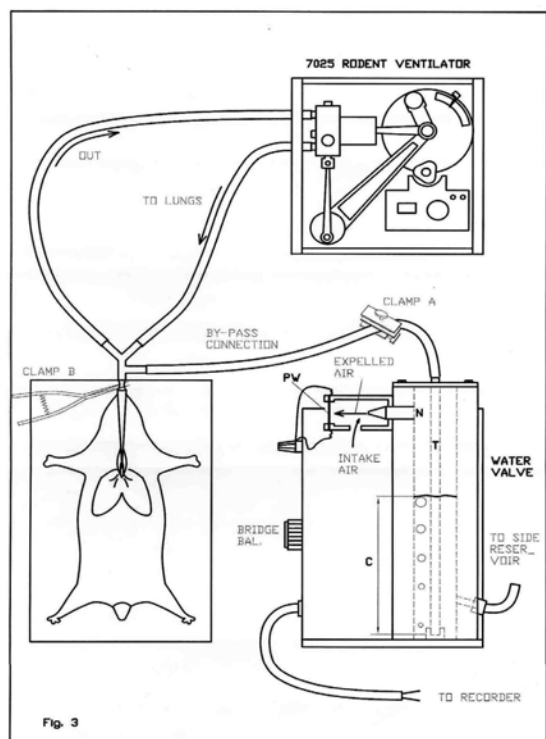

## Sensitivity

The sensitivity of the instrument in comparison with conventional Konzett-Roessler apparatus as illustrated in the table below

Minimum dosage in  $\mu\text{g/Kg}$  giving significant readings

|               | K-R Apparatus | UGO BASILE 7020 |
|---------------|---------------|-----------------|
| Histamine     | 3 - 6         | 0.3 - 0.6       |
| Acetylcholine | 20 - 40       | 3 - 10          |
| Serotonin     | 6 - 15        | 1 - 3           |

## Air Flow Meter

The recording system monitors respiratory dynamics by providing a tracing appearing as a succession of spikes.

When bronchospasm occurs, overpressure displaces the water column inside the T-tube and air bubbles through the water, escaping through an air flow transducer which are generating an electrical signal.

When Bronchodilators are administered, overpressure is reduced to below normal breathing values, as the bronchi exert less aerodynamic resistance to forced inspiration.

The tracing will decrease in amplitude to a marked degree, enabling the action of bronchodilators to be assessed.

## Ordering Information

**7020** Bronchospasm Transducer, complete with following parts:

- 7021** Air Flow Head
- 7022** Set of two Steel Rods
- 7024** Side Reservoir
- 7015** Instruction Manual (on CD)

Ask for details about:

**7025** Rodent Ventilator

## Physical

|                 |            |
|-----------------|------------|
| Weight          | 2.7Kg      |
| Shipping Weight | 5.2Kg      |
| Packing         | 40x39x30cm |

## Bibliography

### Method Paper

- H.Konzett & R. Roessler: Arch. Exp. Path. Pharmacol.: 195, 171, 1940

### Papers which include mention to 7020

- K. Ogino et alia: "Anti-inflammatory Effect of Arginase Inhibitor and Corticosteroid on Airway Allergic Reactions in a Dermatophagoides farinae-induced NC/Nga Mouse Model" *Inflammation* 36 (1): 141-151, 2013
- S.J.S. Flora et alia: "Interactive effect of arsenic and fluoride on cardio-respiratory disorders in male rats: possible role of reactive oxygen species" *BioMetals* 24 (4): 615-628, 2011
- N.R.F. Nascimento et alia: "1,8-Cineole induces relaxation in rat and guinea-pig airway smooth muscle" *J. Pharmacy and Pharmacol.* 61 (3): 361-366, 2009
- A. Floch et alia: "Characterization of NK1 and NK2 tachykinin receptors in guinea-pig and rat bronchopulmonary and vascular systems" *Br. J. Pharmacol.* 111 (3): 759-768, 1994
- W.A.K. Lau et alia: "Methoxyphenamine inhibits basal and histamine-induced nasal congestion in anaesthetized rats" *Br. J. Pharmacol.* 101 (2): 394-398, 1990
- C. Broquet et alia: "Aminoacylates and Aminocarbamates of 2-substituted 4-hydroxymethyl 1,3-dioxolans as Ammonium Salts. A new Series of PAF Antagonists" *Eur. J. Med. Chem.*: 25: 235-240, 1990.

## Gas Anesthesia Systems

Cat. No. 21100

### General

The Ugo Basile New Gas Anesthesia is a compact, modular and reasonably-priced system, intended to match the highest technical requirements of animal labs that do not compromise on quality.

A wide range of options and accessories are available, most of which can be added in a scalable manner, making the system modular and with an excellent value for price!

Typical anesthesia procedures involve an induction phase and a maintenance phase, which require at least:

- Flow-meter and anesthetic Vaporizer
- Induction box and/or mask with breathing circuit
- Scavenger or flow hood (for gas anesthetic removal)

The Ugo Basile New Gas Anesthesia system include all of the above! ... and much more!



**Portable**

**Modular**

- **THE IDEAL MATCH TO UGO BASILE LINE OF VENTILATORS**

### Main Features

- Digital Flowmeter with wide range (up to 16 litres per minute) for multiple animal delivery
- Up to six Animals with one Station
- Manifold for mask/induction-box switch and full range of accessories
- NEW Tec3 Vaporizers (non-refurbished)



## Overview

The unique digital flowmeter, coupled to non-refurbished vaporizers for Isoflurane or Sevoflurane, result in an innovative yet sturdy and reliable system to anesthetize animals of virtually any size and up to 6 animals simultaneously.

An ample selection of modular components and accessories enables the user to customize and expand the anesthesia system upgrading from a **basic** (flowmeter & vaporizer) to a **full system** (with induction boxes, breathing circuits with masks of any size, switch valves, multiple delivery systems active or passive scavengers, etc.)

The blue 4mm thick aluminum rack has a highly resistant paint to protect against stains from aggressive anesthetic liquids & solvents.

Two universal attachment blocks are mounted on the back, to connect the device easily to any rail or mobile floor model anesthesia rigs of sizes 25x8mm up to 35x10mm.

## Digital Flowmeter

The Ugo Basile Gas Anesthesia System includes a unique digital flowmeter.

Its wide flow range (from 0.3 to 16 l/min.) and fine resolution (0.1 l/min.) guarantees enough gas flow to anesthetize up to 6 animals simultaneously!

Small and large animals could be anesthetized with the same system (virtually, from mouse to horse!)



## Nose-cone/Masks with diaphragm

Unlike many rodent masks available on the market, these masks incorporate a latex diaphragm, which holds the rodent nose, keeping the animal in correct position and ensuring a continuous positive flow of fresh oxygen & anesthetic.

The membrane also provides a positive seal reducing the exposure of the user to anesthetic gases.

Available in several sizes:

- Small/Large Mice
- Small/Medium/Large Rats
- Large Rodents/Feline



*The picture shows a mouse nose-cone/mask, connected to an evacuation tubing.*

## Induction Box

The **7900** Induction Box is a conveniently dimensioned (25x13x13cm), cost-effective solution to confine one guinea pig, one rat or several mice.

It incorporates a sliding lid and tubing connectors (vaporizer input and scavenger output).



A larger size, **7910** is also available, dimensioned 44x22x21 cm, for larger animals such as rabbits.

## Dual Diverter Manifold with Humidifier

All of the Ugo Basile Gas Anesthesia Systems come with a pre-installed mounting bracket to fit the Dual Diverter Manifold (as shown in the picture).



The anesthetic gas flow can be diverted toward 2 independent devices (i.e., an induction chamber and a breathing mask).

A simple and efficient humidifier is included with the manifold. It is especially recommended for long-term anesthesia, when dehydration may become an issue.

## Multiple Delivery System

The Multiple Delivery accessory allows the connection of up to six devices to one anesthesia system for simultaneous operation.



Each device (for 2, 3, 4, 5 or 6 animals) has independent flow regulation.

## F/AIR Scavenger

A solution to handling waste anesthetic gases when active evacuation systems are not available, activated charcoal canisters remove approx. 50g of halogenated anesthetic agents from the waste gas stream before being discarded.



## Ordering Information

### ANESTHESIA SYSTEMS

- 21050 Basic Single-Output Anesthesia System** including Digital Flowmeter (for O<sub>2</sub> or Medical Air) and TEC-3 vaporizer for Isoflurane (\*)
- 21100 Single-Output Anesthesia System**, including 21050 (\*), 2 passive scavengers (\*\*), evacuation tubing.
- 21200 Double-Output Anesthesia System**, including 21050 (\*), 4 passive scavengers (\*\*), evac. tubing & dual diverter manifold with humidifier
- 21400 Multiple-Animal Anesthesia System**, including 21050 (\*), 8 passive scavengers (\*\*), evac. tubing and Multiple Delivery System for 4 animals.
- 21600 Multiple-Animal Anesthesia System**, including 21050 (\*), 12 scavengers (\*\*), evac. tubing and Multiple Delivery System for 6 animals.

### ACCESSORIES

#### Delivery Systems (Masks & Induction Boxes)

- PS-0525-A Nose-Cone/Mask Circuit for Small Mice**,
- PS-0305-A Nose-Cone/Mask for Large Mice**, 3cmØ
- PS-0306-A Nose-Cone/Mask for Small Rats**, 4.5cmØ
- PS-0307-A Nose-Cone/Mask, Medium Rats**, 5cmØ
- PS-0308-A Nose-Cone/Mask for Large Rats**, 5.5cmØ

*All masks are complete with diaphragm and inlet connector*

- 7900 Induction Box for small rodents** (rats and mice), dimensioned 25x13x13 (h) cm
- 7910 Large Induction Box**, 40x22x21(h)cm
- 21100-790 Induction Box for small rodents**, airtight model, with latch, 25x13x13 (h) cm

#### Special Systems with N<sub>2</sub>O

- 22100 O<sub>2</sub>/N<sub>2</sub>O Anesthesia System**, with 2 Analog Flowmeters, TEC-3 vaporizer for Isoflurane (\*), passive scavenger (\*\*), evac. tubing.

\* Vaporizers for other anesthetic agents available on request

\*\* Activated Charcoal Canisters

### Multiple-Output Delivery Systems

**PS-0529-02 Dual Diverter Manifold** with humidifier, *see complete model 21200*

**PS 30-459 Multiple-Animal Delivery System**, 6 Flowmeters, *see complete model 21600*

*Multiple delivery systems for 2, 3, 4, and 5 animals available*

### Anesthetic Scavenger and Evacuation

**PS-0581-00 F/air filter** (activated charcoal canister)

**PS-0581-01 F/air filter**, pkg. of 8

**PS-0582 Evac.Tubing for F/air**, 1.8 m with 19 mm male x 22 mm female adaptor

**21100-833 Active Scavenger System**, to remove the anesthetic agent by negative pressure (to be connected to an activated charcoal canister)

### Heating Pads and Surgical Tables

**PS-0491 Rodent Workstation**, with heated surgical table and outlet for anesthetic removal. Includes connection tubing for Gaymar water circulator (not included)

**PS-0766 Water Circulating Pump** Gaymar TP702, adjustable from 10 to 42°C

**PS-0811 Heating Pads** Delta-Phase Isotherm (pkg of 3), 20x20x0.65 cm. Maintains animal body temperature near 37°C up to several hours. Ideal for NMR.

### Other Recommended Accessories

#### Fill Devices

- PS-0950** for Isoflurane
- PS-0951** for Sevoflurane
- PS-0949** for Halothane



#### Physical (21100)

|                 |                     |
|-----------------|---------------------|
| Weight          | 8.5Kg               |
| Dimensions      | 26(w)x18(d)x24(h)cm |
| Shipping Weight | 12Kg                |
| Packing         | 67x42x53cm          |

## Anesthetizing Box

Cat. No. 7900 (rodents) 7910 (rabbits)

### General

Our Induction Boxes are conveniently dimensioned induction boxes, featuring a sliding lid. They are made of Perspex and prove to be particularly useful to confine laboratory animals during anesthetizing.

The **7900**, for small rodents, is **dimensioned 25x13x13(h)cm**; the larger model 7910, for rabbits is **dimensioned 40x22x21(h)cm**;

The transparent acrylics permits the animal to be kept under constant observation.

Two tubing connectors of nickel plated brass are fitted into each end, one located at the top of the box and the other at the bottom.

Any (non-explosive!) gas mixture can be used. In case small quantities of liquid, as ether or chloroform are used, soak a cotton wool flock and place it in a small Becker, in-side the box.

For more demanding application, and higher safety, an airtight model, with latch, is also available, see picture below.



**Our Induction chambers  
are ideal to work with  
our new  
Anesthesia Systems**

**TO CONFINE SMALL  
LABORATORY ANIMALS  
DURING  
ANESTHETIZING**



### Ordering Information

- **7900** Induction box for small rodents  
25x13x13(h)cm, ID 23x12x12(h) cm
- **7910** Induction box for rabbits  
40x22x21(h)cm, ID 38x20x19(h) cm
- **2100-790** Airtight model, with latch, see picture  
25x13x13 (h) cm, ID 21x11x13(h) cm



## Touch-Screen Controller for Conditioning Cages

Cat. No. 40500-001

### General

The **40500-001 Touch-Screen Controller** is a powerful and versatile tool, which will function as main unit in a number of tests:

**Fear Conditioning** (see 46000-100)

**Passive Avoidance (step-through)**

**Passive Avoidance (step Down)**


**Active Avoidance**

**Learned Helplessness**

**Startle Response/PPI** (available soon)

For each test, a specific application software is available for installation; each software is sold separately, so it is easy to customize each controller.

By the application "Launcher UB" installed on the 12" touch-screen, the user chooses the experimental routine among the ones installed.

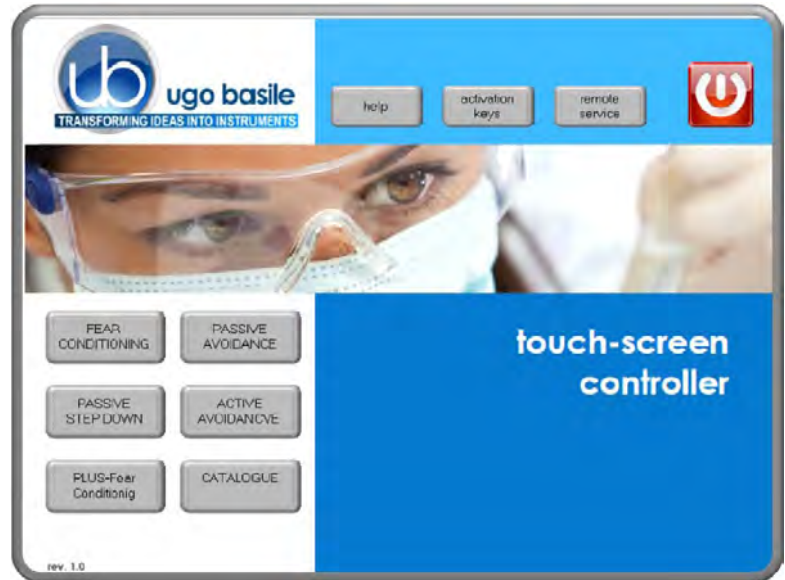

### A SINGLE UNIT TO CONTROL:

- EXPERIMENT SETTINGS (LIGHT, SOUND, ETC.)
  - SHOCK PARAMETERS
  - DATA ACQUISITION, MANAGEMENT & EXPORT
- IN ALL UB CONDITIONING CAGES!**

**Great Versatility**

**Outstanding Adaptability**

**tip!**

this is part of UB conditioning cage project; you buy a single touch-screen controller, and manage all UB cages.  
**Ask for details!**

### Main Features

- The electronic unit encompasses all controls for **up to 4 animal cages**
- The same controller will function as main unit in a number of **behavioral tests**; just purchase the hardware and the application software for the additional test!
- The new "**launcher**" application, makes it easy to manage other UB behavioral cages
- The **remote control** feature will make remote service and software upgrades extremely simple!



## System Description

Different set-ups, depending on animal (rat or mouse), type of behavioral test and number of cages, can be obtained by combining the following elements:

- **Programming/Recording Unit with Shocker**
- **Behavioral Cage/s** (up to 4 with one controller)
- **Expansion Box**, for multiple cage set-up
- **Isolation Cubicle/s Box**, (if required)

## Instrument Description

The **40500-001** Programming/Recording Unit, encompassing all controls, incorporates a constant-current high precision 8-pole shocker, and manages data acquisition: data are stored inside the unit and can be downloaded via the USB key provided as standard, for further processing via Excel, Access, etc.



The Unit, incorporating a 12" touch-screen, manages the test via the dedicated application software, see ordering info.

Up to 4 cages of the same type can be connected to the same Controller; in this case, an expansion box **40500-005** is required for each extra cage.

## Launcher Menu

Beside the buttons to select the desired application software, this menu features the following options:



- **Help** Pressing the "help" button will display the Launcher user manual online
- **Activation Keys** the software activation keys can be entered via a virtual keyboard. Additional software activation keys may be purchased separately
- **Remote Service** remote service is possible thanks to a specific software installed on the Touch-Screen Controller.

## Ordering Information

**40500-001** Touch-Screen Controller & Shocker

### Available Software and Activation Keys

- 40530-010** Activation of **Active Avoidance** Software
- 40550-010** Activation of **Passive Avoidance** (step-through) Software
- 40570-010** Activation of **Passive Avoidance** (step-down) Software
- 46100-010** Activation of **Fear Conditioning** Software

### See also the following datasheets

- 40530** **Passive Avoidance** (step-through)
- 40550** **Passive Avoidance** (step-through)
- 40570** **Passive Avoidance** (step-down)
- 46000** **Fear Conditioning**

## System Specifications

|                          |   |
|--------------------------|---|
| Inputs                   | 4   |
| Input voltage            | TTL input 0-5Vdc opto-isolated  |
| LCD                      | 12" with resistive touch screen   |
| CPU Module Port          | 2 USB Port 2.0<br>1 Ethernet port 10/100Mb<br>1 DVI port for external monitor |
| Peripheral Port          | 4 output for Sound, Shock and light   |
| Power supply             | 12V-2A  |
| Expansion Bus Connection | 2 RJ11 connectors   |
| USB port                 | type B ( only for software connection)  |

### Physical:

|                 |                           |
|-----------------|---------------------------|
| Weight          | 2.7Kg                     |
| Shipping Weight | 4Kg                       |
| Dimensions      | 25(d) x 33(w) x 5.5(h) cm |
| Packing         | 53x41x13cm                |



# New Fear Conditioning System

Series 46100

## General

The Ugo Basile Fear Conditioning Systems 46000 includes all the components to run experiments on mice or rats, according to the paradigms:

- **Contextual Fear Conditioning**
- **Cued Fear Conditioning**

The detection of **Freezing** is automated and based on video analysis. The **shock**, **light** and **sound** parameters are controlled by software (USB) or manually, via the new Electronic Unit, based on touch-screen technology.

## System Configuration

A typical **Basic System** consists of:

- Controller with touch-screen
- Animal box with electrified floor and Context Kit (3 floors, 9 walls)
- Isolation Cubicle, with dual (visible/I.R.) light, speaker and fan

The **complete system** also include:

- Freezing-detection Software
- USB Videocamera




Preinstalled PC can be supplied as optional



**Memory**

**Behaviour**

**NEW MODEL**  
all controls managed  
by a single  
**Electronic Unit**

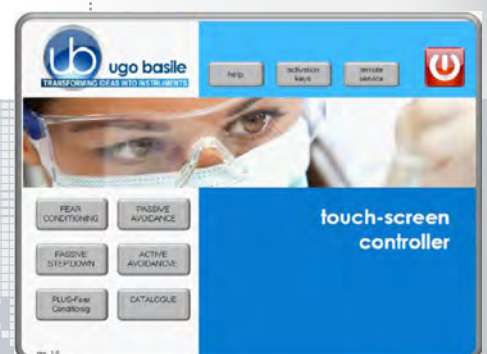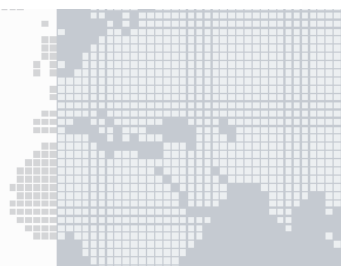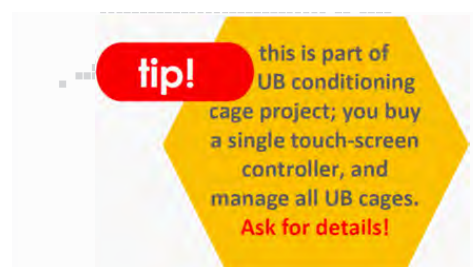

## Main Features

- **AUTOMATIC** detection of **FREEZING** also in **Total Darkness**
- Specific versions for rats or mice
- Multiple Cage Set-up (up to 16 cages, in groups of 4)

## NEW on the 2014 version!

- The new "launcher" application, makes it possible to manage other UB behavioral cages with the same **Touch-Screen Controller 40500-001**; just purchase the hardware and the application software for the additional test!
- **Remote Control feature** will make remote service and software upgrades extremely simple!

## System Components

### Software and IR-CCD camera


The Ugo Basile Fear Conditioning system benefits of a specific version of the Any-maze software, Cat. No. 60000-FC; the software controls the Ugo Basile hardware, automatically detects the *freezing* behavior and analyzes the results across time.



Measured parameters include:

- Total *Freezing* time
- No. and duration of freezing episodes across time

The USB videocamera 47400-030, is sensitive to IR light and allows for *freezing* detection even in total darkness.

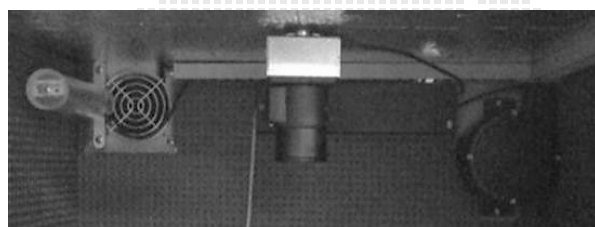

Wide angle lenses and IR filters are included.

### Controller

The **new** FC Controller 46000-100 consolidates all controls in a single compact electronic unit.

On its 12" touch-screen, the researcher sets the following parameters via the user-friendly interface

- **Sound**, in the range 100Hz-40KHz; 1-150dB or white noise. The speaker is included in the Cubicle.
- **Shock**: constant current (from 0.1 to 2.9 mA in 0.1 mA steps). The shock can be controlled via external operation (via 5V TTL signals)
- **Light**

Connections are arranged on the controller back panel:



### Animal Box with Electrified Grid Floor

- **46003** Mouse Box:  
inside dimensions: 17x17x25(h) cm
- **46002** Rat Box  
inside dimensions : 26x26x30(h) cm

### Context Kit


A complete set of removable contexts is provided to alter the colour of the box walls and floor. Each animal box includes a kit with: 3 striped walls, 3 chessboard, 3 grey walls and 3 plastic floors (white, black, grey).

## Isolation Cubicle

The new-design Isolation Cubicle 46000-590 includes:

- a dual (visible and I.R.) LED light
- a loudspeaker
- a noiseless fan

all conveniently positioned inside the sound-attenuating cubicle. Multiple-cage set-ups include an expansion-cubicle with its slave electronics on board

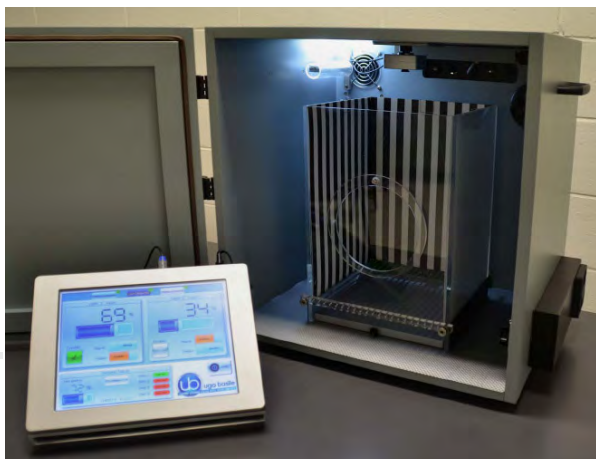

### Preinstalled PC (optional)

The Ugo Basile Fear Conditioning systems can be used with Windows-based laptop or desktop PCs, which can be purchased directly from Ugo Basile, with software & hardware pre-installed, tested, and ready to use.

## Ordering Information

### COMPLETE SYSTEMS (with software and USB camera)

| MOUSE | RAT   |                                |
|-------|-------|--------------------------------|
| 46153 | 46152 | Complete Single-Cage FC System |
| 46253 | 46252 | Complete Two-Cage FC System    |
| 46453 | 46452 | Complete Four-Cage FC System   |

### BASIC SYSTEMS (without software/camera)

| MOUSE | RAT   |                             |
|-------|-------|-----------------------------|
| 46103 | 46102 | Basic Single-Cage FC System |
| 46203 | 46202 | Basic Two-Cage FC System    |
| 46403 | 46402 | Basic Four-Cage FC System   |

Additional Animal Kits, including Cage, expansion cubicle & electronics (no camera) are available:

|           |       |
|-----------|-------|
| 46102-002 | Rat   |
| 46103-003 | Mouse |

All components can be ordered separately

## Bibliography

- S. Yusufisq et alia: "Post-Weaning Social Isolation impairs observational fear conditioning" *Behav. Brain Res.* 242 (1): 142-149, 2013
- A. Sirri *et al.*: "Temporal gene expression profile of the hippocampus following trace fear conditioning". *Brain Research* 1308, 14-23, 2010



## Active Avoidance Set-Up (Automatic Reflex Conditioner)

Cat. No. 40532 Rats

Cat. No. 40533 Mice

### General

The new model of **Active Avoidance Set-Up** has been designed to enable the researcher to perform a wide range of avoidance experiments, each according to a flexible schedule.

Via the **TIMELINE** feature, the user will be able to configure a number of different tests, according to the specific experimental needs, namely the classical shuttle-box tests in its various modes.

Ugo Basile Active Avoidance set-up instrument basically consists of a Controller, and a Cage for either rat or mouse.

The tests are conducted in a cage, divided into two sections by a partition with an intercommunicating opening at floor level.

The tilting floor ensures a simple and reliable detection mechanism to score the animal's movement across the two compartments.

The electronic unit encompasses all controls for up to 4 cages, and a scrambling shocker.



**NEW VERSION**

**Multiple-Cage  
Set-up**

**EFFICIENT, RELIABLE  
INSTRUMENT FOR  
THE CLASSIC ACTIVE  
AVOIDANCE TEST**

**tip!**

this is part of  
UB conditioning  
cage project; you buy  
a single touch-screen  
controller, and  
manage all UB cages.  
**Ask for details!**



### Main Features

- **Maximum flexibility:** configure your own Avoidance-Experiment Schedules via the **timeline** function
- The electronic unit encompasses all controls for **up to 4 animal cages!**
- **Reliable tilting-floor detection mechanism**

### NEW on the 2014 version!

- The new "launcher" application, makes it possible to manage other UB behavioral cages with the same **Touch-Screen Controller 40500-001**; just purchase the hardware and the application software for the additional test!
- **Remote Control feature** will make remote service and software upgrades extremely simple!




## Instrument Description

Different set-ups, depending on animal (rat or mouse) and number of cages, can be obtained by combining the following elements:

- **Programming/Recording Unit with Shocker**
- **Rat Cage** (up to 4 with one controller)  
or
- **Mouse Cage** (up to 4 with one controller)
- **Expansion Box**, for multiple cage set-up

## Programming/Recording Unit

The **40500-001** Programming/Recording Unit, encompassing all controls, incorporates a constant-current high precision 8-pole shocker, and manages data acquisition: data are stored inside the unit and can be downloaded via the USB key provided as standard, for further processing via Excel, Access, etc.

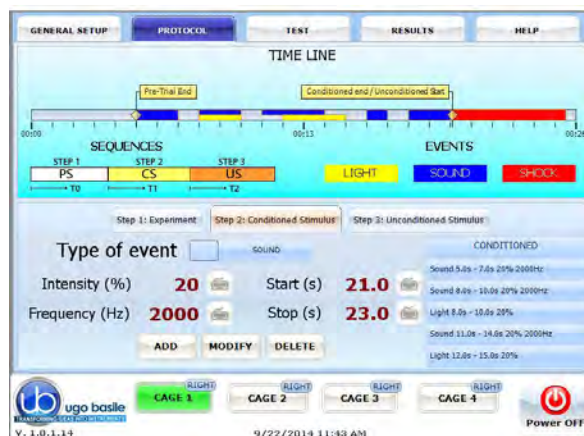

The Unit, incorporating a 12" touch-screen, manages the Passive Avoidance Test via the **40530-010** Software. Up to 4 cages can be connected to the same Controller. If more than one cage is connected, an expansion box **40500-005** is required for each additional cage.

The trials can be configured via the **TIMELINE** feature, entering the setting via the virtual keyboard: trial number, the acoustic/visual stimulus, delay, shock intensity, and timing of the different experimental sequences:

- PS:** pre-stimulus interval (randomizable)
- CS:** conditional stimulus interval
- US:** unconditional stimulus interval.

## Active-Avoidance Cage (shuttle-box)

Two types of cages are available:

- **40532** designed for **Rats**  
dimensioned 57x27x30(h)cm, I.D. 48x20x22(h)cm
- **40533** designed for **Mice**  
dimensioned 47x18x25(h)cm, I.D. 38x9x17(h)cm

Both cages are provided with acoustic and visual conditioning stimulators. The reinforcement consists of an electrical stimulus applied to the floor bars of the cage by an incorporated 8-pole "scrambling" circuit.

The cage is divided into two compartments intercommunicating by an opening at floor level.

When the animal crosses the door, the cage floor tilts, thus operating a reed arrangement, which cuts out all the stimuli or, if the crossing takes place during the pause, records the intertrial crossing.

## Ordering Information

- 40500-001** Programming/Recording Unit & Shocker
- 40530-010** P.A. Software and activation
- 40532** Rat Cage, complete with catch pan
- 40533** Mouse Cage, complete with catch pan
- 40500-005** Expansion Box, for multiple cage set-up

## Specifications :

- Shock Duration in steps of 0.1s
- Shock intensity 0-3mA step 0,1mA
- Light intensity 0-100%, in steps of 5
- Sound intensity 0-100%, in steps of 5
- Sound frequency 100-18.000Hz, in steps of 100Hz
- Light, sound, shock start in seconds, 0,1s precision
- Light, sound, shock stop in seconds, 0,1s precision

## Physical:

- Weight 2.7Kg (40500-001)  
5.3Kg (40532)  
3.4Kg (40533)
- Shipping Weight 4Kg (40500-001)  
9Kg (40532)  
5.8Kg (40533)

## Bibliography

### Papers which quote Ugo Basile A.A. Test (previous model)

- D. Dimitrova, D. Getova: "Effects of Rivastigmine on Learning and Memory Processes in Rats Active Avoidance Test" *Medicine* 4.1, 2014
- G.N. Carmona et alia: "The Dense Core Vesicle Protein IA-2, but not IA-2 $\beta$ , is Required for Active Avoidance Learning" *Neuroscience* 269 (6): 35-42, 2014
- O. Ortiz et alia: "Associative Learning and CA3-CA1 Synaptic Plasticity Are Impaired in D1R Null, Drd1a/ Mice and in Hippocampal siRNA Silenced Drd1a Mice" *J.Neuroscience* 30 (37): 12288-12300, 2010
- J.I. Lemos et alia: "Involvement of the prelimbic prefrontal cortex on cannabidiol-induced attenuation of contextual conditioned fear in rats" *Behav. Brain Res.* 207: 05-111, 2010
- N. Seferos et alia: "Mandibular bone density and calcium content affected by different kind of stress in mice" *J Musculoskelet Neuronal Interact.* 10 (3): 231-236, 2010

## Passive Avoidance Step-Through *New Model*

Cat. No. 40552 Rats

Cat. No. 40553 Mice

### General

**Passive Avoidance Test** is used to assess memory function based on the association formed between a specific environmental context, which the animal learns to avoid, and an aversive stimulus, represented by a mild foot shock.

The tests are conducted in a two-compartment apparatus, where one is dimly lit and preferable to a rodent, and the other is brightly lit.

After the training period, during the test proper, the animal that learned the task will avoid the location previously paired with the aversive stimulus, and show greater latency to enter it.

Ugo Basile Passive Avoidance set-up instrument basically consists of a Controller, and a Cage divided into two compartments by a partition which embodies a sliding door.

The tilting floor ensures a simple and reliable detection mechanism to score the animal's movement across the two compartments.



**Step-Through Cage**

**EFFICIENT, RELIABLE  
INSTRUMENT FOR  
THE CLASSIC PASSIVE  
AVOIDANCE TEST**

**Multiple-Cage Set-up**



### Main Features

- The electronic unit encompasses all controls for **up to 4 animal cages!**
- Silent and automated sliding door to divide the two compartments (no stepping motor!)
- Reliable tilting-floor detection mechanism

### NEW on the 2014 version!

- The new "launcher" application, makes it possible to manage other UB behavioral cages with the same **Touch-Screen Controller 40500-001**; just purchase the hardware and the application software for the additional test!
- **Remote Control feature** will make remote service and software upgrades extremely simple!


## Instrument Description

Different set-ups, depending on animal (rat or mouse) and number of cages, can be obtained by combining the following elements:

- **Programming/Recording Unit with Shocker**
- **Rat Cage** (up to 4 with one controller)  
or
- **Mouse Cage** (up to 4 with one controller)
- **Expansion Box**, for multiple cage set-up

## Programming/Recording Unit

The **40500-001** Programming/Recording Unit, encompassing all controls, incorporates a constant-current high precision 8-pole shocker, and manages data acquisition: data are stored inside the unit and can be downloaded via the USB key provided as standard, for further processing via Excel, Access, etc.

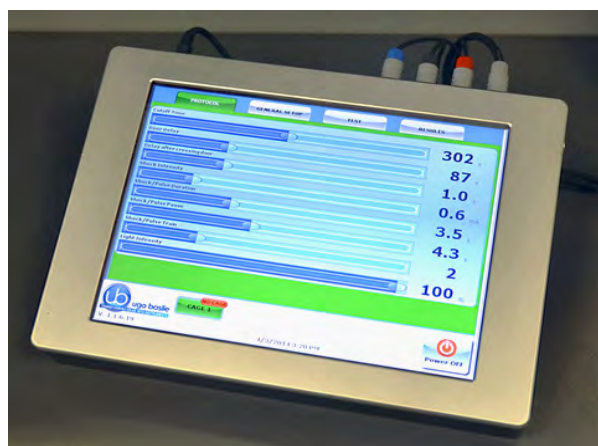

The Unit, incorporating a 12" touch-screen, manages the Passive Avoidance Test via the **40550-010** Software. Up to 4 cages can be connected to the same Controller. If more than one cage is connected, an expansion box **40500-005** is required for each additional cage.

## Passive-Avoidance Cage (step-through)

Two types of cages are available:

- **40552** designed for **Rats**  
dimensioned 57x27x30(h)cm, I.D. 48x20x22(h)cm
- **40553** designed for **Mice**  
dimensioned 47x18x25(h)cm, I.D. 38x9x17(h)cm

The cage is divided into two sections, the **START** and **ESCAPE** compartments. The start compartment is white and **illuminated** by a light fixture (3LED, white-light); the escape compartment is **dark** and its grid floor is connected to the shocker.

The two compartments are divided by a partition which embodies an automatically operated sliding door at floor level. The **door delay** and the **shock parameters** can be preset on the touch-screen of the controller, according to experience or data suggested by the literature.

With the rodent in the START compartment, the START button activates the timer, providing the **opening of the door** after the preset delay.

The opening of the door enables the **latency timer**, which stops at the animal crossing; latency time is displayed in 0.1s steps. The door shuts one second after the crossing, to prevent the animal being upset or hurt by a too close door operation.

## Ordering Information

- 40500-001** Programming/Recording Unit & Shocker
- 40550-010** P.A. Software and activation
- 40552** Rat Cage, complete with catch pan & sliding door assembly
- 40553** Mouse Cage, complete with catch pan & sliding door assembly
- 40500-005** Expansion Box, for multiple cage set-up

## Specifications :

|                 |                              |
|-----------------|------------------------------|
| Latency Time    | 5-digit Read-Out, 0.1s steps |
| Door Delay      | 0-99s, in steps of 1s        |
| Shock Duration  | 0.1-9.9s, in steps of 0.1s   |
| Shock Intensity | 0.1-9.9mA, in steps of 0.1mA |
| CutOff Time     | 1-600s, in steps of 1s       |

## Physical:

|                 |                                      |             |
|-----------------|--------------------------------------|-------------|
| Weight          | 2.7Kg                                | (40500-001) |
|                 | 5.3Kg                                | (40552)     |
|                 | 3.4Kg                                | (40553)     |
| Shipping Weight | 4Kg                                  | (40500-001) |
|                 | 9Kg                                  | (40552)     |
|                 | 5.8Kg                                | (40553)     |
| Packing         | 80x60x44cm (Control Unit & one cage) |             |

## Bibliography

### Papers which quote Ugo Basile P.A. Test (step-through)

- C.I. Navarro-Francés et alia: "Influence of trait anxiety on the effects of acute stress on learning and retention of the passive avoidance task in male and female mice" *Behav. Processes* 105: 6-14, **2014**
- L. Zvejniec et alia: "The cognition-enhancing activity of E1R, a novel positive allosteric modulator of sigma-1 receptors" *Br. J. Pharmacol.* 171(3): 761-771, **2014**
- R.W. Flint et alia: "NMDA receptor antagonism with MK-801 impairs consolidation and reconsolidation of passive avoidance conditioning in adolescent rats: Evidence for a state dependent reconsolidation effect" *Neurobiology of Learning and Memory* 101: 114-119, **2013**
- G. Telegdy et alia: "The action of kisspeptin-13 on passive avoidance learning in mice. Involvement of transmitters" *Behav. Brain Res.* 243: 300-305, **2013**
- V. Capurro et alia: "Pharmacological Characterization of Memoquin, a Multi-Target Compound for the Treatment of Alzheimer's Disease" *PLoS ONE* 8(2): e56870, **2013**
- J. Michaud et alia: "Hematopoietic MyD88-adaptor Protein Acts as a Natural Defense Mechanism for Cognitive Deficits in Alzheimer's Disease" *Stem Cell Reviews and Reports* 8 (3): 898-904, **2012**



## Passive Avoidance Step-Down *New Model*

Cat. No. 40570

### General

The **Passive Avoidance step-down cage**, for **mice or immature rats**, is based on the step-down scheme in which the animal is dropped on an elevated platform which becomes uncomfortable because of vibrations.

The instrument basically consists of an **arena**, shaped as a cage (Cat. No. **47573**) and a control unit with touch-screen

The method is based on the mouse tendency to step-down a small platform, uncomfortable because of vibrations, onto the floor of the testing apparatus, which is electrified.

The animal inhibits its behaviour in order to avoid shock; this is measured by longer latency or refusal to step down. Latency is used to assess memory.


Increase or decrease of the **retention latency** gives an indication of improvement or impairment in memory and learning processes.



### Step-Down Cage

### Multiple-Cage Set-Up

Measures the increase/decrease of retention latency to study memory & learning processes

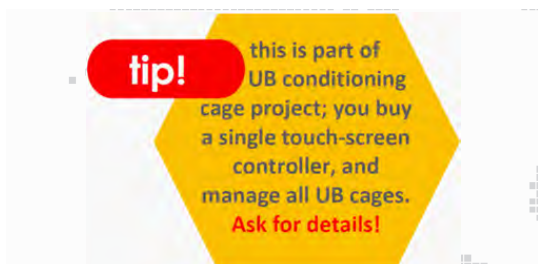

### Main Features

- The electronic unit encompasses all controls for **up to 4 animal cages!**
- Specifically designed for mice or immature rats
- Latency time recorded down to 0.1 seconds

### NEW on the 2014 version!

- The new "launcher" application, makes it possible to manage other UB behavioral cages with the same **Touch-Screen Controller 40500-001**; just purchase the hardware and the application software for the additional test!
- **Remote Control feature** will make remote service and software upgrades extremely simple!



## Instrument Description

Different set-ups, depending on the number of cages, can be obtained by combining the following elements:

- **Programming/Recording Unit with Shocker**
- **Mouse Cage** (up to 4 with one controller)
- **Expansion Box**, for multiple cage set-up

## Programming/Recording Unit

The **40500-001** Programming/Recording Unit, encompassing all controls, incorporates a constant-current high precision 8-pole shocker, and manages data acquisition: data are stored inside the unit and can be downloaded via the USB key provided as standard, for further processing via Excel, Access, etc.



The Unit, incorporating a 12" touch-screen, manages the Passive Avoidance Test via the **40570-010** Software. Up to 4 cages can be connected to the same Controller.

If more than one cage is connected, an expansion box **40500-005** is required for each additional cage.

## Passive Avoidance Cage (step-down)

The cage, dimensioned 28(w)x23(d)x26(h)cm, is provided with a top lid; the cage floor is made of 0.2cm diam. bars, spaced 0.5cm apart, wired to the constant current 8-pole scrambling circuit, located in the control unit.



The detachable circular platform, diam. 7cm, is positioned at the centre of the cage, on a protruding stud fastened to the actuator, the mechanism which energizes the platform vibration.

A larger platform diam. 11cm, is also supplied with the standard package.

## Principle of Operation

When the elevated platform onto which the animal is dropped becomes uncomfortable because of vibrations, the animal steps down to an electrified grid.

When the mouse confronts the electrified grid and re-returns to the platform, the stop command (or pedal switch) is used to halt platform vibration, and stop the latency counter; the touch-screen controller records the latency time in tenths of seconds.

The latency figure remains frozen until a new "session" is started. experimental data are stored inside the controller memory, for further processing.

The vibration intensity is selected from 10 to 100Hz, in 10 steps (10Hz each). The shock intensity can be preset in the range 0 to 3mA, in steps of 0.1mA.

A delay up to 15 seconds can be set in steps of 1s.

## Ordering Information

- 40500-001** Programming/Recording Unit & Shocker
- 40570-010** P.A. Software and activation
- 47573** Mouse Cage, complete 2 platforms
- 40500-005** Expansion Box, for multiple cage set-up

## Specifications

|              |  |
|--------------|--|
| Start        | from the touch screen, or via pedal switch |
| Stop         | from the touch screen, or via pedal switch |
| Vibration    | 10-100Hz, in 10 steps (10Hz each)          |
| Shock        | 0 to 3mA, in 0.1mA steps                   |
| Delay        | 0-15 seconds, in 1s steps.                 |
| Latency Time | 0.1s steps                                 |

## Physical

|                 |                                     |
|-----------------|-------------------------------------|
| Dimensions      | 28(w)x23(d)x26(h)cm (Cage)          |
| Dimensions      | 33(w)x25(d)x5.5(h)cm (Control Unit) |
| Weight          | 8Kg                                 |
| Shipping Weight | 16Kg (approx.)                      |
| Packing         | 80x60x44cm                          |

## Bibliography

### Papers which quote the P.A. Test (step-down)

- A. Mikulecká et alia: "Consequences of early postnatal benzodiazepines exposure in rats. I. Cognitive-like behavior" *Front. Behav. Neuroscience* 8 : 101, 2014
- I.K. Celikyurt et alia: "Effect of harmaline, an endogenous  $\beta$ -carboline, on learning and memory in rats" *Pharmacol. Biochem. & Behavior* 103: 666-671, 2013
- D.S. Dimitrova & D.P. Getova-Spassova: "Effects of Galantamine and Donepezil on Active and Passive Avoidance Tests in Rats With Induced Hypoxia" *J. Pharmacol. Sciences* 101 : 199-204, 2006
- M. Sakaguchi et alia: "Effects of beta-casomorphin-5 on passive avoidance response in mice" *Biosci. Biotechnol. Biochem* 67 (11): 2501-2504, 2003

# Learned Helplessness

Cat. No. 40532 Rats

Cat. No. 40533 Mice

## General

When rodents are exposed to inescapable and unpredictable stress, such as forced swim or inescapable footshock, they often develop deficits in memory and learning tasks (**e.g. Active Avoidance**), and they often show also analgesic reactions (**S.I.A. Stress-Induced Analgesia**).

The **Ugo Basile Set-Up for Learned Helplessness** is based on a sophisticated generator of unpredictable random shocks delivered to the grid floor of a rodent box where no escape is possible.

Electric shocks can be randomized in terms of shock length and interval.

Complex trains can be programmed.

**Up to 4 animals** can be treated simultaneously in 4 independent boxes, controlled by the same electronic unit and software.



The set-up for Learned Helpless is part of the new UB Behavioral Cage program, exploiting the potentiality of a modern controller with touch-screen.



## IDEAL TO STUDY

- Depression & Stress
- Learning & Memory Impairment
- Stress-Induced Analgesia (S.I.A.)

## New Model

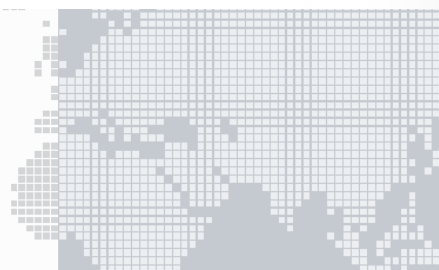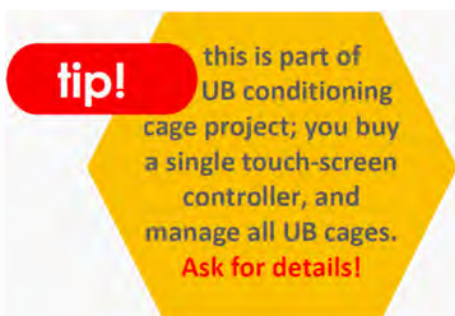

## Main Features

- Randomizable shock patterns
- **Maximum flexibility:** configure your own Experimental Schedules on the touch-screen controller
- The electronic unit encompasses all controls for **up to 4 animal cages!**
- The new "launcher" application, makes it possible to manage other UB behavioral cages with the same **Touch-Screen Controller 40500-001**; just purchase the hardware and the application software for the additional test!
- **Remote Control feature** will make remote service and software upgrades extremely simple!


## System Description

Different set-ups, depending on animal (rat or mouse) and number of cages, can be obtained by combining the following elements:

- **Touch-Screen Controller with Shocker**
- **Rat Cage** (up to 4 with one controller)  
or
- **Mouse Cage** (up to 4 with one controller)
- **Expansion Box**, for multiple cage set-up

## Programming/Recording Unit

The **40500-001** Programming/Recording Unit, encompassing all controls, incorporates a constant-current high precision 8-pole shocker, and manages data acquisition: data are stored inside the unit and can be downloaded via the USB key provided as standard, for further processing via Excel, Access, etc.

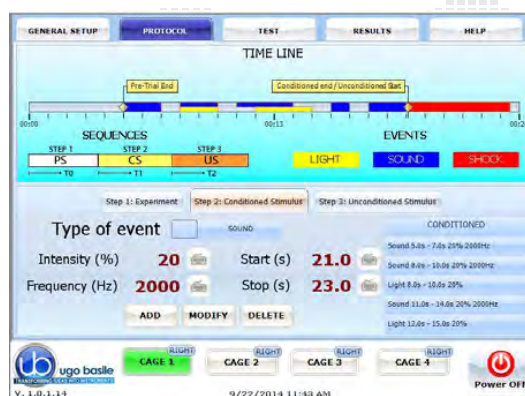

The Unit, incorporating a 12" touch-screen, manages the Learned Helplessness Test via the **40530-010** Software. Up to 4 cages can be connected to the same Controller.

If more than one cage is connected, an expansion box **40500-005** is required for each additional cage.

The trials can be configured on the touch-screen controller, entering the setting via the virtual keyboard: train features, shock and timing of the different experimental sequences.

The system includes a user-friendly reporting software, to collect, visualize and manage data related to the delivered shocks; this is especially important to analyze the randomized shocks and have full control on the performed stimulation.

## Randomizer

The **Touch-Screen controller** is also a sophisticated generator of unpredictable random shocks delivered to the grid floor of the cage.

Electric shocks can be randomized in terms of shock length, interval and complex trains can be programmed. It connects to up to 4 cages.

## Rat and Mouse Cage

The dimensions of **Rat Cage 47502** are 22x22x20(h)cm; **Mouse Cage 47503** is dimensioned 17x17x20 (h) cm.

Both Cages include an electrified floor and a catch pan.

The electrical stimulus is applied to the floor bars of the cage and by an 8-pole "scrambling" circuit incorporated in the touch-screen controller.

All necessary cables and connectors are included to make it a ready-to-use system!

## Ordering Information

- 40500-001** Touch-Screen Controller & Shocker
- 47500-010** Learned-Helplessness Software and activation
- 47502** Rat Cage, complete with electrified floor & catch pan
- 47503** Mouse Cage, complete with electrified floor & catch pan
- 40500-005** Expansion Box, for multiple cage set-up

## Specifications :

Power Requirement 115/230V, 50/60Hz, 30W max.

Shock Parameters : constant current, from 0.1 to 2.9mA in 0.1mA steps

Manual or external operation (via 5V TTL signals), with optional I/O box 46000-150

## Physical:

Weight 2.7Kg (40500-001)  
5.3Kg (47502)  
3.2Kg (47503)

Shipping Weight 4Kg (40500-001)  
9Kg (40552)  
6Kg (40553)

Packing 80x60x44  
(control unit and one cage)

## Bibliography

- Borsini & Cesana: "**Mechanisms of action of flibanserin in the learned helplessness in rats.**" *European Journal of Pharmacology* 433: 81-89, 2001
- Grau et alia: "**Long-term analgesia and activation of the opiate system**" *Science* 213:1409-1411, 1981
- Guilherme dos Santos et alia: "**Antidepressive-like effects of electroacupuncture in rats**" *Physiology & Behavior* 93: 155-159, 2008
- Kademian et alia: "**Biphasic effects of adrenal steroid on learned helplessness behavior by inescapable shock**" *Neuropsychopharmacology* 30: 58-66, 2005



## Conditioned Place Preference Box (CPP)

Cat. No. 42552 for Rat

Cat. No. 42553 for Mouse

### General

The **Ugo Basile Conditioned Place Preference (CPP)** is a 2-compartment box to evaluate the abuse potential of substances and the motivational effects of drugs.

The 2 compartments differ for the wall color and patterns and for the floor patterns and texture.

Both floors and contexts floors are interchangeable so that the visual and tactile difference between the 2 compartments can be easily adjusted by the scientist.

In fact, the CPP box includes the contextual cues required by the experimental paradigm; each box includes:

- 4 interchangeable floors with square and circular patterns
- 3 sets of walls.

The new CPP box has been designed and optimized for visual scoring, or for use with any video-tracking software. See [www.ub.anymaze.com](http://www.ub.anymaze.com).



IDEAL TO STUDY

**Drug Abuse**

**Addiction**

- **Interchangeable floors for tactile stimulation**
- **NEW MODEL with interchangeable CONTEXTS**

### Main Features

- Optimized for Video-Tracking
- Specific models for rats or mice
- Designed for multiple-cage systems
- Interchangeable floors provided for different patterns & texture
- Walls in either compartment can be visually altered, by replacing the context kit



## Rat and Mouse Box


The box **42552** is designed for tests on rats. Its external dimensions are **60x30x30(h)cm**; the box **42553** is similar to the 42502, but its dimensions (**32x15x25(h)cm**) make it suitable for use with mice.

Both boxes have a patterned door in the central wall, 7.5x7.5cm in the rat, 4x6(h)cm in the mouse box.

## Tactile Stimulation: Patterned Floors


One of the major keys to the success of a **CPP** experiment is due to the design of the visual and tactile differences between the 2 compartments.

Ideally the 2 compartments should have clearly distinct contextual cues but should not determine any preference in unconditioned animals.

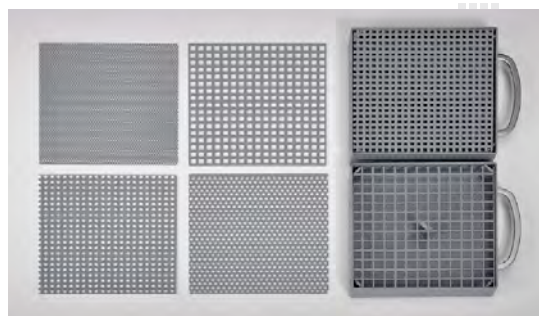

Given the importance of **paw tactile sensitivity** in rodents, while the design of commercially available CPP boxes has traditionally focused only on the wall patterns and colors, the Ugo Basile CPP box includes 4 interchangeable floors with different patterns & texture.

Four sets of floor grids, and 3 sets of replaceable wall contexts (striped, checked, and dotted) are supplied with each box:

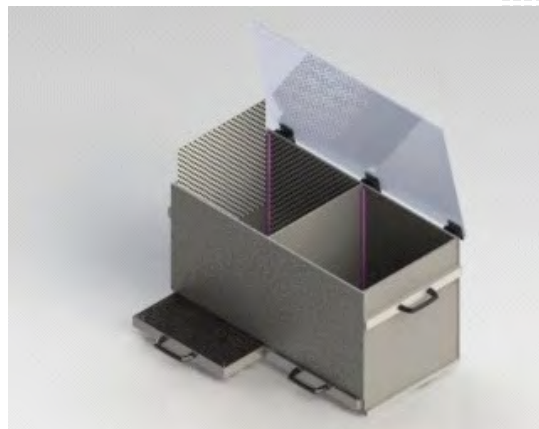

## Rationale and outline of the procedure

The CPP paradigm provides information on the rewarding or aversive effects of visible and tactile contextual cues associated with drugs.

This technique has acquired great popularity in research studies involving addiction, being much easier, if compared to drug self-administration procedures.

First, the animal is conditioned to identify one of the two compartments with the drug experience. Then the time spent in each compartments is measured; preference or aversion to the drug-paired compartment, hence rewarding/aversive properties of drugs, can be easily deduced.

The CPP test only requires the animal to carry out a simple operation (i.e. move from one compartment to the other) to approach or avoid the drug-paired compartment; the animal is expected to spend more time in the drug-paired compartment, if the drug experience produced a positive effect.

## Optimized For Video-Tracking



All floors are grey-colored, to optimize contrast and facilitate tracking of both dark and albino animals.

## Ordering Information

|                   |  |
|-------------------|--|
| <b>42502</b>      | <b>CPP BOX for RAT, including</b>            |
| <b>M-TR 230-F</b> | Floor Drawer (2 pcs.)                        |
| <b>42502-011</b>  | Round 2mm holes, 6mm interax. (2 pcs.)       |
| <b>42502-012</b>  | Round 12mm holes, 16mm interax. (2 pcs.)     |
| <b>42502-014</b>  | Square 6x6mm holes, 9mm interax. (2 pcs.)    |
| <b>42502-013</b>  | Square 10x10mm holes, 12mm interax. (2 pcs.) |
| <b>42552-320</b>  | <b>Wall Context Kit for Rat Cage</b>         |
| Weight            | 22Kg net, 25Kg gross; Packing: 80x60x44cm    |
| <b>42503</b>      | <b>CPP BOX for Mouse, including:</b>         |
| <b>M-TR 230-F</b> | Floor Drawer (2)                             |
| <b>42503-012</b>  | Round 2mm holes, 3mm interax., 2 pcs.        |
| <b>42503-011</b>  | Round 4mm holes, 6mm interax., 2 pcs.        |
| <b>42503-013</b>  | Square 4x4 holes, 7mm interax., 2 pcs.       |
| <b>42503-014</b>  | Square 6x6 holes, 9mm interax., 2 pcs.       |
| <b>42553-320</b>  | <b>Wall Context Kit for Mouse Cage</b>       |
| Weight            | 8Kg net, 10Kg gross; Packing: 36x55x45cm     |

## Acknowledgements & Bibliography

A special thank to Prof. Paola Fadda (Department of Pharmacology, University of Cagliari, Italy) for the initial design of the boxes: her valuable comments and suggestions allowed us to keep the focus on the user needs and opinions.

- L. Fattore et alia: "Baclofen Prevents Drug-Induced Reinstatement of Extinguished Nicotine-Seeking Behaviour and Nicotine Place Preference in Rodents" *Eur. European Neuropsychopharmacol.* (in press 2009)
- M. Scherma et alia: "Inhibition of Anandamide Hydrolysis by Cyclohexyl Carbamic Acid 3'-Carbamoyl-3-yl Ester (URB597) Reverses Abuse-Related Behavioral and Neurochemical Effects of Nicotine in Rats" *J. Pharmacol. and Exper. Therap.* 327:482-490, 2008

## Lickometer - Vogel Test

Cat. No. 45100 Set-up for Rat

Cat. No. 45150 Set-up for Mouse

### General

The **Ugo Basile Lickometer - Vogel Test** is a versatile system that can function as a simple software-based lickometer or as a Drinking-Conflict set-up to assess the anxiolytic effect of drugs.


In the Drinking-Conflict Vogel paradigm, a water deprived animal is exposed to a lickometer and the licking events are coupled to electric shocks.

The animal is in a motivationally conflicting situation, hence his licking behavior is affected by anxiety and anxiolytic drugs.

The Lickometer controller and software can manage up to 5 animal cages for either rat or mouse; one shocker is required for each cage.

The friendly-to-user software, provided as standard, manages the system and experimental configuration, collects and saves the experimental data, and provides a detailed report.

Data are saved as .csv file and .rpt file (a proprietary format which can be opened only within the Lickometer software)

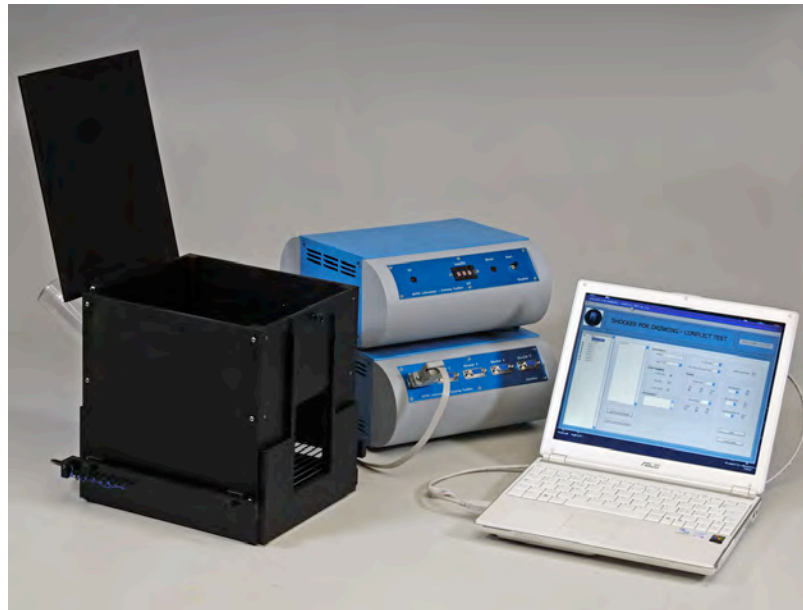

### Specific Models

for Rat

for Mouse

- Vogel Conflict Test
- Lickometer
- Anxiety Testing
- Multiple Chambers

### Main Features

- Up to 5 animal chambers with grid floor, lick sensor, water reservoir
- Software for experiment configuration and data collection
- Two-pole shockers with adjustable shock intensity
- Chambers can be used as a general lickometer

## Rationale of the Test

The Drinking Conflict Vogel test usually consists of three phases:

- Initial wait (triggered by the first licking event)
- Shock phase (the sipper is electrified)
- No-shock phase (no shock is associated to drinking)

For each phase of the experiment, the number and the timing of licking events is recorded and graphically displayed.

The alternation between shock and no-shock phases can be based on TIME or N° OF LICKS, according to the user experimental paradigm.

When no shock is delivered, the system can be simply used as a software-driven lickometer.

The duration of each phase is user-defined for each cage, based either on time or on the animal behaviour (i.e. the sipper is electrified after a defined number of licking events have occurred).

At the end of the test a report will summarize the results; these results can be automatically printed and exported into a spreadsheet.

## System Components

The system is composed of:

- USB-Control Unit for up to 5 cages, including
  - Software
- Drinking Conflict Cage
- 2-Pole Sine-Wave Shocker

## Animal Cages

Drinking-conflict cages are provided with grid floor, electrified sipper and lick sensor. Two sizes are available, for rats and mice.

The rat cage inside dimensions are 35(w)x25(d)x30(h)cm.

The mouse cage is dimensioned 20(w)x24(d)x20(h)cm.



## Lickometer Software


The following parameters, which define the experimental configuration can be set via the software for each cage:

- Trial duration
- Initial Pause
- Time Intervals with and without shock
- Number of licks to deliver a shock etc.



*Experiment configuration*

For each cage, it is possible to assign a specific name to report, operator and animals involved in the experiment; sex and weight of the animals can also be specified.

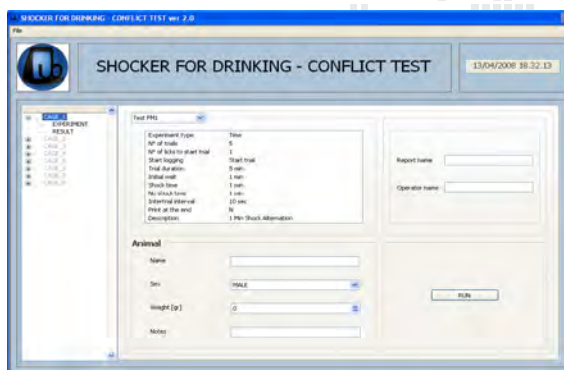

*Cage configuration*

The software collects the experimental data and saves them as .csv file & .rpt file (the latter a proprietary format which can be opened only within the Lickometer software). A complete report file is provided at the end of the experiment; results can be automatically printed and exported into a datasheet.

## Ordering Information

**45100 Lickometer Set-up for RAT, one cage, including:**

- 45100-002 Rat Cage**
- 45100-001** 5-channel Electronic Unit
- 45100-005** Software
- 45100-004** Shocker
- 45100-302** Instruction Manual

**45150 Lickometer Set-up for MOUSE, one cage:**

- 45100-003 Mouse Cage**
- and other components as for 45100**

| Physical        | 45100      | 45150      |
|-----------------|------------|------------|
| Weight          | 8.5Kg      | 7.5Kg      |
| Packing         | 80x60x44cm | 80x60x44cm |
| Shipping Weight | 12Kg       | 10Kg       |

## Bibliography

- P. Ohara et alia: "Evidence for a Role of Connexin 43 in Trigeminal Pain Using RNA Interference In Vivo" J. Neurophysiol 100: 3064-3073, **2008**
- J.P. Vit et alia: "Silencing the Kir4.1 Potassium Channel Subunit in Satellite Glial Cells of the Rat Trigeminal Ganglion Results in Pain-Like Behavior in the Absence of Nerve Injury" J. Neurosci. 28(16): 4161-4171, **2008**

## Sociability Apparatus ( 3-chambered social test )

Cat. No. 46503

### FOR STUDIES ON:


- Autism
- Social Memory & Novelty
- Pair-bonding
- Dominance hierarchies

### General

The new Ugo Basile Sociability Apparatus (3-chambered social test) is a valuable tool to study SOCIAL INTERACTION in mice. It consists of a transparent Perspex cage, a special non reflective, grey colored floor and 2 grid enclosures.

Many authors (e.g. Moy et al. 2004; Nadler et al. 2004) have shown that a 3-chambered box can be used to test:

- Social Novelty Preference
- Sociability
- Dominance

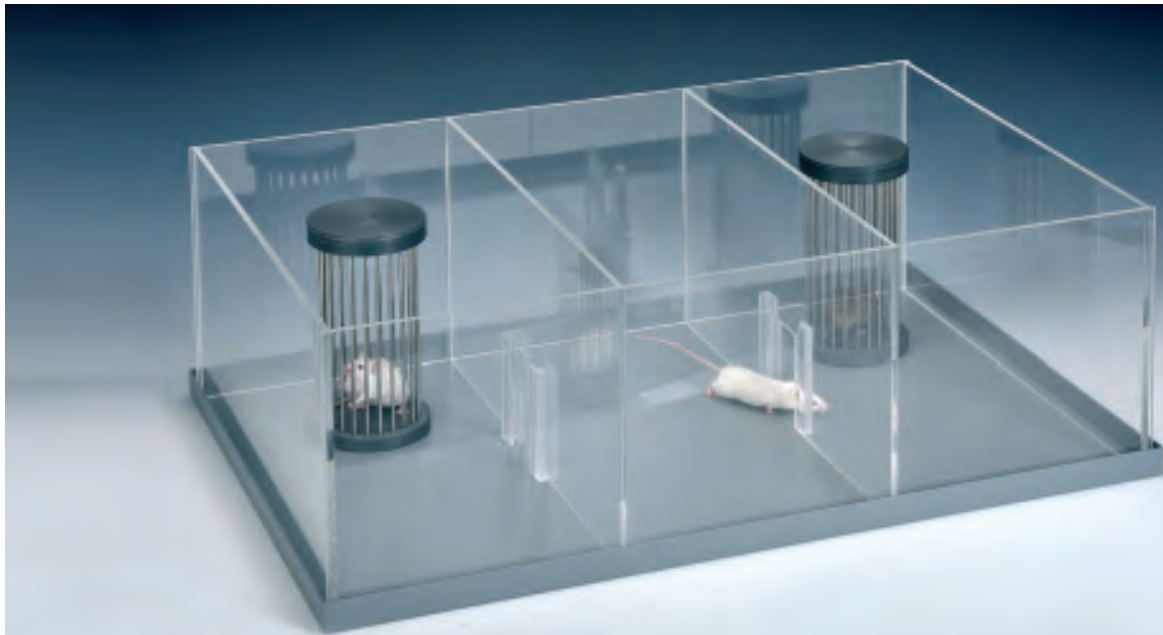

### Main Features

- Works even with the most basic video-tracking software
- Grid Enclosures maximize animals interaction
- Different colours and sizes are available on request



## Rationale and Outline of the Procedure

The Ugo Basile 3-Chambered Apparatus can be used with many different procedures.

In their 2004 paper, Moy and co-authors (see bibliography), describe a typical procedure: after a period of habituation a mouse sociability is determined by measuring the time spent by the freely-moving mouse in the proximity of the grid enclosures containing the first 'stranger' mouse.

A second 'stranger' mouse is then introduced in the box and the preference for the new 'stranger' mouse can be easily assessed.



## 3-Chamber Box & Grid Enclosures

The clear Perspex box gives ideal transparency for visual observation of the experiment or for side positioning of the video-camera.

Dimensions for each of the 3 compartments are 20x40x22(h)cm; two sliding doors (5x8(h)cm), opening on the central compartment, can be closed to confine the animal.

The grey floor gives high contrast with both light and dark animals, allowing for automated video-tracking of the animals. Its special painting also gives a slightly rough surface, pleasant for the animals to walk on.

The grid enclosures allow mice to interact closely; the grid bars have a diameter of 3mm and are spaced 7mm.




The standard enclosures are 15cm tall and have an internal diameter of 7cm. The top and the bottom are made of grey (46503-003) or white (46503-013) PVC.



Different sizes/colors are available on request.

## Optimized for Video-Tracking

The grey floor gives best contrast to both light and dark animals, which is the most critical factor for **all video-tracking softwares** to work properly.

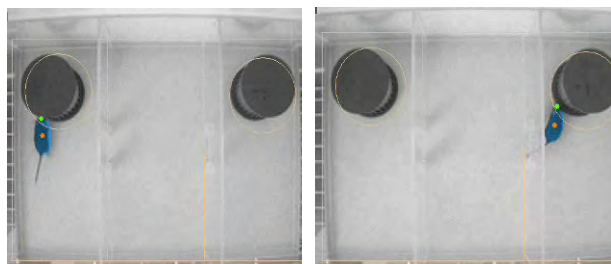

Images and videos, courtesy of Dr. Patrizia D'Adamo (San Raffaele Institute, Milan, Italy)

## Ordering Information

- 46503 Mouse Cage for Sociability**, complete with following components:
- 46503-003** Grid Enclosure for Sociability ("stranger" cage) **grey**, I.D. 7cm, 15cm h (2 pieces)
  - 46503-010** 3-Chambered Transparent Cage, with two doors
  - M-LM 570** Metal Floor, special non-reflective color
- 46513 Mouse Cage for Sociability**, complete with following components:
- 46503-013** Grid Enclosure for Sociability ("stranger" cage) **white**, I.D. 7cm, 15cm h (2 pieces) and other components as for **46503**

## Physical

|                 |               |
|-----------------|---------------|
| Dimensions      | 60x40x22(h)cm |
| Weight          | 9Kg           |
| Shipping Weight | 12Kg          |
| Packing         | 67x42x53cm    |

## Bibliography

- M. J. Kane et alia: "Mice Genetically Depleted of Brain Serotonin Display Social Impairments, Communication Deficits and Repetitive Behaviors: Possible Relevance to Autism" *PLoS ONE* 7(11): e48975, 2012
- M. Yang et alia: "UNIT 8.26 Automated Three-Chambered Social Approach Task for Mice" *Current Protocols in Neuroscience* Published Online: 1 July 2011

## Method Papers

- S.S. Moy et alia: "Sociability and Preference for Social Novelty in Five Inbred Strains: an Approach to Assess Autistic-Like Behavior in Mice" *Genes, Brain and Behavior* 3(5):287-302, 2004
- J.J. Nadler et alia: "Automated Apparatus for Quantitation of Social Approach Behaviors in Mice". *Genes, Brain and Behavior* 3(5): 303-314, 2004.

## HYDRAULIC “ATLANTIS” PLATFORMS

for WATER MAZE experiments

Cat. No. 40100

LIFTING CONTROL

LOWERING CONTROL

NO ELECTRICITY

NO HANDS IN  
THE POOL !

### Why Automated Platforms?

Despite being very effective, the **Morris Water Maze** task has some limitations, related to the platforms normally used having fixed height, which cannot be raised during probe tests. Probe tests run with the use of a **lift platform** give more reliable indications on the presence of true **spatial learning**.

The Ugo Basile Atlantis Platforms are made of clear Perspex and are operated by hydraulic pressure. No electricity is present inside the pool; the electrical parts of the mechanism (i.e. the electro-hydraulic actuators) are safely located outside.



### Main Features

- 4-Platforms with one Controller
- Remote lifting/lowering control
- Manually or PC-Operated
- Consistency of positioning in the 4 quadrants
- No more hands in the pool!
- No Electricity in the pool

## System Description

Up to 4 platform/motor combination connect to the 4-channel control unit.

Each platform is **driven independently**, so that the Water Maze experiment can be completely automated by positioning a platform in each of the 4 quadrants of the pool.

Once the 4 platforms have been positioned in the pool, each is connected to the related external motor, via the connectors conveniently fitted to the water tank ([ask for information about our models](#)); the whole experiment can then be run automatically, via the control unit or external triggers.

## Specifications

- 4 independent channels : manual or TTL mode
- Platform vertical range : 25-35cm
- Vertical travel : 10cm, in 9 steps
- Platform Speed : 10mm/s
- Platform diameter : 10cm

## Manual or Automated Modes

The platforms go up and down in steps of 1 cm, for a total vertical travel of 10 cm.

Different operation modes are possible using Ugo Basile Atlantis platform system: in the **manual** mode the vertical travel is controlled by simply depressing a key.

In the **automated mode** the platforms can be operated by external triggers (TTLs), controlled by any videotracking software.



Each platform can be kept submerged, and raised automatically when the animal swims above it. This protocol allows one to exclude from the test "navigation strategies" in which spatial memory is not involved.

platform up →

↓ platform down



When used as stand-alone tool, without motor/controller, the Atlantis hydraulic **platform 40101-002** can also conveniently replace standard fixed platforms.

## Ordering Information

- 40100** Complete 1-Platform System, including standard components as listed below
- 40400** Complete 1-Platform System, including standard components as listed below

|                  |                      | 40100 | 40400 |
|------------------|----------------------|-------|-------|
| <b>40100-001</b> | 4-Channel Controller | 1     | 1     |
| <b>40101-002</b> | Platform             | 1     | 4     |
| <b>40101-003</b> | Motor                | 1     | 4     |
| <b>40101-320</b> | Connection Cable     | 1     | 4     |
| <b>40101-321</b> | 100ml Syringe        | 1     | 4     |
| <b>40101-322</b> | Stretch of Tube (3m) | 1     | 4     |
| <b>40100-302</b> | Instruction Manual   | 1     | 1     |
| <b>E-WP 008</b>  | Mains Cable          | 1     | 1     |

## Available Accessories

- 40101** Additional platform and motor assembly

[Ask for information about our Water Mazes and ANY-maze videotracking software](#)

| Physical        |    | 40100    | 40400 |
|-----------------|----|----------|-------|
| Weight          | Kg | 11       | 30    |
| Shipping Weight | Kg | 17       | 39    |
| Packing         | cm | 80x60x44 | (x2)  |

## Bibliography

- R.I.W. Spooner et al.: "The Atlantis Platform: A New Design and Further Developments of Buresova's On-demand Platform for the Water Maze" [Learn. Mem.](#) 1: 203-211, 1994
- G. Riedel et al.: "Reversible Neural Inactivation Reveals Hippocampal Participation in Several Memory processes" [Nature Neurosc.](#) 2 (10): 898-905, 1999
- I.Q. Wihshaw et al.: "The Behavior of the Laboratory Rat: A Handbook with Tests" [Oxford Univ. Press, USA](#): 1, 2004



# ANY-maze

## Advanced Videotracking

Cat. No. 60000


### General

ANY-maze is a flexible video tracking system designed to automate testing in behavioural experiments.

Packed with advanced features ANY-maze is one of the most comprehensive video tracking systems available today

### Flexibility

With a single ANY-maze system you can easily automate a range of apparatus, for example, a plus maze, a water maze and a set of 6 locomotor activity boxes.

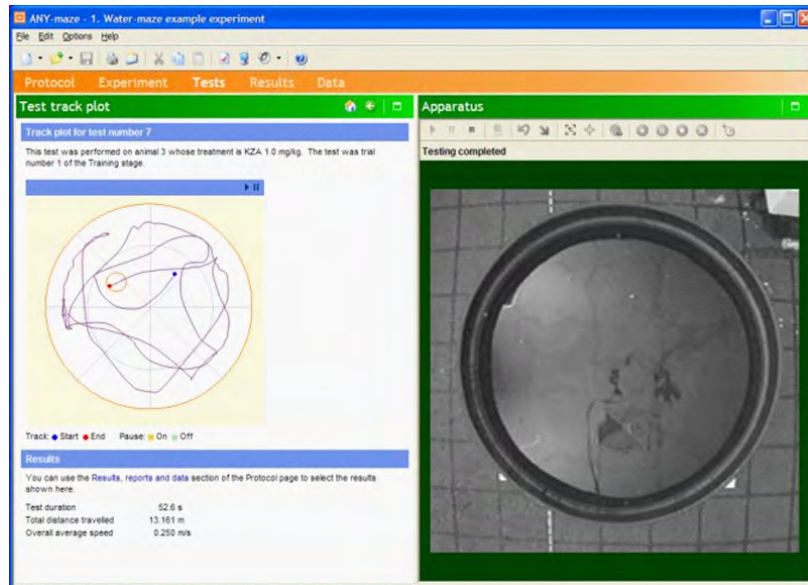

**User-friendly interface  
and flexibility**

**Compatible with most cameras  
and digitizers**



### Take a tour and see for yourself

The quickest way to learn more about ANY-maze is to take a brief introductory tour...



### Download ANY-maze and try it out

Why not try ANY-maze for yourself - you can download the complete system for free!

We've even included some experiment videos so you can see the tracking in action.

### Video tracking your animals in a wide range of behavioural apparatus:

- Morris Water-Mazes
- Elevated Plus Mazes
- O-T- Y-Mazes
- Radial Mazes
- Open Fields
- Home Cages
- Metabolic Cages
- Place Preference Boxes
- Porsolt Forced Swim Tests
- Tail Suspension Tests

## Equipment

ANY-maze's flexible design makes it easy to set up experiments in a wide range of different apparatus: plus maze, water-maze, T-maze, activity boxes, forced swim test, open-field cages, Fear Conditioning, etc.


Our extensive range of **high-quality mazes** & test apparatus have been designed in cooperation with experienced behavioural scientists and are optimized for video tracking, include

What's more, all our mazes and test apparatus can be automated using the standard version of ANY-maze; so you only need a **single piece of software to automate any of these tests**; additionally, some devices, such as our **Fear Conditioning system** are available with low-cost versions of ANY-maze specific to the device.

## Simultaneous Testing

Using ANY-maze you can perform tests in **up to sixteen pieces of apparatus simultaneously**. This provides a great way to increase throughput and also makes it easier to control for environmental variables.

And ANY-maze's versatile camera management means you can use one camera, or many, to view the apparatus. For example, in these place preference boxes four cameras are being used, one on either side of each box.

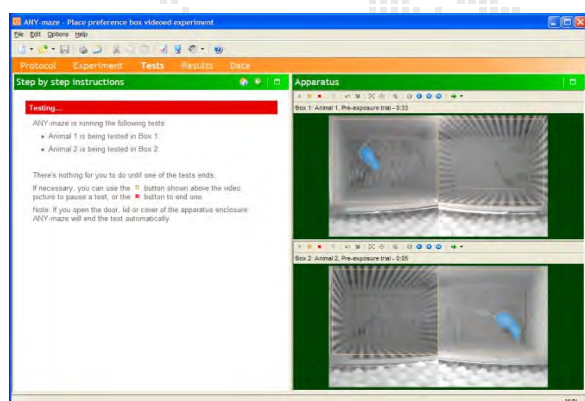

## Cameras & Computers

With such flexibility, how do you determine the computer, cameras etc., that you'll need?

The answer's provided by the **ANY-maze equipment wizard** which quizzes you about all the apparatus you want to automate and then creates a detailed report of the equipment required.

You can use ANY-maze with inexpensive USB web-cams, high quality 'machine-vision' USB cameras, DV camcorders or almost any analogue CCTV camera.

This breadth of support not only makes it easy to find a compatible device (indeed, you may already own one), but also means that the system can meet a range of differing requirements, such as low cost, notebook connectivity, simultaneous capture from multiple cameras, tracking in darkness, etc.



ANY-maze has been designed to work with modern computers running Windows Vista, Windows XP, Windows 7.

However, that doesn't mean you can't use it with older equipment or other versions of Windows check **computer compatibility** on our web site.

## Ordering Information

|                  |  |
|------------------|--|
| <b>60000</b>     | <b>ANY-MAZE LICENSE, including technical support and updates for 1 year</b>  |
| <b>60050</b>     | ANY-Maze, 1 year support extension (*)   |
| <b>60000-FC</b>  | ANY-Maze, Freezing detection only, for Fear Conditioning   |
| <b>47400-030</b> | USB Camera, with 2.1 & 4.3 lenses, visible block filter, cables, and ceiling support   |
| <b>47400-010</b> | Black and white high sensitivity videocamera including varifocal day & night lens & ceiling support  |
| <b>47400-011</b> | Analogue-Digital Converter PCI RTV24, 4 channels, for connection of 47400-010 to desktop PC. Complete kit including cables                               |
| <b>47400-012</b> | Analogue-Digital Converter FireWire, 1 channel, for connection of 47400-010 to laptop PC. Complete kit including FireWire adaptor, cables & power supply |

## ANY-maze License

### How ANY-maze licensing works

- You can download ANY-maze from this site for free and install it on any number of computers.
- You only need a license for copies which will be used for tracking - you can use other, free copies to set up experiments, analyse results, transfer data etc.
- To license a copy of ANY-maze, so you can use it for tracking, you supply us with its serial number and pay the purchase price. We then supply you with your license number which will **permanently** enable the tracking system and will permit updates to be installed for a period of 1 year.

### What's included in the price

- The ANY-maze software itself.
- All updates to the system for a period of 1 year.\*
- Technical support for a period of 1 year.\*

### \* Extended Support

- When you purchase ANY-maze, we supply technical support and all upgrades for a period of 1 year.
- To get support and upgrades after this period you need to purchase an extended support contract.
- When you purchase an extended support contract we will supply you with a new ANY-maze license number, this will permit updates to be installed for a further period of one year.

# Animal Mazes for Video-Tracking

## FOR STUDIES OF:

- Anxiety and Stress
- Memory and Learning
- Spatial Memory
- Activity and Exploration

## General

The Animal Mazes manufactured by Ugo Basile are designed to give optimal results with any Video-Tracking software. This is achieved by:

- *high-contrast colors*: grey, white, black or the NEW Ugo Basile Light-Blue
- *non-reflective colors*: reflections are a common source of error in animal tracking. Let's avoid them!

All maze materials were selected to be *sturdy and easy to clean*, to construct reliable and durable mazes.



## Main Features

- High-contrast, non-reflective colors optimized for Video-Tracking
- Quality materials: light, easy to clean and to store
- Surface texture selected for best rodent's comfort (reasonable rough, "warm" surface)



## Water Maze Pool

The Ugo Basile Water Mazes are water pools specifically manufactured for Morris Water Maze experiments (i.e., not a cattle drinking trough) and include:

- wheels and drain hose

- built-in connectors for Hydraulic Atlantis Platforms (not included)

- customizable colors and dimension on request

- animal platform (fixed height, 10 or 12 cm diameter)

Pools are 60 cm high and 120, 150 or 180 cm diameter.



## Barnes Maze

- Mouse version: 100 cm diameter, 5 cm hole diameter

- Rat version: 130 cm diameter, 10 cm hole diameter



Both versions are 60 cm high and are painted in non-reflective grey or light-blue (white, black or other custom colors are available on request). The animal shelter is included and is magnetically attached to the maze, for quick and easy experiments.

## Elevated Plus-Maze and Zero-Maze



Elevated Plus-Maze



Zero - Maze

These mazes are manufactured from high-tech metal alloy and can be painted in different colors. Dimension (cm):

- Elevated Plus-Maze, Mouse: arm length 35, arm width 5, closed wall height 15, height from the floor 60

- Elevated Plus-Maze, Rat: arm length 50, arm width 10, closed wall height 40, height from the floor 60

- Zero-Maze, Mouse: diameter 55, corridor width 5, wall

height 15, height from the floor 60 cm

## Y-maze, T-maze



These mazes have a metal base painted in non-reflective grey (more colors on request) and plastic arms that can be disassembled and closed with the included doors. Dimension (cm):

Y-maze, Mouse: arms length 35, width 5, wall height 10

Y-maze, Rat: arms length 50, width 10, wall height 20

T-maze, Mouse: stem length 35, arm length 30, width 5, wall height 10

T-maze, Rat: stem length 50, arm length 40, width 10, wall height 20

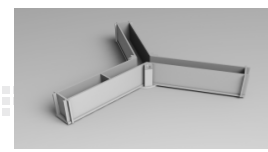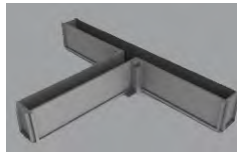

## Open-Field

Open Fields are available in non-reflective grey color, for mice (44 cm) or for rats (100 cm); both versions have detachable walls for ease of storage.

## NEW Mouse Radial Maze

The new Mouse Radial Maze is manufactured from high-tech metal alloy and durable plastics to be as sturdy as possible.

A new automated model, with retractable doors is also available.

Different colors are available, all non-reflective, and arms can be detached, for easy cleaning.

Dimension (cm):

arms length 35, width 5, wall height 10



## Ordering Information

- 40125 Water Maze, 120 cm, for mice
- 40155 Water Maze, 150 cm, for mice and rats
- 40185 Water Maze, 180 cm, for rats
- 40193 Barnes Maze, for mice
- 40192 Barnes Maze, for rats
- 40142 Elevated Plus-Maze, for rats
- 40143 Elevated Plus-Maze, for mice
- 40163 Elevated Zero-Maze, for mice
- 40173 Y-maze, for mice
- 40172 Y-maze, for rats
- 40133 T-maze, for mice
- 40132 T-maze, for rats
- 47432 Open-field, 44 cm, dark walls
- 47433 Open-field, 44 cm, transparent walls
- 47100 Open-field, 100 cm, dark walls
- 47150 Open-field, 100 cm, with 4 partitions

## MULTI-MAZE SYSTEM for Mouse

Cat. No. 41500

**FULLY CONFIGURABLE MULTI-MAZE  
(Y, T, RADIAL MAZE)**

**OPTIMIZED FOR VIDEO-TRACKING  
AND OPTOGENETICS**

The new MOUSE MULTI-MAZE Cat. No. 41500, will help the researcher to conduct fully automated memory experiments such as:

- Assessing spatial memory
- Testing basic working memory
- Discriminating working from reference memory
- Evaluating impairments in the working memory

The electronic unit features USB interface, 8 independent TTL inputs and integration with video-tracking software.

The proprietary sliding doors retract in the maze floor, ensuring unobstructed animal tracking, while guaranteeing smooth, silent, totally automated up and down movements.

All the animal mazes manufactured by Ugo Basile are designed with high-contrast colors and non-reflective coatings, providing optimal results with any video-tracking software.

Ugo Basile mazes are constructed of sturdy and easy to clean materials, making them the most reliable mazes on the market.

Surface texture was selected for best rodent's comfort.

The proprietary doors sliding underneath the floor make it the best solution for videotracking and optogenetics test.



### Main Features

- New proprietary modular system
- Doors slide underneath the floor
- Smooth and silent operation
- Easily cleanable
- Manual or PC-driven operation modes (via TTL or USB connection)
- Interchangeable walls for egocentric or allocentric spatial memory tests, different colors or textures available on request

## System Description

The new **MOUSE MULTI-MAZE** Cat. No. **41500** is a modular system, enabling the user to set-up an electronically controlled maze, by combining one of the different arenas provided, and the required number of arms, in the following configurations:

- **Mouse Y-Maze**
- **Mouse T-Maze**
- **Mouse 8-Arm Radial Maze**

This feature is facilitated by the new concept introduced with this recent version: the door-controlling kinematics of each section is actually integral part of the arm itself, being positioned below the door area, while a control unit, positioned below the central arena, consolidates the motor control board, the interface with the external electronic unit, and the interface with the video-tracking software (ANY-maze, not included).

Each arm may be limited in length, by the use of spacers (provided). The corridor side walls, made of plastic material, are easily removable, for cleaning purposes. Moreover, it will be easy to switch from high profile to low profile wall (optional), according to the research needs.

## System Configurations

### 41503 Y-Maze for Mouse

- 1 **41500-001** Central Control Unit
- 3 **41500-002** Standard Arm, with automated door
- 1 **41153-010** Electronic Unit (8 TTL outputs)



**Arm dimensions:** Length 35cm  
Width 5cm  
Height 12cm

An automated door is provided on each arm, at the central arena end.

### 41504 T-Maze for Mouse

- 1 **41500-001** Central Control Unit
- 3 **41500-002** Standard Arm, with automated door
- 1 **41500-003** "Start" compartment with automated door
- 1 **41153-010** Electronic Unit (8 TTL outputs)

**Arm dimensions:** Length 35cm  
Width 5cm  
Height 12cm



An automated door is provided on each arm, at the central arena end; the additional "start" compartment, with automated door, attached to the end of the stem-arm, completes the T-Maze configuration.

### 41508 8-Arm Radial Maze for Mouse

- 1 **41500-001** Central Control Unit
- 8 **41500-002** Standard Arm, with automated door
- 1 **41153-010** Electronic Unit (8 TTL outputs)

**Arm dimensions:** Length 35cm  
Width 5cm  
Height 12cm

An automated door is provided on each arm.

*The 41504 and 41508 configurations also enable the Y-maze test to be carried out, without any extra accessories.*

## Ordering Information

### Components

- 41500-001** **Central Control Unit**, incorporating motor controls and interface to external unit, provided with 3 central arenas (Y, T, radial)
- 41500-002** **Standard Arm**, provided with automated door, and high profile (12cm h)
- 41500-003** **"Start" Compartment** for T-maze, with automated door Arm
- 41153-010** **Electronic Unit** (8 TTL outputs)
- 41500-310** **ANY-maze interface board**

### Configurations

- 41503** **3-Arm configuration**, for Y-maze test, standard profile walls, automated doors
- 41504** **3-Arm configuration**, and "Start" Compartment, for Y-maze and T-maze test, standard profile walls, automated doors
- 41508** **8-arm Radial Maze configuration**, for Radial and Y-maze test, standard profile walls, automated doors



## Isolated Organ Baths

Cat. No. 4000 / 4050 / 4400

### General

The Isolated Organ Baths have been designed for accurate recording of isometric or isotonic tissue contraction/release.


Research involving effects of electrical stimuli or drugs on isolated organs, uterus, trachea, vessel strips, auricle, can be performed under optimum conditions.

Wash or test solution enters the chamber after passing through the temperature equilibrating coils and the syringe valve. The tissue in the chamber is washed by flushing the chamber through an overflow drain tube, to avoid exposing the tissue to the air.

Water stirring is accomplished by a water jet delivered by a centrifugal pump.

A 200W stainless steel heating element is mounted on the Perspex tank floor. A precise solid state "proportional" thermostat maintains the temperature within the excellent limits of  $\pm 0.1^{\circ}\text{C}$  on all models.

*Note : the Isometric and Isotonic Transducers featured in the picture are not included and can be ordered separately.*

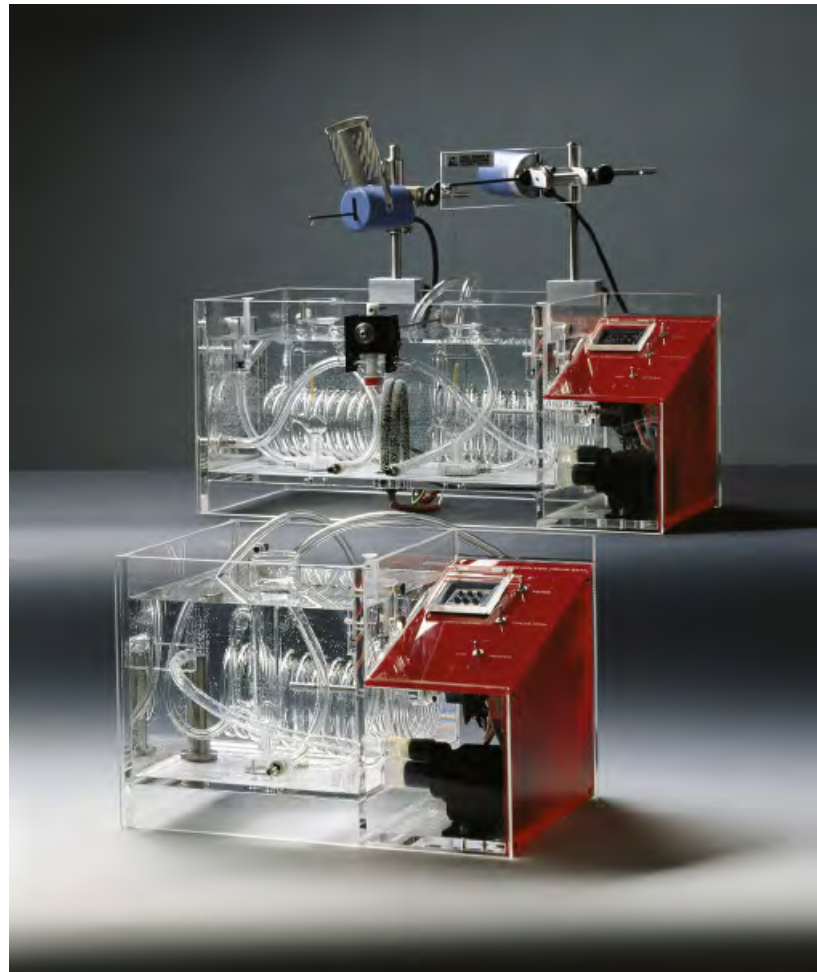

**SAFETY**

**EFFICIENCY**

4000 One Muscle Chamber

4050 Two Muscle Chambers

4400 Four Muscle Chambers

### Main Features

- All components visible through the clear Perspex tank: great for teaching!
- Tissue washing without exposure to air
- Water-jet bath stirring provided by a noiseless vibration-free centrifugal pump
- Easy and quick mounting of tissue

## Bath 4000

The 4000 water bath consists of a clear Perspex tank, 19x19x17cm which contains one tissue chamber, one temperature equilibrating coil, one adjustable support rod on which transducers can be fastened to the tank via the holder provided.

## Bath 4050

This is similar to the one-chamber bath 4000 but the tank is dimensioned (34x19x17cm) to accommodate two muscle chambers and syringe valves, two coils, two adjustable support rods and holders for transducers.

## Bath 4400

The bath 4400 lodges up to 4 preparations; they maintain the features of the 4050 but heating power and dimensions are upgraded accordingly, the tank being 47x29x22cm.

## Tissue Chamber Configuration

The tissue chambers provided with porous frit, available in 5, 10, 20, 30 or 50 ml are standard. Unless otherwise specified, we supply our tissue baths with 10ml muscle chambers.

An accurately positioned glass hook is provided in the chamber to which the thread loop can be easily attached, ensuring the organ being well centered in the chamber.

Tissue chambers are also available provided with an aeration side arm in 20, 30 or 50 ml volume. Tissue chambers without hook are available on request.

## Control Box

The control section of the bath lodges the electronics; the temperature regulator, the temperature sensor & the circulation motor are connected to by connectors enabling quick disconnection for servicing purposes.

The upper panel consolidates all commands and the temperature regulator, with keys to preset water temperature in the range 25-45°C, enabling an accurate temperature setting in 0.1°C steps.



## Recording & Transducers

Ugo Basile offers a complete line of Transducers (Iso-metric 7003, 7004, 7005, 7010 or Isotonic 7006) and a versatile 4-channel digital recorder, the 17304 Data-Capsule-Evo. Ask for details!

## Ordering Information

**4000** **Isolated Organ Bath, One Muscle Chamber**, with circulation pump, heater, thermostat, temperature sensor, complete provided with following standard accessories:

**4005** Temperature Equilibrating Coil

**4100** Muscle Chamber, 10ml, provided with porous frit and hook

**14110** Lead-Screw Positioner for 10 & 13mm rods

**4004** Supporting Rod (10mm diam.)

**4000-302** Instruction Manual

**E-WP 008** Mains Cord

**4050** **Isolated Organ Bath, 2 Muscle Chambers**, as above but all standard accessories multiplied by two, i.e., 2x4005, 2x4100, etc.

**4400** **Isolated Organ Bath, 4 Muscle Chambers**, as above but all standard accessories multiplied by four, i.e., 4x4005, 4x4100, etc.

### Physical:

**4000** Dimensions : 32x20x22cm  
Weight : 4Kg  
Shipping Weight : 10.5Kg  
Packing : 67x42x53cm

**4050** Dimensions : 47x20x22cm  
Weight : 6.5Kg  
Shipping Weight : 11.5Kg  
Packing : 80x60x44cm

**4400** Dimensions : 47x29x22cm  
Weight : 9Kg  
Shipping Weight : 16.5Kg  
Packing : 680x60x44cm

### Power Requirement:

115 or 230V, 50-60Hz

250VA max. for 4000/4050, 400VA max. for 4400

## Bibliography

- N. Bektas et alia: "Effect of phenolic acids on functions of rat aorta, vas deferens and on metabolic changes in streptozotocin-induced diabetes" *Indian J.Pharmacol.* 44 (2): 184-188, 2012
- A. Rizzo et alia: "In vitro effects of L-arginine on spontaneous and Homocysteine-induced contractility of pregnant canine uterus" *Theriogenology* 76 (4): 715-720, 2011
- E. N. Gorbatova et alia: "In Vitro Effects of Pentifin on Some Neurotransmitter Systems in the Brain" *Bull. Exper. Biology & Medicine* 136 (2): 174-175, 2003
- G. Re et alia: "Identification of Functional  $\alpha$ -Adrenoceptor Subtypes in the Bovine Female Genital Tract During Different Phases of the Oestrous Cycle" *Vet. es. Communications* 26 3): 479-494, 2002

## Multiplexing Pulse Booster

Cat. No. 3165

### General

The 3165 Multiplexing Pulse Booster is a useful complement to any stimulator, delivering up to 800mA of constant current to up to four in-vitro preparations (e.g., smooth muscles) at the same time.

The Multiplexing Pulse Booster has been designed to obviate the inconveniences connected to the use of single-channel stimulators, that lack the independent output connections and the individual adjustment capability to deliver pulses of preset intensity to more than one preparation.

It is especially useful when "field electrodes" and other low impedance stimulation arrangements are used.

Bear in mind that the one-channel stimulator can be conveniently replaced by a data acquisition system, as our **17304 DataCapsule-Evo!**



**Four in-vitro preparations  
can be driven by a single  
one-channel stimulator**

### Main Features

- High Power constant current: up to 800 mA
- Independent Isolated Circuits to eliminate interference
- Unipolar or Bipolar Mode
- Adequate Voltage (45V) enabling stimulation by field electrodes of most in-vitro preparations
- Continuous Monitoring of the preset current flowing through each preparation



## Instrument Description

### The 3165 features:

- High Power, digitally adjustable constant current: up to 800 mA
- Adequate Voltage (45V) which enables stimulation by field electrodes of most in-vitro preparations described in the literature
- Unipolar or Bipolar Mode
- Independent Isolated Circuits to eliminate interference
- Continuous Monitoring of the preset current flowing through each preparation

The current level of each channel is set via its individual 3-digit thumb-wheel switch.

The current output is adjustable in each channel to equal or different levels in the range 0-799mA in 1mA steps.

**These current levels are independent of the Stimulator output voltage.**

The pulse mode, either unipolar or bipolar, can be selected on one or more channels.

## Optional Timer

An optional **Timer (Cat. 3175)** can be supplied, housed in its individual mini-box, to enable the Pulse Booster to deliver pulse trains, when the Stimulator lacks this feature.

This timer is provided with both train and pause-between-trains duration adjustments. Both adjustment time-scales span the interval 0-999 seconds in 1 second steps.

A standard field electrode pair (Cat. 3160) can be supplied. Special electrodes can be designed and manufactured on request.

Please ask for details!

## Connection to Digital Recorder

A one-channel stimulator can be conveniently replaced by a data acquisition system, as - for example - our **17304 DataCapsule-Evo!**

## Ordering Information

**3165**      **MULTIPLEXING PULSE BOOSTER,**  
complete

**3165-302**      Instruction Manual

**E-PE 015**      Connection Cable

**E-WP 008**      Power Cord

### Optional

**3175**      Timer for 3165

### PHYSICAL

**Power Requirement**      115/230 V, 50/60 Hz, 30W

**Dimensions**      26(w)x30(d)x12(h)

**Weight**      4.4Kg

**Shipping Weight**      6.5Kg approx.

**Packing**      46x38x27cm

## Bibliography

- D. Currò: "**Voltage-gated calcium channels involved in the inhibitory motor responses and vasoactive intestinal polypeptide release in the rat gastric fundus**" *Eur. J. Pharmacol.* 628: 207-213, 2010
- F. Guagnini et alia: "**Tolerance to cannabinoid response on the myenteric plexus of guinea-pig ileum and human small intestinal strips**" *Br. J. Pharmacol.* 148, 1165-1173, 2006
- F. Borrelli et alia: "**Effect of Boswellia serrata on intestinal motility in rodents: inhibition of diarrhoea without constipation**" *B. J. Pharmacol.* 148, 553-560, 2006
- M.G. Matera et alia: "**Immune Sensitization of Equine Bronchus: Glutathione, IL-1 $\beta$  Expression and Tissue Responsiveness**" *Respir. Res.* 6 (1): 104, 2005
- S. Tambaro et alia: "**Evaluation of Tamsulosin and Alfuzosin Activity in the Rat Vas Deferens: Relevance to Ejaculation Delays**" *J. Pharmacol. Exper. Therap.* 312: 710-717, 2005
- S. Ruiu et alia: "**Synthesis and Characterization of NESS 0327: a Novel Putative Antagonist of the CB1 Cannabinoid Receptor**" *J. Pharmacol. Exper. Therap.*, 2003
- D. Licheri et alia: "**Long-Term Voluntary Ethanol Consumption Induces Impairment of the Mechanical Performance in the Papillary Muscle of Sardinian Alcohol-Preferring Rats**" *Alcohol and Alcoholism* 36 (1): 44-47, 2001

# Superfusion System

Cat. No. 14900

## General

Neurotransmitter release is the major step of neurotransmission. Abnormalities in neurotransmitter release have been proposed to be involved in many pathological conditions.

Therefore, understanding the physiological mechanisms of transmitter release and how the process can be modified by pathological states is essential to develop therapeutically useful pharmacological agents.

**UGO BASILE 14900 Superfusion System has been especially designed to perform release studies from synaptosomes**, although brain slices can be employed as well.

On the other hand, presynaptic nerve terminals are the sites where release specifically occurs; therefore superfusion of synaptosomes is best suited to explore presynaptic events.

Superfused synaptosomes are the preparation of choice to study release-regulating presynaptic receptors and to explore the intimate mechanisms of neurotransmitter release.



## RAITERI'S METHOD

## Synaptosomes Release Studies

## Main Features

- Specifically designed to perform release studies from synaptosomes
- Brain slices can be employed as well
- More than 300 full papers using superfused synaptosomes have been published

## Introduction

UGO BASILE **14900 Superfusion System** is a semi-automated version of that originally developed in Raiteri's laboratory, where about 200 papers have been published exploiting the technique.

We have developed this Superfusion System in order to make commercially available an instrument in which the original design of the superfusion chambers has remained intact.

The 14900 Superfusion System consists of 12 parallel open superfusion chambers with 12 upper reservoirs, all thermoregulated by a water-jacket. Prewarmed oxygenated media of the desired composition can be concomitantly delivered from the reservoirs to the superfusion chambers.

Synaptosomes are accommodated as very thin layers on microporous filters placed on glass filter supports.



Superfusion is provided by a multi-channel peristaltic pump and superfusate samples are directly collected into scintillation vials.

## Physical

|                   |  |
|-------------------|--|
| Weight            | 34Kg (complete assembly)   |
| Shipping Weight   | 48Kg   |
| Dimensions        | 14900-001: 38(w)x30(d)x13(h)cm<br>14900-002: 46(w)x28(d)x60(h)cm |
| Packing           | 1 box 80x60x44cm<br>1 box 62x65x84                               |
| Power Requirement | 115 or 230V, 50/60Hz, 100W max.                                  |

## Bibliography

### Method Paper:

- M. Raiteri, F. Angelini, G. Levi: "A simple apparatus for studying the release of neurotransmitters from synaptosomes" *Eur. J. Pharmacol.* 25: 411-414, 1974

### Papers quoting 14900:

- A. Pittaluga et alia: "Effects of the neoclerodane Hardwickiic acid on the presynaptic opioid receptors which modulate noradrenaline and dopamine release in mouse central nervous system" *Neurochemistry Intl.* 62 (4): 354-359, 2013
- S. Zucchini et alia: "Increased excitability in tat-transgenic mice: Role of tat in HIV-related neurological disorders" *Neurobiology of Disease*: available only 2013
- F. Giribaldi et alia: "Group I metabotropic glutamate autoreceptors induce abnormal glutamate exocytosis in a mouse model of amyotrophic lateral sclerosis" *Neuropharmacology* 66: 253-263, 2013
- J. Marrocco et alia: "Anxiety-Like Behavior of Prenatally Stressed Rats Is Associated with a Selective Reduction of Glutamate Release in the Ventral Hippocampus" *J. neuroscience* 32 (48): 17143-17154, 2012
- C. Romei et alia: "The GABAB receptor antagonists CGP35348 and CGP52432 inhibit glycine exocytosis: Study with GABAB1- and GABAB2-deficient mice" *Pharmacological Res.* 61: 547-552, 2010
- M. Grilli et alia: "Specific Inhibitory Effect of Amyloid- $\beta$  on Presynaptic Muscarinic Receptor Subtypes Modulating Neurotransmitter Release in the Rat Nucleus Accumbens" *Neuroscience* 167: 482-489, 2010
- G. Bonanno et alia: "Release of [3H]D-aspartate induced by K<sup>+</sup>-stimulation is increased in the cervical spinal cord of the wobbler mouse: a model of motor neuron disease" *Neurochemistry Intl.* 55: 302-306, 2009
- M. Grilli et alia: "Release-enhancing pre-synaptic muscarinic and nicotinic receptors co-exist and interact on dopaminergic nerve endings of rat nucleus accumbens" *J. Neurochemistry* 105 (6): 2205-2213, 2008

*In addition, more than 300 full papers using superfused synaptosomes have been published*

## Ordering Information

**14900 SUPERFUSION SYSTEM (Raiteri's method)**, standard package, including:-

|                  |  |
|------------------|--|
| <b>14900-001</b> | Electronic Unit  |
| <b>14900-002</b> | Superfusion Bath Complete Assembly, including upper & lower chambers, valves, set of tubes, etc. |
| <b>14900-004</b> | Suction Pump   |
| <b>14900-302</b> | Instruction Manual   |
| <b>14900-328</b> | Set of Phials  |
| <b>14900-338</b> | Set of Filters   |
| <b>14900-325</b> | Phial Rack   |
| <b>14900-302</b> | Drain Pan  |
| <b>E-WP008</b>   | Mains Cord   |

### Optional:

|                  |  |
|------------------|--|
| <b>14900-003</b> | Water Circulator/Heater                              |
| <b>14900-015</b> | Masterflex Peristaltic Pump, 12 channels, expandable |
| <b>14900-024</b> | Masterflex Peristaltic Pump, 24 channels             |



## Isometric Transducers

Cat. No. 7003 / 7004 / 7005 & 7010

### General

The three models 7003-7004-7005 cover the range from 0 to 50 g (see table on the facing page). The high sensitivity 7010 is designed for the mg range.

The force exerted on a hollow carbon fibre beam is converted into proportional electric signal via strain-gauges, conveniently wired in Wheatstone bridge circuit.

### Model Selection

Ugo Basile transducers are of robust construction and can withstand forces of up to 5-10 times the rated value.

It is possible to use 7003 which is generally used for trachea rings or artery strips, where forces of 5-10 grams are involved, by operating at minimum amplifier sensitivity; however, the cantilever will deflect with a load of the mentioned magnitude

Generally speaking, it is advisable to use a stiff transducer, operating at high amplifier sensitivity, and use the most sensitive transducer only when



The picture shows an **Isometric Transducer** (right) & an **Isotonic Transducer** (left), see separate datasheet

### Also available from Ugo Basile:

- Tissue Baths, 1, 2, 4-chambers
- Digital Recorder DataCapsule-Evo
- Electrodes & Stimulators

### Main Features

- Ugo Basile transducers have been designed for precise measurement of force in muscular preparations under isometric conditions
- An Isometric Transducer measures changes in force at constant length whereas an Isotonic Transducer is basically a displacement meter under constant load

## Isometric Transducer Specifications

| Model                              | 7010              | 7003              | 7004              | 7005              |
|------------------------------------|-------------------|-------------------|-------------------|-------------------|
| <b>Electrical</b>                  |                   |                   |                   |                   |
| Excitation Voltage (max.)          | 6V                | 6V                | 6V                | 6V                |
| Excitation Voltage (typical)       | 3V                | 3V                | 3V                | 3V                |
| Sensitivity ( $\mu V$ per g per V) | 110               | 70                | 25                | 10                |
| Non linearity & Hysteresis         | +/-3%             | +/-3%             | +/-3%             | +/-3%             |
| <b>Mechanical</b>                  |                   |                   |                   |                   |
| <b>Force Range</b>                 | <b>0-800 mg</b>   | <b>0-2g</b>       | <b>0-10g</b>      | <b>0-50g</b>      |
| Overload Rating                    | 5g                | 20g               | 50g               | 200g              |
| Moment of Inertia                  | 7gcm <sup>2</sup> | 7gcm <sup>2</sup> | 7gcm <sup>2</sup> | 7gcm <sup>2</sup> |
| Lever Arm Displacement             | 0.5 mm/g          | 0.3 mm/g          | 0.1 mm/g          | 0.06 mm/g         |
| <b>Physical</b>                    |                   |                   |                   |                   |
| Weight                             | 270g              | 270g              | 270g              | 270g              |
| Shipping Weight                    | 900g              | 900g              | 900g              | 900g              |
| Packing                            | 29x26x29cm        |                   |                   |                   |

## Compatibility

**Before ordering, check the connection compatibility** of your amplifier/recording system.

The Isometric & Isotonic Transducers are normally supplied with a connector designed for Ugo Basile Data-Capsule-Evo Recorder (see datasheet).

If the customer wishes to make use of other recording apparatus, the transducers can be supplied with appropriate connectors on request: we will be glad to provide transducer with different connectors, if available, or to provide wiring information and instruction.

## Ordering Information

- 7003** Isometric Force Transducer , type DY1
- 7004** Isometric Force Transducer , type DY2
- 7005** Isometric Force Transducer , type DY3
- 7010** High-Sensitivity Transducer , type DY0

## Bibliography

### Isometric Transducers 7003, 7004, 7005

- H. Ellers et alia: "Pungent General Anesthetics Activate Transient Receptor Potential-A1 to Produce Hyperalgesia and Neurogenic Bronchoconstriction" *Anesthesiology* 112: 1452-63, 2010
- A. Rizzo et alia: "Effects of rosiglitazone, a PPAR- $\alpha$  agonist, on the contractility of bovine uterus in vitro" *J. vet. Pharmacol. Therap.* 32, 548-551, 2009
- L. Natale et alia: "Interleukins 1 Beta and 6 Induce Functional Alteration of Rat Colonic Motility: An In Vitro Study" *Eur. J. Clin. Investigation* 33 (8): 704-712, 2003
- D. Mitolo-Chieppa et alia: "Involvement of  $\kappa$ -Opioid Receptors in Peripheral Response to Nerve Stimulation in  $\kappa$ -Opioid Receptor Knockout Mice" *Autonomic & Autacoid Pharmacology* 22:4: 233-239, 2002
- M.R. Accomazzo et alia: "Leukotriene D4-Induced Activation of Smooth-Muscle Cells From Human Bronchi Is Partly Ca<sup>2+</sup>-Independent" *Am. J. Respir. Crit. Care Med.* 163:1: 266-272, 2001
- M. Shalev et alia: "Stimulation of P2y Purinoceptors Induces, Via Nitric Oxide Production, endothelium-Dependent Relaxation of Human isolated Corpus Cavernosum" *J. Urol.* 161: 955-959, 1999
- M.C. Breschi et alia: "Effects of Noise Stress on EFS-Mediated Cholinergic and Inhibitory NANC Responses in Tracheae from Normal and Sensitized Guinea-Pigs" *J. Autonomic Pharmacol.* 17:6: 353-363, 1997
- M.K. Sim et alia: "Presence of an Endothelial Esterase in the Rat Aorta: Effects on the Actions of Ester and Non-Ester Muscarinic Antagonists" *Endothelium* 1: 109-114, 1993

### High-Sensitivity Transducer 7010

- L.W. Tait et alia: "Hagfish natriuretic peptide changes urine flow rates and vascular tensions in a hagfish" *Comparative Biochemistry and Physiology C* (150) 45-49, 2010
- G. Foldi et alia et alia: "Activity of sap from Croton lechleri on rat vascular and gastric smooth muscles" *Phytomedicine* 16: 768-775, 2009

# Isotonic Transducer

Cat. No. 7006

## General

The 7006 Isotonic Transducer basically consists of a carbonfibre lever arm which pivots on the shaft of a Hall-effect rotary motion transducer of original design.

The arm is balanced by an adjustable counterweight of tungsten alloy.

**It is possible to carry out experiments on extremely small muscle fibres**, which can be held under a tension of as little as 100-200 mg so that minimal force and consequent displacement alterations can be recorded.

The lever arm balancing is provided by a tungsten alloy counterweight which can be shifted by turning its knurled section.

This load is monitored by the counterweight rim moving along a scale calibrated in grams.



The picture shows an **Isotonic Transducer** (left) & an **Isometric Transducer** (right), see separate datasheet

## Also available from Ugo Basile:

- Tissue Baths, 1, 2, 4-chambers
- Digital Recorder DataCapsule-Evo
- Electrodes & Stimulators

## Main Features

- Ugo Basile Isotonic Transducer is specially designed for investigating isotonic contractions in isolated organs, particularly smooth muscle, amphibian hearts, etc.
- An Isotonic Transducer is basically a displacement meter under constant load, whereas an Isometric transducer measures changes in force at constant length



## Isotonic Transducer Specifications

|                           |  |
|---------------------------|--|
| <b>Voltage Output</b>     | 300µV per mm displacement of lever arm tip |
| <b>Linearity</b>          | ± 2% to ± 15 ° rotation                    |
| <b>Excitation Voltage</b> | 6 ÷ 15V                                    |
| <b>Excitation Current</b> | 20mA (constant in the range 6 ÷ 15V)       |
| <b>Operating Range</b>    | ± 15° about the centre                     |
| <b>Lever Arm Length</b>   | 10cm                                       |
| <b>Lever Arm Travel</b>   | 6cm  |
| <b>Breakaway Torque</b>   | less than 0.1g x cm                        |
| <b>Moment of Inertia</b>  | 35 g x cm <sup>2</sup>                     |
| <b>Overall Dimensions</b> | 16.5x5.5x11cm (excl. removable handle)     |
| <b>Weight</b>             | 0.35Kg                                     |
| <b>Shipping Weight</b>    | 1.60Kg                                     |
| <b>Packing</b>            | 29x26x29cm                                 |

## Compatibility

**Before ordering, check the connection compatibility** of your amplifier/recording system.

The Isometric & Isotonic Transducers are normally supplied with a connector designed for Ugo Basile Data-Capsule-Evo Recorder (see datasheet).

If the customer wishes to make use of other recording apparatus, the transducers can be supplied with appropriate connectors on request: we will be glad to provide transducer with different connectors, if available, or to provide wiring information and instruction.

## Ordering Information

**7006** Isotonic Transducer

## Bibliography

- O. E. Kiroglu et alia: "The effects of thiol modulators on nitrergic nerve- and S-nitrosothiols-induced relaxation in duodenum" *J. of Basic and Clinical Physiol. & Pharmacol.* 0 (0): 1-8, 2013
- M. Bucci et alia: "Cross-talk between toll-like receptor 4 (TLR4) and proteinase-activated receptor 2 (PAR2) is involved in vascular function" *Br. J. Pharmacol.* 168 (2): 411-420, 2013
- C. Jelen et alia: "Bone scaffolds with homogeneous and discrete gradient mechanical properties" *Materials Science & Engineering: C* 33 (1): 28-36, 2013
- M. Volta et alia: "Pharmacological profile and antiparkinsonian properties of the novel nociceptin/orphanin FQ receptor antagonist 1-[1-cyclooctylmethyl-5-(1-hydroxy-1-methyl-ethyl)-1,2,3,6-tetrahydropyridin-4-yl]-3-ethyl-1,3-dihydro-benzimidazol-2-one (GF-4)" *Peptides* 31:1194-1204, 2010
- P.U. Ertug: "Protective effect of quercetin, a polyphenolic compound, on mouse corpus cavernosum" *Fundamental & Clinical Pharmacology* 24: 223-232, 2010
- O. Desire et alia: "Antispasmodic and antioxidant activities of fractions and bioactive constituent davidigenin isolated from *Mascarenhasia arborescens*" *J. Ethnopharmacology. J. Pharmacol.* Accepted: May 2010, 2004
- D. Currò et alia: "Voltage-gated calcium channels involved in the inhibitory motor responses and asoactive intestinal polypeptide release in the rat gastric fundus" *Eur. J. Pharmacol.* 628: 207-213, 2010
- C. Belloli et alia: "Adrenergic Regulation of Vascular Smooth Muscle Tone in Calf Digital Artery" *J. Vet. Pharmacol. Therap.* 27:4: 247-254, 2004
- F. Carpi et alia: "Electromechanical Characterisation of Dielectric Elastomer Planar Actuators: Comparative Evaluation of Different Electrode materials and Different Counterloads" *Sensors and Actuators.* 107: 85-95, 2003
- P. Tucci et alia: "Cyclo-Oxygenase- and Capsaicin-Sensitive Afferent Fibres Affect Beta-Adrenoceptor-Evoked Response in the Rat Urinary Bladder" *Pharmacology* 64: 57-62, 2002
- P.C. Moser et alia: "SL65.0155, A Novel 5-Hydroxytryptamine4 Receptor Partial Agonist with Potent Cognition-Enhancing Properties" *J. Pharmacol. Exper. Therap.* 302:2: 731-741, 2002
- P. Tucci et alia: "Effects of Natural Tachykinins on Ovine Lower Urinary Tract Smooth Muscle" *J. Autonomic Pharmacol.* 21:2: 79-84, 2001
- C.M.Q. Jesus-Morais et alia: "Yangambin, a Lignan Obtained from *Ocotea duckei*, Differentiates Putative PAF Receptor Subtypes in the Gastrointestinal Tract of Rats" *Planta Med.* 66:4: 211-216, 2000

# DataCapsule-Evo Digital Recorder

Cat. No. 17304


## General

The new DataCapsule-Evo, powered by iWorx, is a new general purpose, four channel data acquisition system that provides high resolution and sensitivity over conventional recorders.

This versatile digital recorder is a unique system, in that each channel is independent, each having its own analog-to-digital converter.

In addition, each channel input is equipped with the appropriate filters and amplification required.

Set-up is plug-and play easy, with connection to PC or MAC computers via USB interface: data acquisition is accomplished via the versatile LabScribe Software provided with the system.

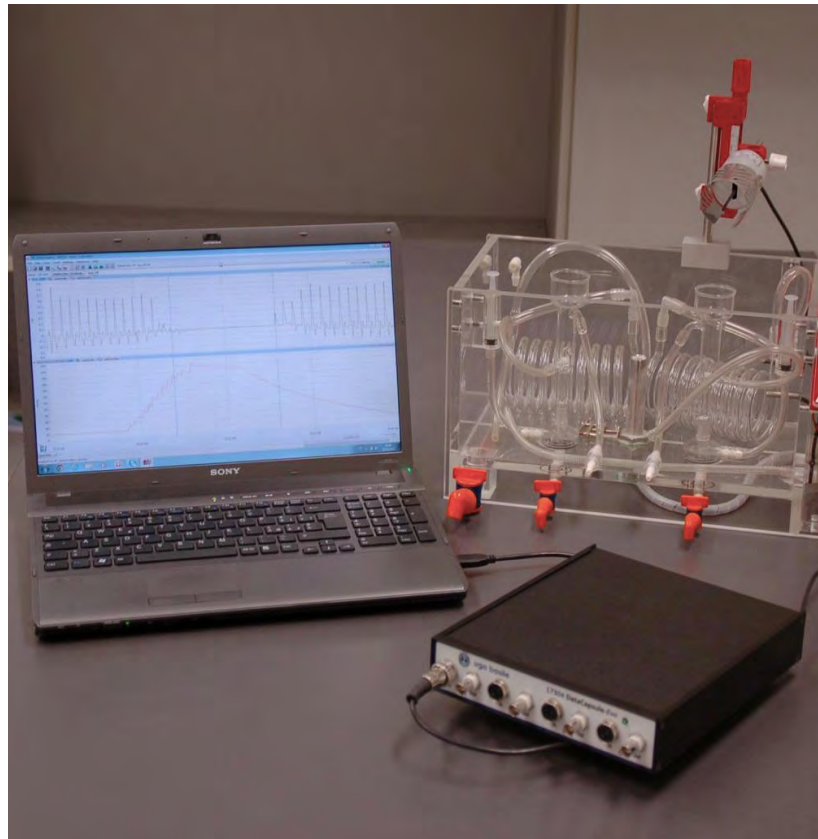

**LabScribe2™**  
Software on-board

**5KHz Sampling  
Speed**

**4-CHANNEL DATA  
ACQUISITION SYSTEM**

**Each channel is  
independent**

## Main Features

- USB connection to PC and MAC
- Connectors for most transducers
- DIN inputs and BNC outputs
- Input trigger to start recording
- High resolution and sensitivity
- Each channel equipped with programmable stimulator

## Connections and Specifications

All four channels have two connectors, an 8-pin DIN (input) and a BNC (output). Cable adaptors allow connection of a variety of sensors and transducers.

Each channel is independent, having its own 24 bit analog-to-digital converter.

The maximum sampling speed is 5KHz on all channels simultaneously.

## Analog Outputs

Each channel of the 17304 is equipped with an independently programmable 16bit, +/-10V stimulator. All stimulus parameters may be controlled via the LabScribe Software, from the stimulator tool bar.

Eight digital outputs are available to control external devices; programming the output lines is point-and-click easy.

## Resolution and Noise

The 17304 features high resolution, combined with an exceptionally low noise (as low as 32 uV).

## Software and Data Management

The DataCapsule-Evo setup is plug-and-play easy with connection to PC or MAC computers via the popular USB interface.

Recorded data are managed by the versatile LabScribe2 Software, featuring optimized scaling of displayed data: time base or y-axis scaling can also be zoomed in or out with a single click of the mouse.

Keyboard input from the user may be time locked to the data; annotations may be positioned in the data, just as you would write on chart paper!

Twenty-four off-line calculations are also supported, including Max-Min, Slope at a Point, and Mean.

Any view of the data can be exported to the disk as a text file or graphic. This capability is ideal for post calculation in programs like Excel™ or MatLab™.

Of course data from any window in the program can always be printed.

## DataCapsule-Evo Specifications

|                    |   |
|--------------------|---|
| Software           | : LabScribe2 onboard  |
| Display            | : Real-time, user-definable Screen Time independent of Sample Rate, user-Definable Units, AugoScale, Full Scale or User-Defined Scale |
| Power Requirements | : 115/230V – 50/60 Hz, 22W max  |

## Analog Input

|                    |  |
|--------------------|--|
| Analog Inputs      | : 4 independent amplified channels   |
| Input Resolutions  | : 24 bits differential   |
| Sample Speed       | : 1, 2, 5, 10, 20, 50, 100, 200, 500, 1K, 2K, 5K samples/second                        |
| Input Range        | : $\pm 10$ V   |
| Excitation Voltage | : $\pm 5$ VDC @ 50mA per channel   |
| Input Impedance    | : 667 Gig $\Omega$ typ./182 Gig $\Omega$ Min   |
| Trigger Mode       | : external trigger /TTL or Contact Closure), Threshold Trigger from Data, User Trigger |
| System Noise       | : 32 microvolts  |

## Analog Output

|                    |  |
|--------------------|--|
| Analog Outputs     | : 4  |
| Output Resolution  | : 16 bits                                      |
| Output Range       | : $\pm 10$ VDC                                 |
| Stimulator Modes   | : Pulse, Train, Constant, Step, Ramp, Triangle |
| Time Step          | : 0.04ms, 0.4ms, 4ms                           |
| Pulse Width (max.) | : 1.2s, 12s, 120s                              |
| Frequency          | : 12.5KHz, 1.25KHz, 125KHz                     |

## Physical

|                     |                               |
|---------------------|-------------------------------|
| Power Requirements  | : 115/230V - 50/60Hz, 22W max |
| Instrument Size     | : 22(w)x26(d)x5(h)cm          |
| Net/Shipping Weight | : 2.0Kg / 4.0Kg               |

## Ordering Information

**17304 DataCapsule-Evo**  
**4-Channel Digital Recorder**, standard package, including LabScribe2™ Software

## Transducers

**The DataCapsule can be connected to a variety of transducers.**

Among the ones offered by Ugo Basile:

|                |                                       |
|----------------|---------------------------------------|
| <b>7003-F</b>  | Isometric Force Transducer, type DY1  |
| <b>7004-F</b>  | Isometric Force Transducer, type DY2  |
| <b>7005-F</b>  | Isometric Force Transducer, type DY3  |
| <b>7010-F</b>  | High Sensitivity Transducer, type DY0 |
| <b>7006-F</b>  | Isotonic Transducer                   |
| <b>17844-F</b> | Pressure Transducer                   |

# ECT Unit

Cat. No. 57800

## General

The ECT apparatus is specially designed for neurochemical and neuropharmacological research.

A constant current output is used, which ensures reproducible results and accurate determination of the EC threshold while also pinpointing any variations in the threshold, brought about by drugs having a specific action on the cortex and subcortical regions.

The shock parameters have been selected after consulting the most recent literature, to supply the most suitable range when operating with mice and rats.

Consistent reproducible current levels are produced by feedback circuitry that adjust for variance in impedance of the contact from animal to animal.

The Electroconvulsive Device is supplied with auricular (ear lobe) electrodes.



**DESIGNED FOR  
INDUCING  
CONVULSIONS IN  
RESEARCH ANIMALS**

**FOR NEUROCHEMICAL  
&  
NEUROPHARMACOLOGICAL  
RESEARCH**

## Particularly useful for:-

- General screening of potentially neurotropic substances
- Evaluating the depressant or stimulating action of drugs on the CNS
- Endocrinological investigations on the relationship between the nervous system and the hypophysis



## General

Consistent reproducible current levels are produced by feedback circuitry that adjust for variance in impedance of the contact from animal to animal.

The impedance of the animal can be previously measured and displayed, and a warning signal flashes if the impedance is too great to deliver the desired current level.

The special output circuit enables any type of electrode to be used.

The **auricular electrodes 57800-002**, supplied with the standard package, allow a single operator to deliver shock to a number of animals in a short time.



The above picture features **Corneal Electrodes Cat. 57800-003**, which can be provided as **optional**.

Different types of electrodes can be manufactured on request.

## Specifications



Rectangular Positive

Pulse : by H.V. transformer  
 Constant Current : controlled by a feedback network  
 Pulse Rise&Fall Time : 20µs  
 Pulse Width (ms) : 0.1 to 0.9 in 0.1ms steps  $\pm 1\%$   
 Frequency (pulses/s) : 1-299 in 1 pulse/s steps  $\pm 1\%$   
 Shock Duration : 0.1 to 9.9 in 0.1s steps  $\pm 1\%$   
 Pulse Voltage : 2.5KV max.  
 Current Range : 0-99mA in 1mA steps  $\pm 2\%$   
 Output Resistance : min 00hm - max. 25KOhm (at max. current)  
 KOhm Display : 0-199 KOhm - 1KOhm resolution  
 Power Requirements : 115/230V - 50/60Hz - 70VA

**WARNING: due to HIGH VOLTAGE involved, the operator should always wear rubber gloves when handling the electrodes.**

## Bipolar Inverter 57800-010

An optional Biphasic Unit may be placed between the animal and the Electroconvulsive Device to invert every second pulse. Maximum frequency in this case becomes 100 Hz.

## ECT Monitor 57800-015

When connection to an oscilloscope or data acquisition system, this useful accessory is required to guarantee a simple and safe way to monitor the ECT output.



The risk of damage to the ECT Unit due to accidental wrong connections is avoided when using the ECT Monitor.

## Ordering Information

**57800 ECT Unit**, standard package including:-

- 57800-001** Pulse Generator
- 57800-002** Set of Auricular Electrodes
- 57800-302** Instruction Manual (on CD)
- E-WP 008** Mains Cord

## Accessories and Spares

- 57800-003** Set of Corneal Electrodes
- 57800-320** Set of 4 Felt Pads for Auricular Electrodes
- 57800-010** Bipolar Inverter
- 57800-015** ECT Monitor

## Physical

Instrument Size 27(W)x37(D)x13(H)cm  
 Weight 3.4Kg  
 Packing 45x34x26cm  
 Shipping Weight 5Kg

## Bibliography

- K.M. Ryan et alia: "Electroconvulsive Stimulation Alters Levels of BDNF-associated microRNAs" *Neuroscience Letters* available May 2013
- N.K. Thomas et alia: "Triheptanoin in Acute Mouse Seizure Models" *Epilepsy Research* 99(3): 312-317, 2012
- S.A. Epps et alia: "Seizure Susceptibility and Epileptogenesis in a Rat Model of Epilepsy and Depression Co-Morbidity" *Neuropsychopharmacol.* 37: 2756-2763, 2012
- M. Gasior et alia: "Anticonvulsant and Proconvulsant Actions of 2-deoxy-D-glucose" *Epilepsia* 1-10, 2010
- A.L. Hartman et alia: "Efficacy of the Ketogenic Diet in the 6-Hz Seizure Test" *Epilepsia* 49(2): 334-339, 2008

## Lesion Making Device

Cat. No. 53500

### General

This compact, **solid state DC Lesion Maker** has been designed for the production of localized lesions in small animals, when direct current (DC) is preferred to RF.

It features a regulated power supply combined with a constant DC current generator which operates on either continuous or timed mode.

The Lesion Making Device provides constant DC current in mA from 10  $\mu$ A to 99 mA. The pulse duration may be timed by the instrument between 1 and 99 seconds, or manually controlled.

The current generator is protected against short circuit, preventing the electronics to get damaged due to the electrodes coming accidentally in contact with each other.

Particular emphasis has been placed in the design of a good circuit output/ground insulation; this feature also minimizes spurious current field lines across the tissue, outside the pattern preset by the operator.



**New Model!**

**A precision instrument, which provides constant DC current in mA**

### Main Features

- Violation warning circuit
- Current Range : from 10  $\mu$ A to 99 mA
- 3 modes of Operation
- Digital setting of constant current and time duration
- Pulse Duration : timed between 1 and 99 seconds

## Controls



The instrument controls are all placed on the top panel; the parameter are set by two thumb-wheel switches:-

- **Current output adjustment**, in the range 10 $\mu$ A to 99mA
- **Pulse duration** from 0.1 to 99 seconds.

The mode of operation can be selected via a 3-position switch:-

- **Continuous**: the current flows through the preparation in a continuous mode
- **Stand-By**: the instrument is ready to operate but the output stage is not energized
- **Pre-set Duration**: the current flow is timed according to the setting

Three binding posts are located at the upper right of the Lesion Maker: either the red (+) and the black (-) can be connected to the lesion making electrode.

The other binding post is usually connected to a pad electrode with electrolyte on the preparation. Either red (+) or black (-) may be grounded via the green binding post.

## Led Indicators

Three LED indicators are embodied on the top panel:-

- **POWER** (green) which lights when the unit is ON
- **MONIT.** (red) which monitors the presence of lesion current
- **VIOL.** (yellow) which indicates when the current does not correspond to the setting

## Electrodes

Usual needle electrodes, prepared by the researcher according to his/her experimental needs can be used in conjunction with the 53500 Lesion Making Device.

We have the capability and will to manufacture electrodes based on the customer's request.

## Ordering Information

**53500** Lesion Making Device  
standard package, including:-

**53500-310** Set of 3 output plugs  
**53500-302** Instruction Manual  
**E-WP 008** Mains Cord

## Technical Specifications

|                    |  |
|--------------------|--|
| Current Range      | from 10 $\mu$ A to 99 mA                               |
| Pulse Duration     | timed between 1&99 seconds or manually controlled      |
| Compliance Voltage | 200 V DC   |
| Max. Electrode R   | 20M $\Omega$ (10 $\mu$ A) down to 2K $\Omega$ (100 mA) |
| Mains Supply       | 115 or 230V / 50-60 Hz                                 |
| Power Consumption  | 20 W max.  |

## Physical

|                 |               |
|-----------------|---------------|
| Dimensions      | 25x15x11 cm   |
| Weight          | 1.5Kg         |
| Shipping Weight | 2.8Kg approx. |
| Packing         | 45x34x26cm    |

## Bibliography

- S. Stroobants et alia: "**PIncreased gait variability in mice with small cerebellar cortex lesions and normal rotarod performance**" Behav. Brain Res. 241: 32-37, 2013
- L.B. Cruz et alia: "**PraEffect of the bone marrow cell transplantation on elevated plus-maze performance in hippocampal-injured mice**" Behav. Brain Res. available online Apr. 2013
- M.B. Gomes et alia: "**Glucose levels Observed in Daily Clinical Practice induce Endothelial Dysfunction in the Rabbit Macro- and Microcirculation**" Fund. & Clin. Pharmacol. 18 (3), 2004
- C. Hamani et alia: "**Bilateral Anterior Thalamic Nucleus Lesions and High-frequency Stimulation Are Protective against Pilocarpine-induced Seizures and Status Epilepticus**" Neurosurgery, 54 (1): 191-197, 2004
- T. Lee and J.J. Kim: "**Differential Effects of Cerebellar, Amygdalar, and Hippocampal Lesions on Classical Eyeblink Conditioning in Rats**" J. Neuro-science 24 (13): 3242-3250, 2004
- K.C. Bicego and L.G.S. Branco: "**Discrete Electrolytic lesion of the Preoptic Area Prevents LPS-Induced Behavioral Fever in Toads**" J. Exper. Biol. 205: 3513-3518, 2002

## Stereotaxic Instruments by Stoelting

Cat. No. 51600

### General

The Lab Standard™ Stereotaxic Instrument, manufactured by Stoelting, is ideal for researchers in need of a versatile, reliable instrument for stereotaxic procedures with small animals.

Precision alignment when using the Lab Standard™ ensures accurate placement of electrodes, micropipettes, and other devices.

The time-proven 'U'-Frame design concept, sturdy construction, and adaptability to most model species make this the best choice for a stereotaxic instrument.



**SLEEK, COMPACT DESIGN**

**ACCESSORIES AVAILABLE FOR USE  
WITH A WIDE VARIETY OF SMALL  
ANIMALS**



### Classic and Proven U-Frame Design

- Large, easy to read vernier scales. Scales are laser engraved — accurate to 100 microns
- Triple lead screws for fast positioning 80 mm of vertical, lateral and anterior-posterior travel
- Absolute lock at 90 degrees (vertical) Brass bushings in manipulator arm permit electrical grounding



Stoelting's Lab Standard™ offers several advantages over competing instruments:

### Easily Read Scales

All scales are oriented to be read easily from the open end of the 'U'. This is the position from which most scientists prefer to work. The numerals on the scales are larger, and therefore more easily read. The scale lines are laser engraved, to allow finest possible permanent marking of scales on all 3 axes. Precise alignment with facing vernier scales gives accurate resolution to 0.1mm.

### Smooth Movements

The Lab Standard's™ exclusive, triple lead screws allow the fastest positioning possible, consistent with lining up the scales easily at a given coordinate.

### Versatility of Positioning

The manipulator arm controls medio-lateral and vertical positioning via lead screws, and antero-posterior movement via dovetail slide movement, with 80 mm of travel possible in each direction. A Universal Joint allows the investigator to change the angle of the probe up to 90° in either the antero-posterior or medio-lateral planes. The improved locking mechanism on the Lab Standard™ will hold any angle position without slippage. And of course, it also provides an absolute lock at 90° vertical.

In addition, a swing joint allows the investigator to conveniently swing the manipulator arm and probe out of the way for performing a procedure — then reliably return the probe to the same point.

### Convenient for Electrophysiology

Integral brass bushings in the manipulator arm allow grounding directly to the closest metal on the manipulator arm — even the probe holder.

### Selection of Accessories

Species adaptors are available to fit rat, cat/monkey, dog/monkey, mouse, guinea pig and small bird. Probe holders and species adaptors for 'U' frame stereotaxic instruments from other manufacturers are generally compatible with the Lab Standard™ frame.

### Ordering Information

- 51600** Lab Standard w/18 Degree Earbars
- 51650** Lab Standard w/45 Earbars
- 51653** Dual Lab Standard Stereotaxic w/45 deg. Ear Bars
- 51603** Dual Lab Standard Stereotaxic w/18 deg. Ear Bars
- 51601** Lab Standard without Manipulator Arms

## INFUSION PUMPS

*by KDS*

Cat. No.5000

**SO ADVANCED THEY'RE SIMPLE !!**

### General

Ugo Basile presents an entirely new generation of micro-processor controlled syringe pumps. They are designed specifically for applications requiring high metering precision at low, pulse free flow rates.

KDS pumps, manufactured by KD Scientific Inc., U.S.A., provide a unique combination of sophisticated features and advanced microstepping motor-drive technology. The result? KDS pumps routinely perform many of the tasks that other pumps make you do manually. So you have more time for what's really important: your research.

KDS pumps are engineered by the designer of the best selling laboratory syringe pump, to ensure you of years of unsurpassed accuracy and reliability. In addition, you'll find they are extremely simple to set-up and use.

And surprisingly affordable.



### Setup is as easy as:

- Select syringe from displayed table
- Enter volume to be dispensed
- Enter flow rate, then press "start" button. It's that fast...and that simple!

## Common to all models

- A simple menu-driven set up without printed look-up tables **performs rate and volume control and automatic shut-off**. Just set the volume you want dispensed. Volume is tracked continuously on the LC display. Then, when the preset volume has been dispensed, the pump shuts off automatically.
- An **alphanumeric display helps eliminate reading errors**. Their easy-to-read display provides real-time readings using both parameters and values for clearer, mistake-free readings.
- You can control KDS pumps in many different ways. Built-in TTL and RS-232C interfaces permit easy external control.

## Operation

1. Find the syringe you use from the displayed table. Enter its code number.
2. Enter the volume to be dispensed
3. Enter the flow rate, then press the "start" button. It's that fast and simple! Your settings are permanently stored in memory – there's no need to re-enter them each day

## Ordering Information

| Cat. No. | Mode                    | N. of Syringes        | Dim. cm    | Weight Kg. |
|----------|-------------------------|-----------------------|------------|------------|
| KDS 100  | Infusion                | 1                     | 23x15.3x14 | 2.00       |
| KDS 101  | Infusion                | 2                     | 23x15x14   | 2.00       |
| KDS 120  | Push/pull               | 1+1                   | 23x15x14   | 2.00       |
| KDS 200  | Infusion                | 2                     | 28x23x14   | 4.00       |
| KDS 210  | Infusion/<br>Withdrawal | 2                     | 28x23x14   | 4.00       |
| KDS 220  | Infusion                | Multiple              | 28x30.5x14 | 4.25       |
| KDS 230  | Infusion/<br>Withdrawal | Multiple              | 28x30.5x14 | 4.25       |
| KDS 250  | Infusion                | 4<br>(different size) | 28x23x15.3 | 4.00       |
| KDS 260  | Push/pull               | 2+2                   | 28x23x14   | 4.25       |
| KDS 310  | Nano Pump               | 1                     | 2 items    | 2.00       |

## Flow Rates

### Models KDS 100 & KDS 120

| Syringe | Minimum  | Maximum    |
|---------|----------|------------|
| 10 µl   | 0.1 µl/h | 126.5 µl/h |
| 25 µl   | 0.1 µl/h | 318.8 µl/h |
| 50 µl   | 0.2 µl/h | 625 µl/h   |

|        |           |            |
|--------|-----------|------------|
| 100 µ  | 1.0 µl/h  | 1274 µl/h  |
| 250 µ  | 2.0 µl/h  | 3164 µl/h  |
| 500 µ  | 3.0 µl/h  | 6359 µl/h  |
| 1ml    | 0.01 ml/h | 13,2 ml/h  |
| 2,5 ml | 0.02 ml/h | 31,7 ml/h  |
| 3 ml   | 0.02 ml/h | 44.9 ml/h  |
| 5 ml   | 0.03 ml/h | 87.0 ml/h  |
| 10 ml  | 0.1 ml/h  | 125.0 ml/h |
| 20 ml  | 0.1 ml/h  | 219.0 ml/h |
| 30 ml  | 0.1 ml/h  | 282.0 ml/h |
| 60 ml  | 0.2 ml/h  | 426.0 ml/h |

### Model KDS 101

| Syringe | Minimum      | Maximum      |
|---------|--------------|--------------|
| 10 µl   | 0.001 µl/min | 0.350 µl/min |
| 25 µl   | 0.001 µl/min | 0.884 µl/min |
| 50 µl   | 0.001 µl/min | 1.759 µl/min |
| 100 µl  | 0.001 µl/min | 3.526 µl/min |
| 250 µl  | 0.01 µl/min  | 8.78 µl/min  |
| 500 µl  | 0.01 µl/min  | 17.65 µl/min |
| 1 ml    | 0.1 µl/min   | 35.2 µl/min  |
| 3 ml    | 0.1 µl/min   | 122.5 µl/min |
| 5 ml    | 0.1 µl/min   | 176.2 µl/min |
| 10 ml   | 0.001 µl/min | 0.351 µl/min |
| 20 ml   | 0.001 µl/min | 0.602 µl/min |
| 30 ml   | 0.001 µl/min | 0.773 µl/min |
| 60 ml   | 0.001 µl/min | 1.175 µl/min |

### Models KDS 200/220, KDS 210/230, KDS 250/260

| Syringe | Minimum    | Maximum      |
|---------|------------|--------------|
| 10 µl   | 0.001 µl/h | 21.1 µl/min  |
| 25 µl   | 0.003 µl/h | 53.15 µl/min |
| 50 µl   | 0.005 µl/h | 105.8 µl/min |
| 100 µl  | 0.009 µl/h | 212.6 µl/min |
| 250 µl  | 0.021 µl/h | 527.6 µl/min |
| 500 µl  | 0.042 µl/h | 1060 µl/min  |
| 1 ml    | 0.083 µl/h | 2119 µl/min  |
| 3 ml    | 0.288 µl/h | 7360 µl/min  |
| 5 ml    | 0.414 µl/h | 634 ml/h     |
| 10 ml   | 0.828 µl/h | 1270 ml/h    |
| 20 ml   | 1.414 µl/h | 2171 ml/h    |
| 30 ml   | 1.817 µl/h | 2789 ml/h    |
| 60 ml   | 2.757 µl/h | 4234 ml/h    |
| 140 ml  | 5.746 µl/h | 8834 ml/h    |

## Blood Pressure Recorder (non-invasive)

Cat. No. 58500 for Rats

Cat. No. 58600 for Mice

Cat. No. 58550 for Rats & Mice

### General

The BP RECORDER 58500 combines three main systems

- pressure generation-pressure monitoring system
- a pulse amplifier and
- a thermal-array analog & digital recording unit

with two auxiliary systems

- pulse rate measuring and recording
- microprocessor controlled functions to self diagnosis, calibration, signal filtering, signal storage.

### Instrument Description

Pressure is transmitted to the tail cuff; as soon the cuff pressure exceeds the diastolic pressure and starts to narrow the tail artery, the amplitude of the recorder pulse wave gradually decreases until the artery is completely constricted (ischemic), the graph becoming a straight line.

This point indicates the maximum internal pressure of the artery (**systolic pressure**) on the paper grid, on which the **actual pressure** of the system is **digitally printed in 10 mm Hg steps**.

At the end of the recording a second pressure measurement can be started, with decreasing pressure. The systolic pressure is indicated, this time, by the return of the pulse tracing.

The animal **pulse rate** can be assessed in real time by a pulse rate counter which picks the signal from the pulse transducer.



### INDIRECT MEASURING & RECORDING OF THE SYSTOLIC AND DIASTOLIC PRESSURE IN UNANAESTHETIZED RATS & MICE

### Main Features

- graphic printer
- graphic display
- analog output to digital recorders
- pulse transducers of superior performances
- analogue & digital recording of all experiment phases
- reliable pressure generator, providing smooth, stepless pressure build-up



## Animal Restrainers

A convenient animal restrainer is provided with the standard package. Our models are particularly suitable, being purposefully designed for this task, as they feature:-

- a conical "muzzle" to confine the animal head
- availability in 4 different diameters for rat and one for mouse, to fit various animal sizes
- telescope-adjustable length
- a quick fit/release back lid with an ample U-shaped tail slot
- convenient ventilation slots and selection of heat conductive materials, to guarantee body heat dissipation.

## Optional Rat Heater / Scanner

The **58000-845 Heating Box for Rats** is a compact temperature controlled "cupboard", inside dimension 57(w)x47(d)x20(h) cm, to lodge and prewarm 5 rats, each in its individual holder; **58000-840**, designed for mice, has the same dimensions, but it accommodates 6 mouse holders.



The **58000-850 Rat Scanner** is also available, combining the pre-warming features, with an electrical/pneumatical switch which enables connection of up to 5 rodents, tail cuff and pulse pick-up positioned on their tail, to scan their blood pressure in sequence.

Both Rat Scanner and Heating Boxes come complete with holders of selectable diameter.

## Ordering Information

- 58500** BP RECORDER, with accessories for **RAT**: 8mm pulse pick-up, 13mm cuff, 50mm holder
- 58600** BP RECORDER, with accessories for **MOUSE**: 3mm pulse pick-up, 6mm cuff, 30mm holder
- 58550** BP RECORDER, with accessories for **RAT & MOUSE**

Each BP Recorder includes as standard: dedicated software 52050-08, serial cable & USB adaptor, paper roll.

### Available Pulse Pick-Ups

- 58000-503** Pulse Pick-up for Mouse, diam. 3 mm
- 58000-504** Pulse Pick-up for Mouse, diam. 4 mm
- 58000-505** Pulse Pick-up for Rat, diam. 5 mm
- 58000-506** Pulse Pick-up for Rat, diam. 6 mm

- 58000-507** Pulse Pick-up for Rat, diam. 7 mm
- 58000-508** Pulse Pick-up for Rat, diam. 8 mm
- 58000-509** Pulse Pick-up for Rat, diam. 9 mm

### Available Tail Cuffs

- 58000-606** Tail Cuff for Mouse, diam. 6 mm
- 58000-609** Tail Cuff for Rat, diam. 9 mm
- 58000-611** Tail Cuff for Rat, diam. 11 mm
- 58000-613** Tail Cuff for Rat, diam. 13 mm

### Available Holders

- 58000-343** Mouse Holder, 30 mm I.D.
- 58000-344** Rat Holder, 40 mm I.D.
- 58000-345** Rat Holder, 50 mm I.D.
- 58000-346** Rat Holder, 60 mm I.D.
- 58000-348** Rat Holder, 80 mm I.D.

### Optional

- 58000-840** Mouse Heater, compl. with 6 mouse holders
- 58000-845** Rat Heater, complete with 5 rat holders of selectable I.D.\*
- 58000-850** Rat Scanner, complete with 5 rat holders of selectable I.D.\* °

\* if diameter is not specified, the 50mm size will be supplied

° pressure cuffs & pulse pick-ups are not included, and should be ordered separately

## Specifications

|                    |                              |
|--------------------|------------------------------|
| Pressure Range     | 50 to 290 mm Hg              |
| Power Requirements | 115 or 230 V, 50/60 Hz, 25 W |
| Weight (net)       | Kg 10.6                      |
| Shipping Weight    | Kg 15.0 approx.              |
| Dimensions         | 35x35x17(h)cm                |
| Packing dimensions | 80x60x44cm                   |

## Bibliography

- M. Gerold & H. Tschirky "Measurement of Blood Pressure in Unanaesthetized Rats" *Arzneimittelforschung* 18: 1285-287, 1968
- M. Gerold & H. Fünfschilling: "Abhängigkeit der Indirekten Blutdruckmessung an Ratten von der Größe der Kompressionsmanchetten" *Arzneimittelforschung* 21: 2071-2074, 1971.

## Papers quoting Ugo Basile Model

- L. Testai et alia: "The activation of mitochondrial BK potassium channels contributes to the protective effects of naringenin against myocardial ischemia/reperfusion injury" *Biochemical Pharmacol.* 85(11): 1634-1643, 2013
- A. Kolosov et alia: "Intravenous Injection of Leconotide, an Omega Conotoxin: Synergistic Antihyperalgesic Effects with Morphine in a Rat Model of Bone Cancer Pain" *Pain Medicine*: 12(6): 923-941, 2011
- M.A. Gouda et alia: "Synthesis and anti-hypertensive activity of novel sulphadimidine derivatives" *Med. Chem Res.*: 21(11): 3902-3906, 2011
- J. Tchekalarova et alia: "Diurnal rhythms of spontaneous recurrent seizures and behavioral alterations of Wistar and spontaneously hypertensive rats in the kainate model of epilepsy" *Epilepsy & Behavior* 17: 23-32, 2010
- C. Bolego et alia: "Selective estrogen receptor- $\alpha$  agonist provides widespread heart and vascular protection with enhanced endothelial progenitor cell mobilization in the absence of uterotrophic action" *FASEB Journal*: fj.09-139220, pub. online 2010

## Blood Pressure Transducer (invasive)

Cat. No. 17844

**Easy to fill**

**High accuracy**

**Robust, reusable transducer**

### Typical Applications

- Arterial or venous blood pressure measurement
- Connects to Data Acquisition Systems or Chart Recorders
- Urodynamic measurement
- Intrauterine Pressure Measurement
- Intracranial Pressure Measurement
- Catheterization
- Intensive Care Unit




### Main Features

- MPG Klasse II b, CE 0470
- Gold plated for easier cleaning
- Only wiping cleaning necessary
- Disinfection/Sterilisation with VIRKON (10 to 30 min) possible
- Short adapter cable with transducer + separate monitor cable
- Dome with "Snap-on" coupling
- Very high frequency response
- High overload protection (10.000 mm/Hg)
- Dome dry-coupled to the transducer

## Specifications

|   |  |
|---|--|
| <b>Pressure Range</b>                     | -20...+300mmHg                         |
| <b>Overpressure max.</b>                  | 10 000mmHg                             |
| <b>Sensitivity</b>                        | 50 $\mu$ V/V/cmHg                      |
| <b>Resonance Frequency</b>                | 300Hz typical<br>(Transducer and Dome) |
| <b>Electrical Excitation max.</b>         | 15V DC or AC                           |
| <b>Input Resistance<br/>(Input)</b>       | 7000Ohm                                |
| <b>Output Resistance<br/>(Output)</b>     | 10000Ohm                               |
| <b>Non-Linearity &amp;<br/>Hysteresis</b> | max. 0.5% FS                           |
| <b>Zero Balance</b>                       | max. $\pm$ 30mm/Hg                     |
| <b>Thermal Sensitivity Shift</b>          | 0.15% / °C                             |
| <b>Thermal Zero Shift</b>                 | max. 0.25mm/Hg/°C                      |
| <b>Operating Temperature<br/>Range</b>    | +10...+50°C                            |
| <b>Storage Semperature<br/>Range</b>      | -20...+70°C                            |
| <b>Insulation Resistance</b>              | min. 103MOhm                           |
| <b>Leakage Current</b>                    | max. 1.5 $\mu$ A at 250V-50Hz          |
| <b>High Voltage Resistance</b>            | 10KV between Dome and<br>Transducer    |
| <b>Length of Adapter Cable</b>            | ca. 30cm                               |
| <b>Length of Monitor cable</b>            | ca.250cm                               |
| <b>Connector</b>                          | see "compatibility"                    |

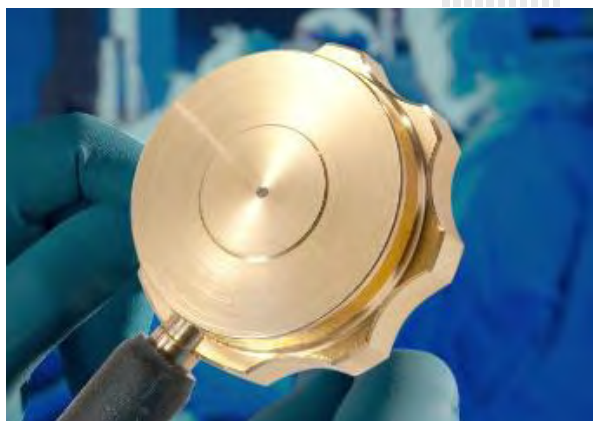

## Compatibility


**Before ordering, check the connection compatibility** of your amplifier/recording System.

The Pressure Transducers are normally supplied with a connector (type -F) designed for Ugo Basile **DataCapsule-Evo Recorder** (see datasheet).

If the customer wishes to make use of other recording apparatus, the transducers can be supplied with appropriate connector on request: we will be glad to provide transducers with different connectors, if available, or to provide wiring information and instruction.

## Domes

The 17844 is provided with a dome provided with with stop cock. The dome has wings, for easy fitting on the transducer.

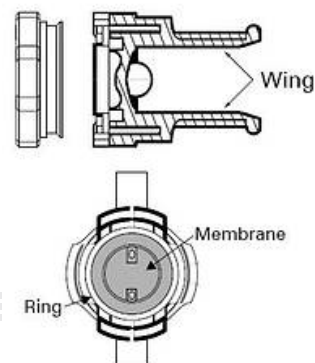

The dome should be filled bubbleless at max. pressure of 50mmHg

## Ordering Information

**17844** Pressure Transducer "Sensoror", type SP-844, complete with one dome 17844-001 lodged in its plastic case.

### Accessories

**17844-001** Clear Polycarbonate Dome (with Luer-Lock Fitting), complete with 3-way stopcock

**17844-002** Set of 10 Clear Polycarbonate Dome (with Luer-Lock Fitting), complete with 3-way stopcock



### Physical

|                     |   |
|---------------------|---|
| Weight              | 0.024Kg (without cable)<br>0.2Kg (with cable) |
| Shipping Weight     | 0.4Kg   |
| Shipping Dimensions | 46x38x27cm                                    |

# MouseOx *Plus*

## Pulse Oximeter for Mice and Rats

### General

The MouseOx® is the world's first and only patented **non-invasive** vital signs monitor, for small laboratory animals; specifically designed for mice, it can be used on larger rodents too!

The MouseOx and The MouseOx Plus® are being used by over 1,500 researchers and veterinarians from Universities, Pharmaceutical Companies, and CRO.

It is fully controlled by PC with a **user-friendly interface**.

The new **MouseOx Plus®** uses the same technology as the original MouseOx® but also includes significant improvements:

- the enhanced signal processing ability improves response to the motion of conscious subjects; the pulse signal is maintained and quickly re-acquired following significant movement.
- the modular software design allows the end user to purchase only the functionality that is needed.
- measurement of core body temperature is now available
- the optional Multiplexer makes it possible to monitor up to 16 animals (or 8 animals with temperature), using 1 MouseOx Plus.



Anesthetized Subjects

Conscious Subjects

MRI Compatible

**SMALL  
ANIMAL  
VITAL SIGNS  
MONITOR**

### Main Features

- Simple non-invasive sensor clips for mice and rats
- Monitor data in real time while recording
- USB plug-and-play, user-friendly interface
- High accuracy at heart rates up to 900 BPM
- Works on neonates through adults



## General

The **MouseOx Plus Small Animal Vital Signs Monitor** provides the following measurements:

- Arterial Oxygen Saturation
- Heart Rate
- Breath Rate
- Temperature (optional)
- Pulse Distention
- Breath Distention

MouseOx Plus® works with both mice and rats; there are 16 variations of the MouseOx® sensor available to accommodate various sensor placement options on mice and rats, ranging in size from neonatal mice to rats over 500gm. The subject must have a heart rate of at least 90 BPM and no greater than 900 BPM.

The MouseOx® oxygen saturation measurement has only been validated with mice and rats, but the instrument is being used in many research projects on subjects other than mice and rats. Some examples include Guinea pigs, hamsters, rabbits and small primates such as marmosets.

## Cardiopulmonary Data Recorder

When used as a Cardiopulmonary Data Recorder, the MouseOx Plus provides:

- Quick Check of Vital Signs
- Real-time Changes in Heart Rate, Breath Rate & O2 Saturation
- Oxygen Saturation During Hypoxemia
- Analog Data Output

## Surgery Monitor:

When used as a Surgery Monitor MouseOx Plus:

- Prevents Hypoxia During Surgery
- Titrates Mechanical Ventilation
- Ensures Proper Depth of Anesthesia
- Titrates Supplemental Oxygen

## Features:


- Immediate responding, beat-by-beat measurements
- High accuracy at heart rates up to 900 BPM and breath rates up to 600 BrPM
- Drawing of blood is not required for any reason
- Simple non-invasive sensor clip enables quick and easy attachment to the subject
- USB plug-and-play technology easily turns your Windows based PC into a low cost physiologic monitor
- Monitor data in real-time, while recording to a file
- Experiment event markers allow the user to mark important events in the data file
- Enhanced signal processing ability improves response to the motion of conscious subjects

## Standard Software and Options

The **Standard software** includes basic monitoring and parameter alarms for all of the vital signs provided by the MouseOx Plus; it is included with all MouseOx Plus systems and is intended for basic monitoring applications.



The **Premium Monitoring & Recording Software** includes trending charts, real time recording options, file markers for noting important events, and a quick averaging diagnostic feature for spot-checking

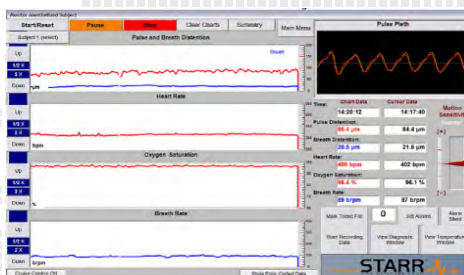

The **Conscious Applications Software** includes enhanced filters and control algorithms to allow the MouseOx® Plus to monitor conscious unrestrained subjects, and provides a subject activity measurement

**MRI software** allows for the use of the MRI sensor.

## Ordering Information

- 015000** MouseOxPlus System, Operation 110V \*
- 015001** MouseOxPlus System, Operation 230V \*
- 015007** Premium Monitoring & Recording Software
- 015017** MRI Module, including Software, 20' Sensor with 15' Copper Wire and 5' Fiber Optic, 2 Mouse Thigh Clips, 2 Rat Foot Clips
- 015002** Conscious Applications Module

### Sensors

\* Two sensors, selectable when ordering, are included free of charge with each MouseOxPlus System:

**CollarClip™** available in XS, S, M, L, XL, 2XL size

**ThroatClip™** available in XS, S, M, L, XL, 2XL size

**Mouse Thigh sensor, Rat Foot Sensor**

### Physical

|                 |              |
|-----------------|--------------|
| Dimensions      | 16x12x4(h)cm |
| Weight          | 2Kg          |
| Shipping Weight | 5Kg approx.  |
| Packing         | 50x39x17cm   |

**NOTE:** Manufacturer's warranty for MouseOx & accessories is limited to 12 months.

# Metabolic Cages

Cat. No. 41800 / 1 / 2 / 3

## USB Connection

### DESIGNED TO MEASURE:

- FEEDING BEHAVIOUR
- ACTIVITY (OPTIONAL)
- EXCRETORY FUNCTIONS

## General

The recently available animal models (for example, obese and diabetic mice) exhibit symptoms similar to those in humans.

Model organisms are closely monitored, revealing differences, which can be correlated with those of the human disorders in fundamental parameters, as feeding/drinking (quantity & frequency of food/drink uptake), activity (with optional I.R. motion detectors) and excretion (the latter assessed by volume or weight, see following subheading).



## For all types of investigations on METABOLISM, including:

- drugs which produce anorexia
- addiction/aversion to particular substances
- thirst arousing and quenching mechanism
- feeding habits and their modification brought about by environmental conditions

## Basic Cage Design

These carefully engineered metabolic cages are manufactured by TECNIPLAST GAZZADA, see separate datasheet.

All components below the cage floor are removable without upsetting the test animal.

## Feeding and Drinking Analysis

The basic Metabolic Cages are upgraded with the addition of miniature scales, which accurately record ingestion of food.

The trough, shaped as an open box, is made of smooth gnaw-proof plastics. It glides into a receptacle made of stainless steel, fastened to its scale-pan, which senses the load of the pellet food and hence monitors the quantity and frequency of food uptake. The crumbs the animal produces are collected in the front compartment, for a precise evaluation of the food consumption. The water bottles are provided with a spout and rest on a support stud fastened to the scale pan for monitoring quantity and frequency of water uptake. Provision is made to collect any dripping for a more precise water consumption evaluation.

## Activity Detection

The coordinate ambulatory activity and "rearing" of the rodent on test can be measured via the optional motion detector Cat. No. 41700-043, consisting of two facing arrays of emitters and receivers which record beam breaks as the animal moves



## Electronics & Data Acquisition

Each feeding and drinking cage comes with feeding and drinking scales and incorporates a preamplifier module (Cat. 41800-010), which directly connects to the PC USB port for computer processing.

The consumed food and liquid and the optional activity can be recorded directly into a computer at preplanned intervals.

The **Cage Monitoring System (CMS) Software 51800** manages up to 8 cages. The software acquires data and provides results related to partial and total food/liquid consumption and to activity.

**Please refer to the software manual for additional information.**

### Ordering Information

#### BASIC METABOLIC CAGES

|                  |  |
|------------------|--|
| <b>41700-002</b> | Metabolic Cage for rats up to 150 grams  |
| <b>41700-003</b> | Metabolic Cage for mice                  |
| <b>41700-004</b> | Metabolic Cage for rats 150 to 300 grams |
| <b>41700-005</b> | Metabolic Cage for rats over 300 grams   |

#### METABOLIC CAGES WITH FOOD & DRINK RECORDING PROVISION

|              |                                |
|--------------|--------------------------------|
| <b>41800</b> | Cage for rats up to 150 grams  |
| <b>41801</b> | Cage for mice                  |
| <b>41802</b> | Cage for rats 150 to 300 grams |
| <b>41803</b> | Cage for rats over 300 grams   |

Each cage is provided with:-

|                  |  |
|------------------|--|
| <b>E-WP008</b>   | Mains Cable, Europe (or E-WP 008-1 U.K. / E-WP 008-2 U.S.) |
| <b>52010-323</b> | USB connector  |
| <b>E-AU 042</b>  | Individual Power Supply                                    |

#### ACCESSORIES

|                  |   |
|------------------|---|
| <b>41700-043</b> | Combination vertical/horizontal sensors for activity detection in metabolic cages |
| <b>41800-302</b> | Instruction Manual for Hardware   |

#### SOFTWARE

|                  |   |
|------------------|---|
| <b>51800</b>     | Data Acquisition Software for up to 8 cages. For recording of food/liquid consumption and activity in cages series 41800. |
| <b>51800-302</b> | Instruction Manual for the software   |

# Feeding & Activity Analyser

Cat. No. 47552 / 47553 / 47555


## General

The recently available animal models (for example, obese and diabetic mice) exhibit symptoms similar to those in humans. Model organisms are closely monitored, revealing differences, which can be correlated with those of the human disorders in fundamental parameters, as feeding/drinking (quantity & frequency of food/drink uptake), activity (with optional I.R. motion detectors).

The Analyser basically consists of an Animal Cage, which can be provided with optional activity detector and an Electronic Unit.

## USB Connection

TO RECORD THE FEEDING BEHAVIOUR AND ACTIVITY (OPTIONAL) IN RODENTS AND THEIR ALTERATION BROUGHT ABOUT BY A NUMBER OF FACTORS.

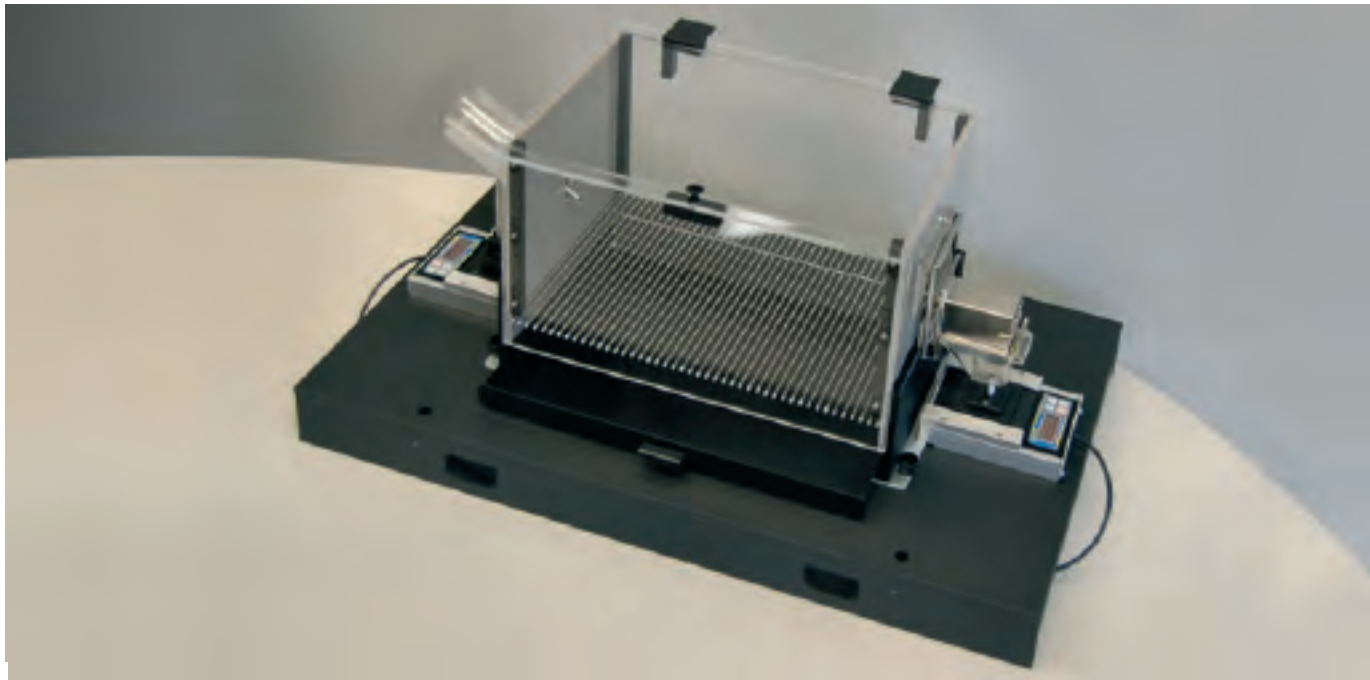

This system has proved to be of great value to quantify at presettable intervals the solid & liquid intake in investigations about:-

- the drugs which produce anorexia
- the thirst arousing and quenching mechanism
- the addiction/aversion to particular substances
- the feeding habits and its modification brought about by environmental conditions.



## Animal Cage

Two types of cages are available: **47552 designed for Rats** and **47453 for Mice**. Both cages have transparent walls and lid to allow the animal observation.

The cages, which are provided with a catch pan, can be easily detached from its Base structure for cleaning or servicing.

They incorporate load cells, which sense the load of food and liquid, thus recording their ingestion, with 0.1 g accuracy and monitoring the frequency of food/liquid uptake.


The **trough**, shaped as an open box, is made of smooth gnaw-proof plastics. It glides into a receptacle made of stainless steel, fastened to its scale-pan. The crumbs the animal produces are collected in the front compartment, for a precise evaluation of the food consumption.

The **water bottles** are provided with a spout and rest on a support stud fastened to the scale pan for monitoring quantity and frequency of water uptake. Provision is made to collect any dripping for a more precise water consumption evaluation.

The inside dimensions of the rat cage are 36x23x18(h) cm, of the mouse cage 26x20x17(h) cm.

## Activity Detection

The coordinate ambulatory activity and "rearing" of

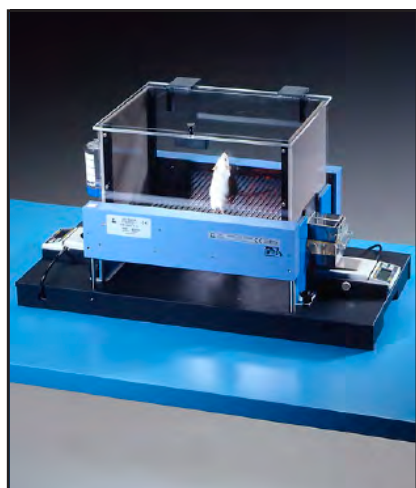

the rodent on test can be measured via the optional motion detector Cat. No. **7435** (horizontal) and **7436** (vertical-rearing), each consisting of two facing blocks of I.R. arrays of emitters and sensors.

For cage **47453** the two sensors are replaced by Cat. **41700-043**, combining both horizontal and vertical sensors.

## FEEDING & DRINKING ANALYSIS

Each cage comes with feeding and drinking scales and incorporates a preamplifier module (Cat. 41800-010), which directly connects to the PC USB port for computer processing.

The consumed food and liquid and the optional activity can be recorded directly into a computer at preplanned intervals.

The Cage Monitoring System (CMS) Software 51800 manages up to 8 cages. The software acquires data and provides results related to partial and total food/liquid consumption and to activity. A PC with USB port is required.

Please refer to the software manual for additional information.

## Mouse Gas Metabolism Cage

The 47555 Gas Metabolism Cage for Mice is an airtight box which has ports to allow monitoring a mouse's O<sub>2</sub> consumption, CO<sub>2</sub> production, VO<sub>2</sub> and VCO<sub>2</sub> via an external Metabolic Monitor (to be purchased separately).

The cages are complete with scales to record food and liquid intake via the 51800 software. Infrared activity arrays 41700-043 may also be attached to monitor the animal's locomotion. The cage complies with the IACUC space requirements for mice.

### Ordering Information

**47552-00 MONITORING CAGE for Food & Drink for RATS up to 150 grams**, complete with food & liquid intake detectors and amplifier module 41800-010.

**Dimensions of the cage: cm 35.4 x 23 x 18 cm**

**47552-002** SAME FOR RATS 150 TO 300 grams

**47552-002** SAME FOR RATS OVER 300 grams

**47453 MOUSE MONITORING CAGE for Food & Drink**, complete with food & liquid intake detectors and preamplifier module 41800-010.

**Dimensions cage: cm 28.6 x 20 x 15.5 cm**

**47555** GAS METABOLISM CAGE FOR MICE

**Each cage is provided with:-**

|           |                         |
|-----------|-------------------------|
| E-WP008   | Mains Cable, Europe     |
| 52010-323 | USB connector           |
| E-AU 042  | Individual Power Supply |

### ACCESSORIES

**7435** Set of horizontal activity sensors for 47552

**7436** Set of vertical (rearing) sensors for 47552

**41700-043** Combination vertical/horizontal activity sensors for 47553/47555

**47552-302** Instruction Manual for the Hardware

### SOFTWARE

**51800** Data Acquisition Software for up to 8 cages.

For recording of food/liquid consumption and activity in cages series 41800.

**51800-302** Instruction Manual for the software

## New Microwave Brain Fixation System

Cat. MMW-05 (5kW)

### General

In neurochemical studies of the brain, it is of great importance to measure accurately neurochemical events *in vivo*.

However, it is difficult to perform reproducible measurement of these events because rapid post-mortem changes occur in the brain concentrations of metabolites and neurotransmitters.

With the NEW Microwave Brain Fixation System by Muromachi, a living mouse or rat is positioned inside the applicator and, in less than 1 second, the microwave beam stops all brain chemistry at the level present in the living animal.

**Measuring brain chemistry *in-vivo* is possible!**



**THE FASTEST AND MOST EFFECTIVE KNOWN METHOD OF HALTING BRAIN CHEMICAL ACTIVITY**

**brain fixation occurs in 1 second**

**activity of degrading enzymes is blocked**

### Prior to analysis of:

- Phosphorylated proteins
- Acetylcholine, Serotonin, Endorphins
- Prostaglandins, Catecholamines
- C-AMP, C-GMP, GABA, DOPA

### NEW features:

- Improved usability - touch screen
- Air-cooled (no water circulation)
- CE-certified
- Absolute safety - negligible leakage

Various techniques have been developed to **prevent post-mortem changes**. One of the more common method is cooling or freezing by immersion of the decapitated head in liquid Nitrogen or cooled Freon to **inactivate enzymes** involved in the metabolism of these compounds.

**Cooling is not fully effective in preventing post-mortem changes** as the time required to freeze deep structure of the brain may range from 10 - 90 seconds; post mortem changes will occur during this period.

An alternate method is microwave heating to inactivate enzymes.

**The microwave method has several advantages over cooling or freezing:**


- The enzymes in the whole brain can be completely inactivated in a very short time
- The brain can be dissected easily and reproducibly at room temperature

Microwave fixation system must satisfy the following criteria:

1. elevate the temperature of brain up to 75-90°C as rapidly as possible, by effectively focusing microwave energy on the animal head
2. give the same results from animal to animal
3. be easily and safely used by personnel not experienced in microwave

## Instrument Description

Thanks to Patented Microwave Focus Applicators, microwaves are channeled and focused by the wave guide from above the head, rather than in front. The entire animal head is placed in a uniform microwave field. Movements of the head do not change the field strength or microwave distribution.

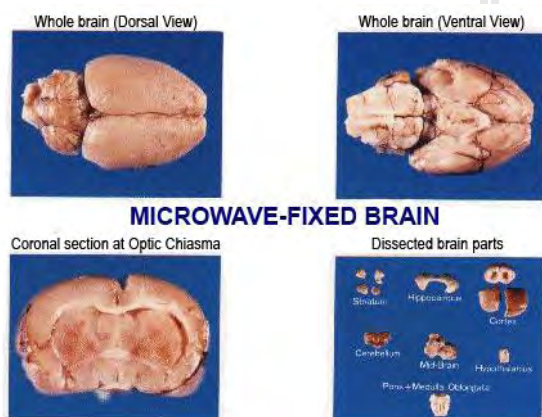


Uniform fixation of the whole brain is thus achieved, without "over cooking" and damage to the hypothalamus.

Muromachi's unique applicators provide protection to the researcher and also compensate for individual differences between animals, giving more reproducible results.

The Muromachi Microwave Fixation Systems are safely designed, so that the microwave leakage will not exceed 1 mW/cm<sup>2</sup> (well below the safety standards).

The Microwave Fixation System comes with specific applicator heads and water-jacketed animal holders:

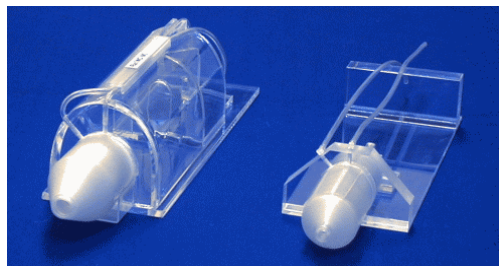

## Ordering Information

**MMW-05** **Microwave Fixation System 5KW**, including 1 Applicator head and 1 animal holder, to be selected

### Applicator heads

**TAW-174P** for mouse holder  
**TAW-424SP** for rat holder WJR-S  
**TAW-424MP** for rat holder WJR-M & L

### Water-Jacketed Animal Holders

**WJM-24** for mice 15-20g  
**WJM-28** for Mice 20-40g  
**WJR-S** for Rats 150-250g  
**WJR-M** for Rats 250-400g  
**WJR-L** for Rats 400-500g

**TAW-174P** for mouse holder  
**TAW-424SP** for rat holder WJR-S  
**TAW-424MP** for rat holder WJR-M & L

### Physical

Power 190/240VAC 30A or 380-440VAC 20A 3-phase  
 Dimensions 75(w)x55(d)x128(h)cm  
 Weight 103Kg  
 Shipping weight 195Kg  
 Packing 81x100x132cm

## Bibliography

- B. Sahin et alia: "Evaluation of neuronal phosphoproteins as effectors of caffeine and mediators of striatal adenosine A2A receptor signaling" *Brain Research* 1120: 1-14, 2000
- P. Svenningsson et alia: "DARPP-32 mediates serotonergic neurotransmission in the forebrain" *PNAS* 99 (5, 2002
- G.L. Caporaso et alia: "Drugs of abuse modulate the phosphorylation of ARPP-21, a cyclic AMP-regulated phosphoprotein enriched in the basal ganglia" *Neuropharmacology* 39:1637-1644, 2000
- A. Nishi et al.: "Amplification of dopaminergic signaling by a positive feedback loop" *PNAS Early Edition* 1-6, 2000

**1963 2013**



For the past 5 decades we have provided scientists with the unmatched tools  
necessary to transform their ideas into meaningful research and results  
We look forward to working with you and to **another 50 years.**



latest revision

05/02/2015